JAY Positioning Supports **EN**

Important Consumer Information.

NOTE: This manual contains important instructions that must be passed on to the user of this product. Please do not remove this manual before delivery to the end user.

SUPPLIER: This manual must be given to the user of this product.

USER: Before using this product, read this entire manual and save for future reference.

WARRANTY REGISTRATION: To validate the warranty on this product, please complete and return the postcard attached to the enclosed instruction and warranty booklet.

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
Bruksanvisning

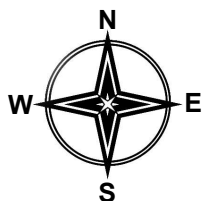ENGLISH

WARNINGS	4
ASSESSMENT AND SIZING	5
INSTALLATION	6
WARRANTY	21

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SHIPPING:

Check all parts are free from shipping damage. In case of damage DO NOT use. Contact Sunrise Medical or your supplier for further instructions.

WARNING!

The products described in this manual were designed to be fitted and installed exclusively by:

- An authorized Sunrise Medical supplier.
- A qualified and competent healthcare professional who has undergone appropriate certified training by either Sunrise Medical or their approved supplier.

THE JAY POSITIONING SUPPORTS INTENDED USE:

The Jay Positioning Supports are intended to provide additional support for a wide variety of wheelchair users.

WARNING!

Any application other than positioning a person in a wheelchair where its failure could result in injury is a misuse of these products and it is not safe.

Maximum User Weight.

Because the harnesses can be installed onto a wide range of wheelchairs, the harness maximum user weight is equal to the maximum user weight of the wheelchair to which it is fitted.

Intended Lifetime.

The intended lifetime of this product is 5 years.

WARNING!

Positioning belts and trunk supports should only be used to help support the user's posture.

DANGER!

- Risk of chest compression or suffocation!
- Do not use any trunk support or chest harness without stabilizing the pelvis.
- Ensure that the user cannot slide down out of the seating system.

NOTE: The use of a JAY contoured cushion may help prevent the user from sliding in the seat.

WARNING!

- Any belt, harness or support strap must not interfere with breathing.
- Any belt, harness or support strap may alter the person's seating position.
- It is crucial that users at risk of skin breakdown are regularly checked for skin integrity.
- Employ an appropriate pressure relief regime.
- A risk assessment is recommended to determine If the user can perform pressure relief activities independently, if not a caregiver/clinician must carry out this task.

WARNINGS FOR TRANSPORTATION

- Jay Positioning belts, harnesses and shoulder straps were designed for positioning, not to be used as transportation safety belt!
- JAY POSITIONING BELTS, HARNESSSES, AND SHOULDER STRAPS may be used while in transit, but are not crash tested or approved to provide protection during a crash.
- When traveling in a motor vehicle, only use these positioning supports in combination with restraints that have been transit tested in accordance with SAE J2249 and/or ANSI/RESNA WC/20.

HYGIENE MEASURES WHEN BEING RE-USED:

Prior to the wheelchair being re-used, it must be carefully prepared. All surfaces which come into contact with the user must be treated with a disinfection spray.

To do this, you must use a disinfectant as authorised/ recommended in your country, for rapid alcohol-based disinfection for medical products and medical devices, which must be disinfected quickly.

Please be aware of the manufacturer's instructions for the disinfectant you are using.

In general, a complete disinfection cannot be guaranteed on seams. We therefore recommend that you dispose of seat and back slings to avoid micro-bacterial contamination with active agents according to your local infection protection law.

CONDITIONS FOR SAFE STORAGE

- Products should be stored in a dry place: at 20 – 75% Relative Humidity.
- Protect from direct sunlight and dust for long periods.
- Away from direct heat and at a temperature between 5°C (40°F) to 40°C (104°F)
- Position and stack the packaged products according to the "this side up" arrows on the packaging.
- If a package is high and narrow and therefore potentially unstable, ensure it is well secured on the pallet or in the racking, to prevent it from falling down.
- Do not stack anything on top of the packaged products when the "do not stack" mark on the packaging is shown.
- Product should never be stored outside and subject to the elements.

Word	Definition
 DANGER!	Advice to the user of potential risk of serious injury or death if the advice is not followed.
 WARNING!	Advice to the user of a risk of injury if the advice is not followed.
 CAUTION!	Advice to user that potential damage to equipment may occur if the advice is not followed.
NOTE:	General advice or best practice.
	Reference To Additional Documentation


ASSESSMENT:**⚠ WARNING!**

An appropriate clinical assessment by a clinician or qualified wheelchair professional is highly recommended to guarantee safety and correct clinical usage.

CLIENT MEASUREMENTS:

Fig.1.0 shows a series of standard measurements.

It is the clinician/wheelchair professionals responsibility to make the appropriate measurements relevant to the product to be used. This is essential to ensure client safety and that the correct size is ordered.

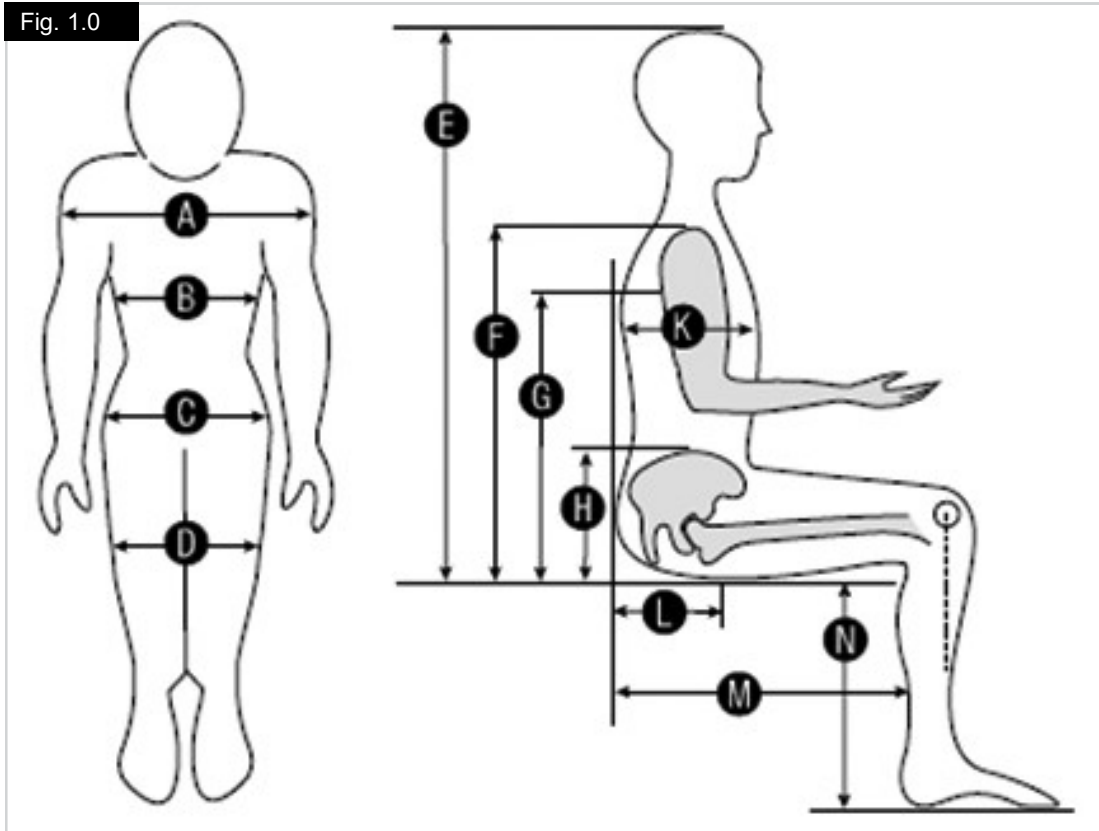Fig. 1.0

- A. Shoulder Width.
- B. Chest Width.
- C. Hip Width.
- D. Width at the Knee.
- E. Seat to Top of Head.
- F. Seat to Top of Shoulder.
- G. Seat to Axilla.
- H. Seat to Psois.
- I. Not Used
- J. Not Used
- K. Chest Depth.
- L. Back to Anterior of Seat Bones, (Ischial Tuberosity)
- M. Seat Depth.
- N. Seat to Footplate.

INSTALLATION:

When assessing posture in terms of seating, always begin with the pelvis.

PELVIC POSITIONING:

Pelvic Area Shaded, (Fig.2.0).

Fig. 2.0**Two Point Pelvic Positioning Belts.**

- Two Point Belt Push Button, (Fig.3.0).
- Two Point Belt Side Release, (Fig.4.0).

Four Point Pelvic Positioning Belts.

- Four Point Belt Push Button, (Fig.5.0).
- V Contour Belt Push Button, (Fig.6.0).

J2 DC CUSHION COMPONENTS**Fig. 3.0****Fig. 4.0****Fig. 5.0****Fig. 6.0**

PELVIC BELT INSTALLATION:

Pelvic belts can be installed in different ways according to user's seating system, wheelchair type and hardware used.

Belt Anchors.

All JAY positioning belts include belt anchors for attachment, (Fig.7.0).

⚠ WARNING!

The products described below must be fitted and installed exclusively by:

- An authorized Sunrise Medical supplier.
- A qualified and competent healthcare professional who has undergone appropriate certified training by either Sunrise Medical or their approved supplier.
- Never place a pelvic belt, or positioning support across the waist or abdomen.

Initial Considerations.

Angles of attachment.

The angle that the pelvic belt is mounted to the wheelchair has a direct impact on the pull angle on the pelvis.

A 60° attachment angle to the seat rail positions the pelvic belt anterior and inferior to the axis, which reduces the tendency of the user with a posterior pelvic tilt to slide under the belt.

A 45° attachment angle on the back post encourages a posterior pelvic tilt, which assists in positioning users with an anterior pelvic tilt.

INSTALLATION PROCEDURE

Two Point Pelvic Belt, (Fig.9.0).

1. Locate desired attachment point.
2. Attach the belt hardware to the wheelchair frame.
3. Consider the user's position to determine the ideal angle for attachment and position the belt on the centre of the pelvis.
4. Thread the webbing end through the belt mount.

Fig. 7.0

Fig. 8.0

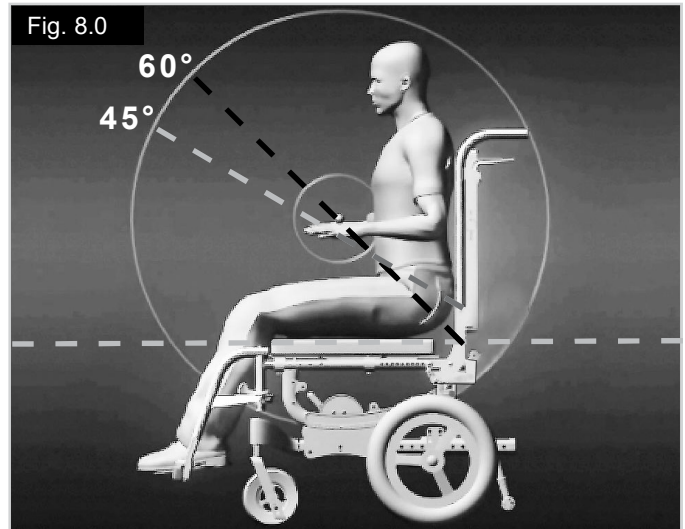
Fig. 9.0

Four Point Pelvic Belt, (Fig.10.0).

- A. Main Straps, (Fig.11-12).
 - B. Secondary Straps, (Fig.11-12).
3. Follow the installation steps described in the previous Two Point Pelvic Belt paragraph, to attach the main straps of the pelvic belt.
 4. Attach the secondary straps to the chair so that they anchor the primary straps.
 5. Tighten secondary straps to prevent movement of the primary straps.

SAFETY:

For both types of belt, when the installation is finished, carry out the following procedure to ensure that the user is appropriately and safely positioned:

1. Check if the user is sitting comfortably.
2. Check if there are any areas of irritation. This usually manifests as a reddening of the skin especially in soft tissue.
3. Check that the user does not slide down through the belt, if they do, adjust the belt as required. If the user still slides under the belt and the belt is a two point model, asses if there is a need for the four point belt, this will provide additional security and adjustment.
4. Check if the pelvic belt mounting hardware is interfering with any other positioning accessories/devices on the wheelchair and reposition it if necessary.


ANTERIOR TRUNK SUPPORTS INSTALLATION:

* Measurements = mm and (inches). Measurements may vary slightly due to binding finish

Fig. 13.0

Classic Style.

Classic Dimensions*								
	A	B	C	D	E	F	G	H
Early Intervention	248 (9.75)	76 (3)	51 (2)	133 (5.25)	76 (3)	76 (3)	25 (1)	165 (6.5)
X-S	318 (12.5)	89 (3.5)	51 (2)	191 (7.5)	89 (3.5)	76 (3)	51 (2)	191 (7.5)
S	381 (15)	114 (4.5)	51 (2)	229 (9)	102 (4)	102 (4)	51 (2)	216 (8.5)
M	482 (19)	140 (5.5)	51 (2)	280 (11)	127 (5)	127 (5)	76 (3)	241 (9.5)
L	572 (22.5)	165 (6.5)	51 (2)	330 (13)	152 (6)	152 (6)	89 (3.5)	267 (10.5)
XL	572 (22.5)	165 (6.5)	64 (2.5)	356 (14)	178 (7)	178 (7)	114 (4.5)	292 (11.5)

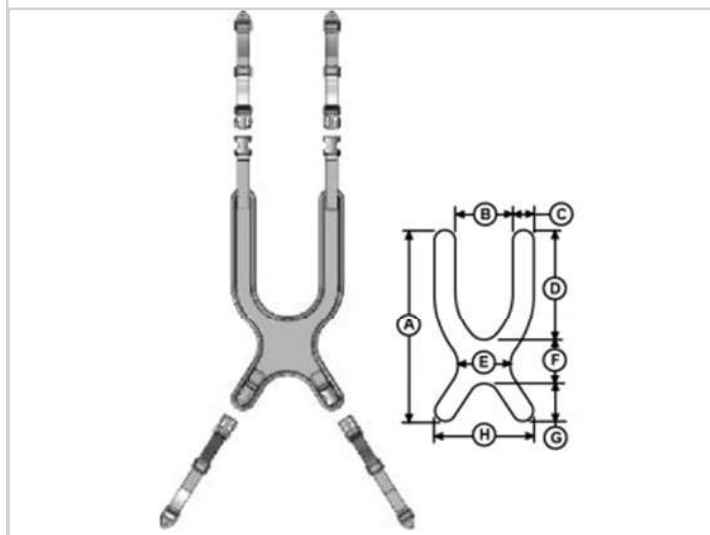**Strap Positioner, (Fig.14.0).**

This accessory is sold separately.

Fig. 14.0


* Measurements = mm and (inches). Measurements may vary slightly due to binding finish

Fig. 15.0

Contour Style.

Contour Dimensions*								
	A	B	C	D	E	F	G	H
Early Intervention	248 (9.75)	76 (3)	45 (1.75)	133 (5.25)	25 (1)	89 (3.5)	38 (1.5)	165 (6.5)
X-S	318 (12.5)	89 (3.5)	51 (2)	191 (7.5)	25 (1)	89 (3.5)	38 (1.5)	191 (7.5)
S	381 (15)	114 (4.5)	51 (2)	229 (9)	25 (1)	102 (4)	51 (2)	216 (8.5)
M	470 (18.5)	140 (5.5)	51 (2)	280 (11)	25 (1)	127 (5)	64 (2.5)	241 (9.5)
L	572 (22.5)	165 (6.5)	51 (2)	330 (13)	38 (1.5)	152 (6)	89 (3.5)	267 (10.5)
XL	572 (22.5)	165 (6.5)	64 (2.5)	356 (14)	38 (1.5)	178 (7)	114 (4.5)	292 (11.5)

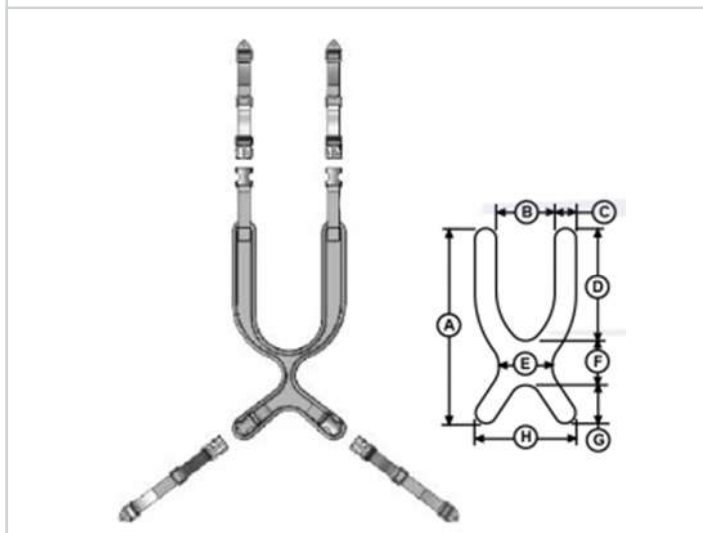
* Measurements = mm and (inches). Measurements may vary slightly due to binding finish

Fig. 16.0

Back Pack Style.

Back Pack Dimensions*		
	A	B
S 305 (12)	305 (12)	51 (2)
L 457 (18)	457 (18)	51 (2)

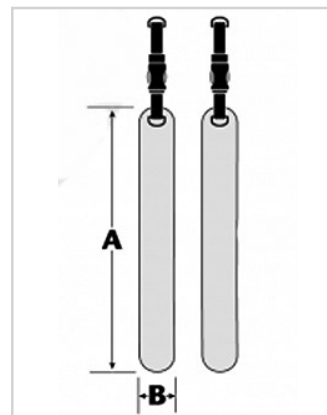
* Measurements = mm and (inches). Measurements may vary slightly due to binding finish

Fig. 17.0

Centre Opening Style.

Centre Opening Dimensions*			
	A	B	C
Early Intervention	248 (9.75)	178 (7)	127 (5)
X-S	280 (11)	203 (8)	147 (5.75)
S	311 (12.25)	229 (9)	165 (6.5)
M	362 (14.25)	280 (11)	248 (9.75)
L	394 (15.5)	292 (11.5)	248 (9.75)
XL	445 (17.5)	324 (12.75)	254 (10.0)

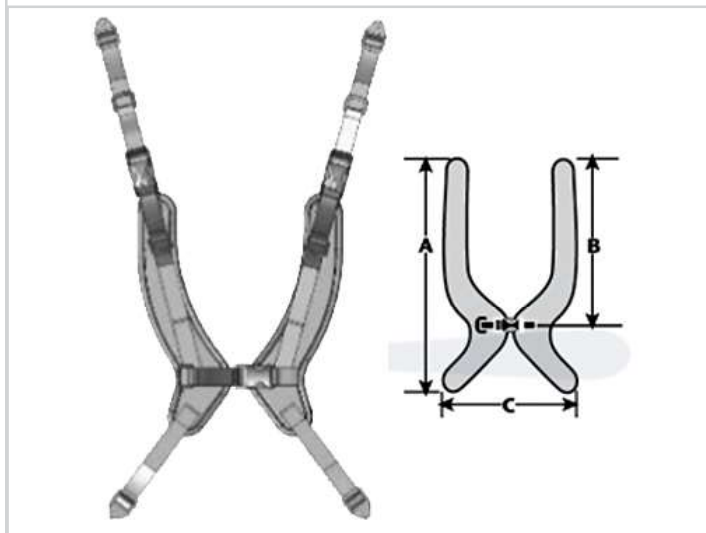Technical Notes:


- Adjustable buckle closure provides quick user access.
- Available in structured construction only.
- Four point attachment provides independent tension adjustment.

* Measurements = mm and (inches). Measurements may vary slightly due to binding finish

Fig. 18.0

Zipper Opening Style.

Zipper Opening Dimensions*								
	A	B	C	D	E	F	G	H
S	381 (15.0)	114 (4.5)	51 (2)	229 (9)	102 (4)	102 (4)	51 (2)	216 (8.5)
M	482 (19.0)	140 (5.5)	51 (2)	280 (11)	127 (5)	127 (5)	76 (3)	241 (9.5)
L	572 (22.5)	165 (6.5)	51 (2)	330 (13)	152 (6)	152 (6)	89 (3.5)	267 (10.5)
XL	572 (25.5)	165 (6.5)	64 (2.5)	356 (14)	178 (7)	178 (7)	114 (4.5)	292 (11.5)

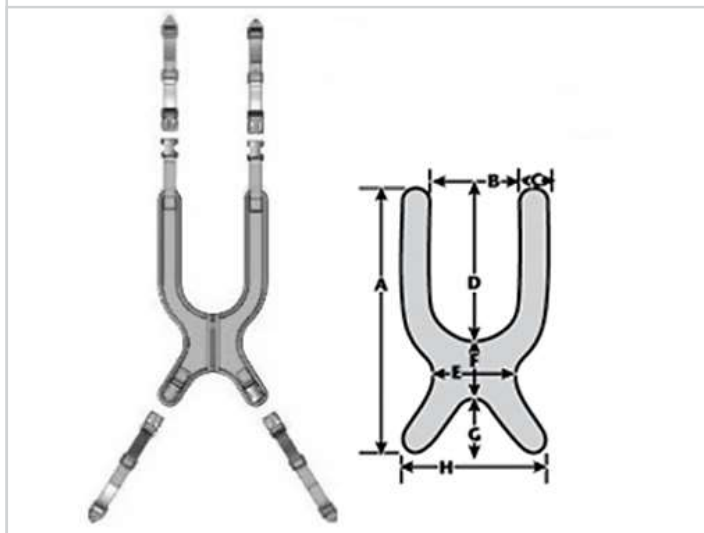Technical Notes:


- Zipper closure provides quick user access.
- Breathable fabric and machine washable.
- Dynamic construction provides support and allows limited range of movement.

* Measurements = mm and (inches). Measurements may vary slightly due to binding finish

Fig. 19.0

Trunk Strap Style.

Trunk Strap Dimensions*		
	A	B
X-S	76 (3)	203 (8)
S	76 (3)	254 (10.0)
M	76 (3)	305 (12)
L	102 (4)	356 (14)
XL	102 (4)	406 (16)

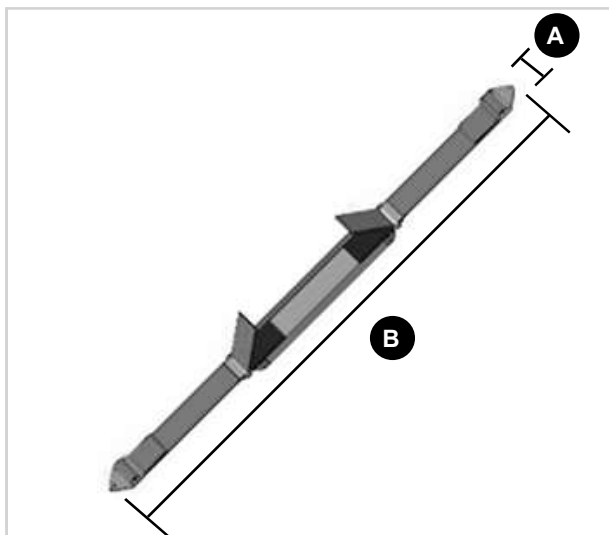Technical Notes:



- Mounts to backrest with strap anchors.
- Anterior Velcro and "D" ring closure provide quick user access.

* Measurements = mm and (inches). Measurements may vary slightly due to binding finish

Fig. 20.0

"Y" Style Shoulder Straps.

"Y" Style Dimensions*			
	A	B	C
Early Intervention	381 (15.0)	178 (7)	45 (1.75)
S	482 (19.0)	203 (8)	51 (2)
M	559 (22)	229 (9)	57 (2.25)
L	635 (25)	280 (11)	64 (2.5)

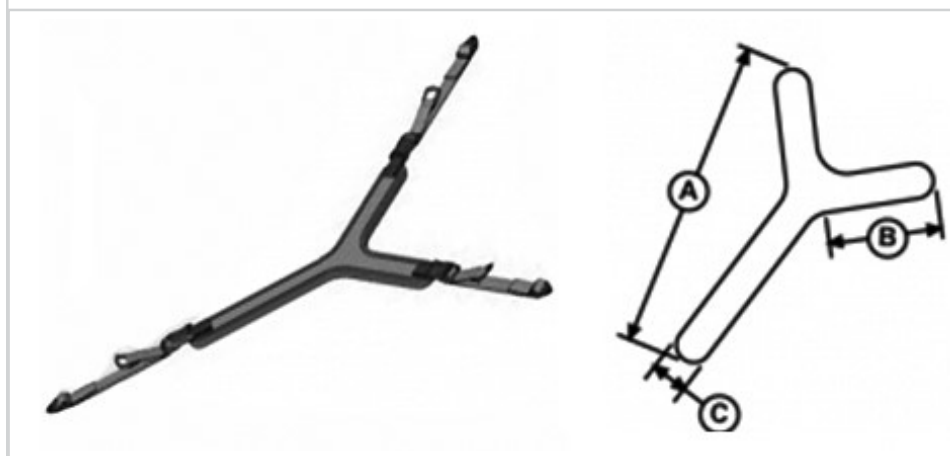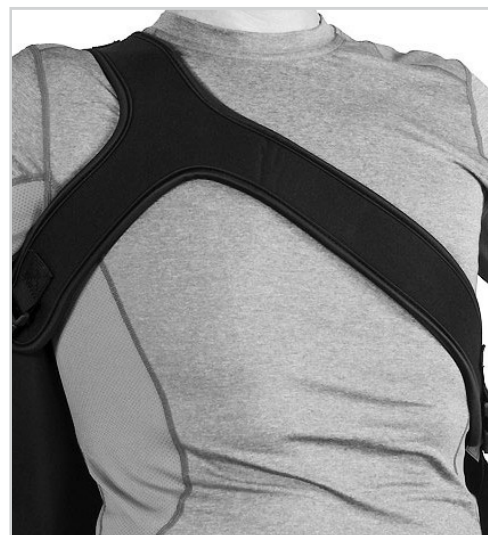INSTALLATION REQUIREMENTS

To be able to mount any anterior trunk support / shoulder harness you have to ensure that the following are in place:

- A wheelchair with a solid backrest set approximately at shoulder height.
- An installed pelvic belt.

! WARNING!

The products described below must be fitted and installed exclusively by:

- An authorized Sunrise Medical supplier.
- A qualified and competent healthcare professional who has undergone appropriate certified training by either Sunrise Medical or their approved supplier.
- Never use anterior trunk supports without stabilizing the pelvis.
- Always use a properly installed pelvic belt.

Installation steps

1. Fitting and installation should be performed with the user sitting in the wheelchair.
2. If the user is not using any pelvic belt: Seat the user in the wheelchair and install the pelvic belt according to the previous instructions of this manual, (pages 9 and 10). Define which shoulder harness type is suitable for the user.
3. The following types are available:
 - i. Classic Style – for optimal distribution of support and maximum stability.
 - ii. Contour Style - for women or individuals who want to limit the amount of contact area between the harness and the body.
 - iii. Zipper Open - to aid in frequent transfers and ease of removal.
 - iv. Adjustable Stretch, Centre Opening - to aid in frequent transfers and provide a more precise fit.
 - v. Trunk Strap - to increase range of motion when shoulder rotation and reaching are required.
 - vi. Y-Style - for increased range of motion for one arm while providing optimal support.
 - vii. Backpack Style – for maximum shoulder retraction and minimal chest contact.
4. Place the desired anterior trunk support /shoulder harness on chest with upper straps over the backrest.
5. Position straps between neck and shoulder joint.

Top Straps

- i. Adjust each shoulder strap for comfortable fit according to the user's shoulder width and sitting position.
- ii. Attach to the solid backrest using the fasteners/strap positioners, (sold separately from harness).

Lower Straps

- iii. Identify optimal position for the lower straps on the backrest without interfering with other positioning devices (lateral supports, etc.).
- iv. Install end mounting between the top of the pelvis and the seat, (Fig.21.0).
- v. Attach using the fasteners supplied.

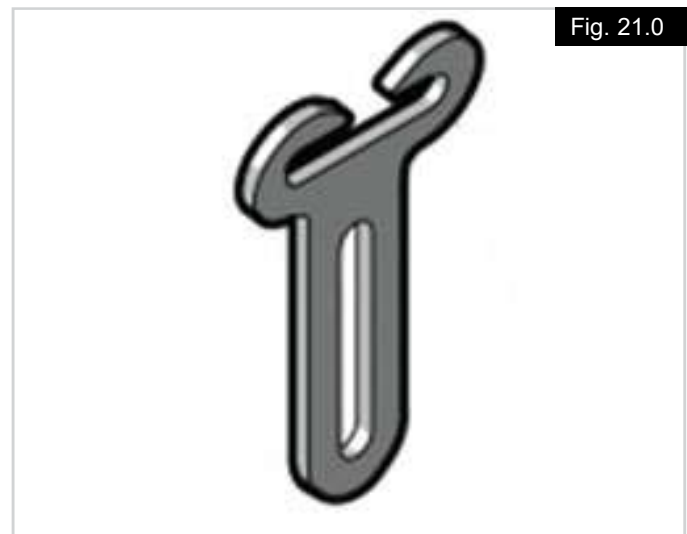
Fig. 21.0

SAFETY CHECK

Once the installation is finished, carry out the following checks to ensure that the user is appropriately and safely positioned.

If possible, ask the user to lean forward and to the sides to test the fit.

Check that:

1. The buckles and adjustment straps are operating properly.
 2. The user sits comfortably in the seating and the harness.
 3. The harness is not too low or too high. The support must not rub or cut into the users neck.
 4. The mounting hardware is not interfering with any other positioning accessories/devices on the wheelchair.
 5. Adjust the harness if necessary.
- Harness Too High, (Fig.22.0).
 - Harness Too Low, (Fig.23.0).
 - Harness Fitted Correctly, (Fig.24.0).

Fig. 22.0

Fig. 23.0

Fig. 24.0

INSTALLATION OF Y-STYLE HARNESS

The Y-Style harness is typically used when increased range of motion for one arm is needed, (one shoulder rotation).

The JAY Y-Style positioning harness provides control of trunk rotation and independent shoulder retraction while minimizing interferences with feeding tubes.

Installation steps.

1. Place the upper part of the Y-Strap above and below the rotated shoulder and install belt anchors to the backrest's rear face, (Fig.25.0).
2. Stretch the lower strap diagonally across the trunk and attach the lower part of the Y-Strap to the backrest's rear face, (Fig.26.0).
3. Using the tension straps on the top portion of the Y-Strap, pull the rotated shoulder back until a proper and comfortable position for the user is achieved, (Fig.27.0).

Fig. 25.0

Fig. 26.0

Fig. 27.0

LOWER BODY POSITIONING SUPPORTS

Shoe holder With Tendon Relief, (Fig.28.0 & Fig.30.0).

Shoe Holders help maintain positioning of the feet and ankles, broadening the seating “footprint”, and improving overall posture.

* Measurements = mm and (inches).

Fig.30	A*	B*
S	89 (3.5)	133 (5.25)
M	102 (4.0)	171 (6.75)
L	121 (4.75)	203 (8.0)
XL	127 (5.0)	254 (10.0)

Foot and Ankle Positioner, (Fig.29.0 & Fig.31.0).

Foot and ankle positioners keep the feet protected on the foot plates without the restrained feeling of a rigid shoe holder.

The foot and ankle positioner has three-point ankle attachments which form a cup for the heel while the side straps wrap around ankle.

Ankle positioners are fitted with side-release buckles that aid access for easy adjustment and release for transfers.

* Measurements = mm and (inches).

Fig.31	A*	B*	C*	Ankle Circumference
S	127 (5.0)	133 (5.25)	114 (4.5)	165 - 229 (6.5 - 9.0)
M	127 (5.0)	171 (6.75)	114 (4.5)	203 - 254 (8.0 - 10.0)
L	152 (6.0)	203 (8.0)	140 (5.5)	241 - 305 (9.5 - 12.0)
XL	165 (6.5)	254 (10.0)	152 (6.0)	279 - 343 (11.0 - 13.5)

Fig. 28.0

Fig. 29.0

Fig. 30.0


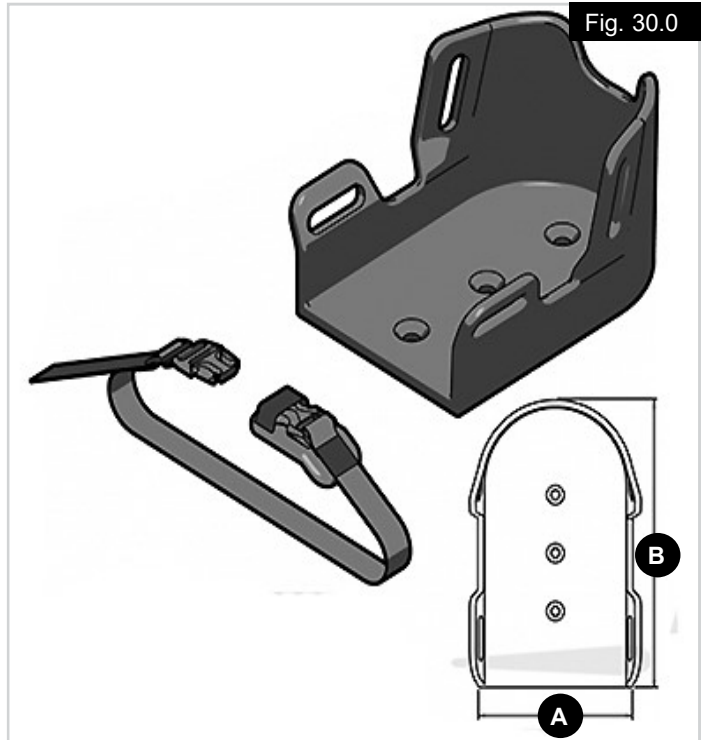
Fig. 31.0

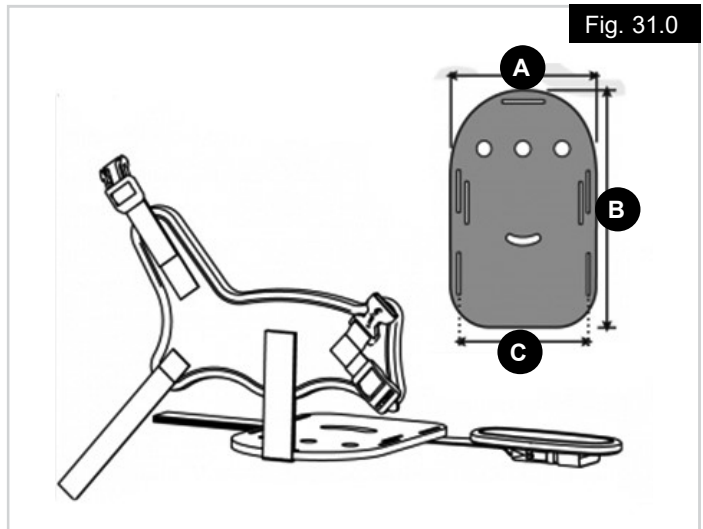
SHOE HOLDER WITH PADDED STRAPS

The shoe holder with padded straps keep the feet protected on the foot plates, but also provide an added degree of comfort.

The shoe holder with padded straps includes two, two-point foot support that provides extra security and comfort.

The padded straps are fitted with side-release buckles that aid access for easy adjustment and release for transfers, (Fig.32.0).

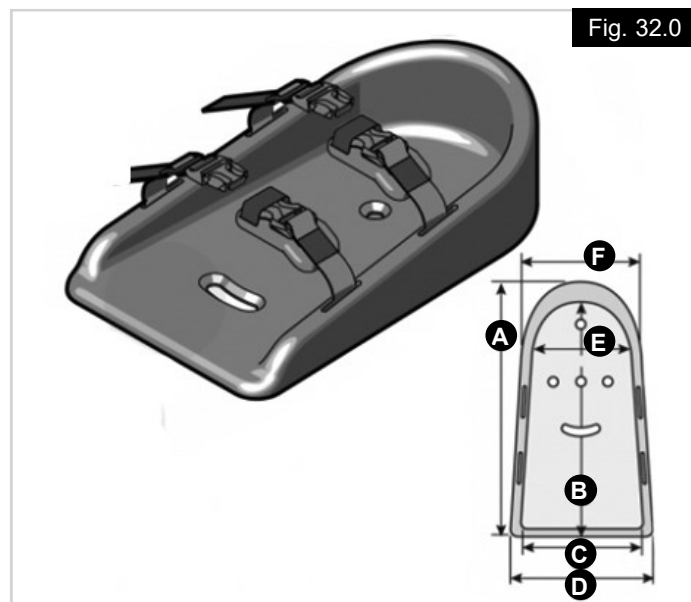* Measurements = mm and (inches).

Fig.32	A*	B*	C*	D*	E*	F*
S	89 (5.75)	133 (5.25)	83 (3.25)	111 (4.375)	83 (3.25)	108 (4.25)
M	102 (7.75)	171 (7.25)	92 (3.625)	121 (4.75)	83 (3.25)	108 (4.25)
L	121 (9.75)	203 (9.25)	114 (4.5)	146 (5.75)	102 (4.0)	133 (5.25)
XL	127 (12.25)	254 (11.75)	152 (6.0)	178 (7.0)	133 (5.25)	159 (6.25)

SHOE HOLDER FLEXIBLE

The flexible shoe holder with fold over padded cover keeps the feet protected on the foot plates. The flexible shoe cover also provides an added degree of comfort and flexibility.

The flexible shoe holder includes two, two-point fold over straps that provide extra security and comfort and flexibility.

The straps are fitted with side-release buckles that aid access for easy adjustment and release for transfers, (Fig.33.0).

FOOT BOXES

Foot Boxes reduce the risk of skin breakdown and/or client injury that may result from periods of high tone or continual contact with the hangers and foot plates.

Single Foot Box Dimensions

* Measurements = mm and (inches).

A	B
203 (8.0)	203 (8.0)
305 (12.0)	203 (8.0)

Full Foot Box & Divider Dimensions

* Measurements = mm and (inches).

A	B
203 (8.0)	406 (16.0)
203 (8.0)	457 (18.0)
203 (8.0)	508 (20.0)
203 (8.0)	559 (22.0)
305 (12.0)	406 (16.0)
305 (12.0)	457 (18.0)
305 (12.0)	508 (20.0)
305 (12.0)	559 (22.0)

Full Foot Box Dimensions

* Measurements = mm and (inches).

A	B
203 (8.0)	305 (12.0)
203 (8.0)	356 (14.0)
203 (8.0)	406 (16.0)
203 (8.0)	457 (18.0)
203 (8.0)	508 (20.0)
203 (8.0)	559 (22.0)
305 (12.0)	406 (16.0)
305 (12.0)	457 (18.0)
305 (12.0)	508 (20.0)
305 (12.0)	559 (22.0)

Fig. 34.0

Fig. 35.0

Fig. 36.0

SINGLE FOOT BOX INSTALLATION

Mounting a Single Foot box onto foot plate:

1. Fasten the Single Foot box onto the foot plate (A) using screws (B), as shown in Fig.37.0.
2. If the existing foot plate holes cannot accommodate the holes in the Single Foot box, please use the drill guide, shown in Fig.39.0.
3. Place Velcro strap (C) around leg rest extension bar (D).

FULL FOOT BOX INSTALLATION (with or without divider)

Mounting a Full Foot box onto foot plate

1. Fasten the Full Foot box onto the foot plate (A) using screws (B), as shown in Fig.38.0.
2. If the existing foot plate holes cannot accommodate the holes in the Full Foot box, please use the drill guide, shown in Fig.40.0.
3. Place Velcro strap (C) around leg rest extension bar (D).





DRILL GUIDES

Single Foot Plate, Fig.39.0.

- A. 102 mm (4.0 Inches).
- B. Diameter = 7.1 mm (9/32 Inch).

Double (Full), Foot Plate, Fig.40.0.

- A. 102 mm (4.0 Inches).
 - B. Diameter = 7.1 mm (9/32 Inch).
- 1 = Left Foot Plate, 2 = Right Foot Plate.

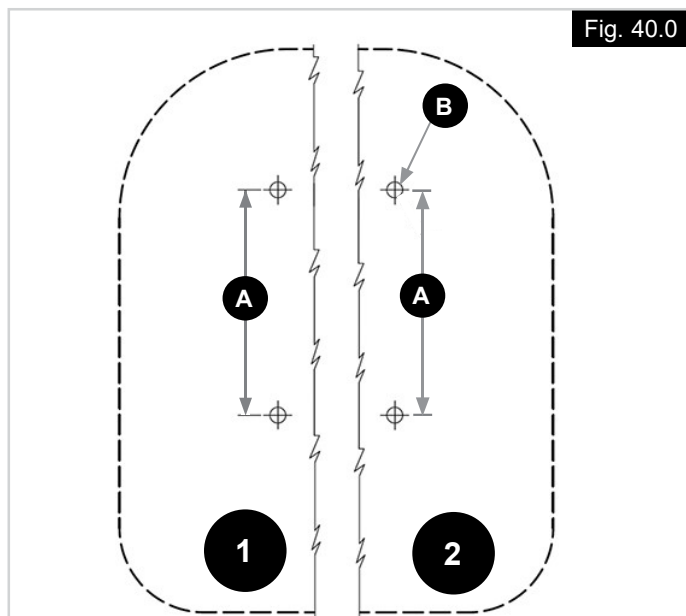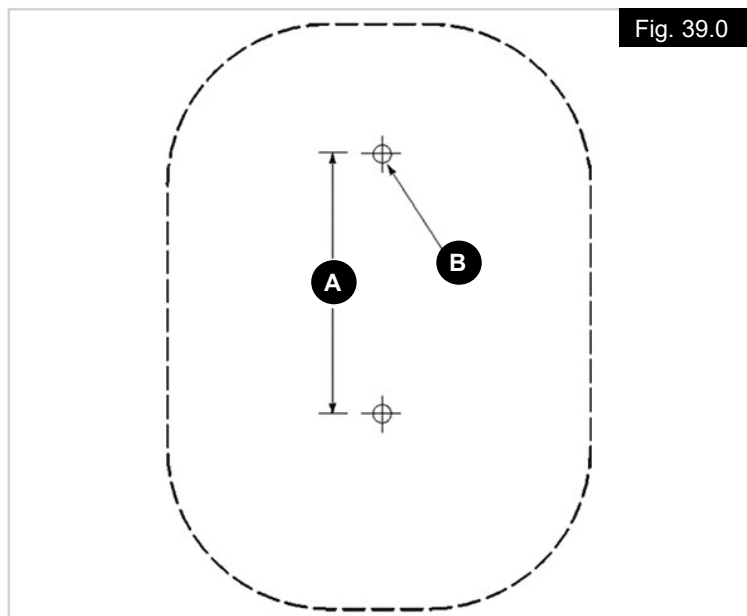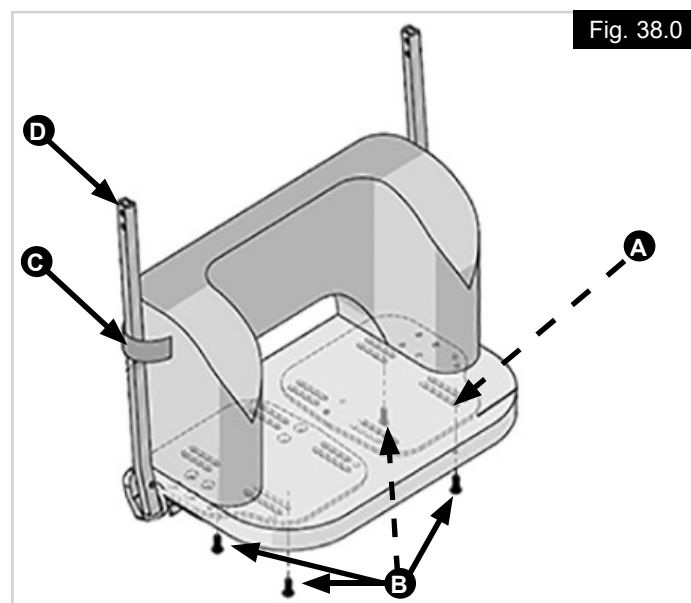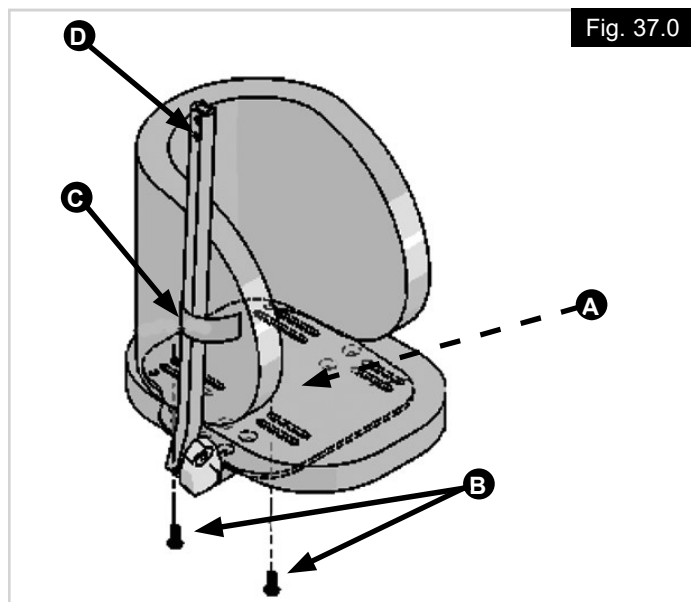24 MONTH LIMITED WARRANTY

Each Product is carefully inspected and tested to provide peak performance. Every Product is guaranteed to be free from defects in materials and workmanship for a period of 24 months from the date of purchase, provided normal use. Should a defect in materials or workmanship occur within 24 months from the original date of purchase, Sunrise Medical will, at its option, repair or replace it without charge. This warranty does not apply to punctures, tears or burns, nor to any removable cover.

Removable covers are guaranteed to be free from defects in materials and workmanship for a period of six (6) months from the date of purchase, provided normal use. Should a defect in materials or workmanship occur within six (6) months from the original date of purchase, Sunrise Medical will, at its option, repair or replace it without charge.

Claims and repairs should be processed through the nearest authorised supplier. Except for express warranties made herein, all other warranties, including implied warranties of merchantability and warranties of fitness for a particular purpose are excluded.

Type:	Product Name/ SKU Number		This symbol means Medical Device
	Lot		Manufacturer's address
	Date of Manufacture		European Authorised Representative
	CE Mark		UK Responsible Person
	UKCA Mark		Swiss Representative's address
	Consult instructions for use		Importer's address

The management system of SUNRISE MEDICAL is certified to EN ISO 13485 and ISO 14001.

We at SUNRISE MEDICAL have been awarded the ISO-13485 certificate, which affirms the quality of our products at every stage, from R & D to production. This product complies with the standards set forth in EU and UK regulations.

Indications

The varieties of fitting variants, as well as the modular design, mean that it can be used by those who cannot walk or have limited mobility because of:

- Joint contractures
- Muscle Dystrophy
- Neurological disease pattern
- Spasticity
- Assymetric Postures
- Hypo- Hyper tonic
- Cerebral Palsy
- MS
- SCI

Contraindication

No contraindications are known or available.

NOTE: General user advice. Not following these instructions may result in physical injury, damage to the product or damage to the environment!

A notice to the user and/or patient that any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

As the manufacturer, SUNRISE MEDICAL, declares that this product conforms to the Medical Device Regulation (2017/745).

As the manufacturer, SUNRISE MEDICAL, declares that the product conforms to the UK Medical Devices Regulation 2002 No. 618.

B4Me special adaptations

Sunrise Medical strongly recommends that in order to ensure that your B4Me product operates, and performs as intended by the manufacturer; all the user information supplied with your B4Me product is read and understood, before the product is first used. Sunrise Medical also recommends that the user information is not discarded after reading it, but it is kept safely stored for future reference.

Medical Device Combinations

It may be possible to combine this Medical device with one or more other Medical Device or other product. Information on which combinations are possible can be found at www.Sunrisemedical.co.uk. All combinations listed have been validated to meet the General Safety and Performance Requirements, section 14.1 of the Medical Device Regulation 2017/745.

Guidance on the combination, such as mounting, can be found at www.SunriseMedical.co.uk.

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Sunrise Medical B.V.
Sunrise Medical HCM B.V.
Sunrise Medical S.A.S