



# Technical / News Bulletin

UK

BULLETIN NO.:

TB 0456

DATE ISSUED

04/02/2015

## BULLETIN DETAILS


<b>Product Group:</b>	Manual Wheelchairs		
<b>Model Affected:</b>	All manual wheelchairs with the option "desk sideguard"		
<b>Effective Date:</b>	Immediate	<b>Serial No. Effective From:</b>	N/A
<b>Bulletin Subject:</b>	Play reduction in desk sideguard		

### Description of issue

Previously, the armpad of the desk sideguard had too much play.

### Solution – Refer to image below

The nylon screw (2) has been used to minimise the play. It is also adjustable to remove play in the field

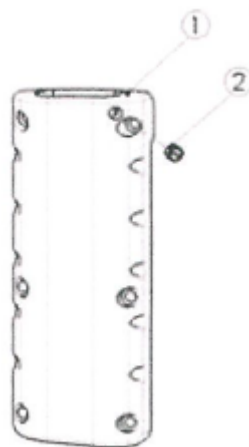

### Benefit to customer

1. Improved quality
2. No play in armrest
3. User feels more secure.

### Literature Affected:

<b>Sunparts:</b>	Running change. No updates required
<b>Order Forms:</b>	No changes required
<b>Website:</b>	No changes required
<b>Owners Manual:</b>	No changes required

If you have any questions regarding this bulletin, please do not hesitate to contact:

Customer Services	Customer Support	Technical Service Centre	E-mail
<b>Telephone:</b>	<b>+44 (0)845 605 66 88</b> Option 2	<b>+44 (0)845 605 66 88</b> Option 3	enquiries@sunmed.co.uk
<b>Fax:</b>	<b>+44 (0)845 605 66 89</b>	<b>+44 (0)845 605 66 89</b>	help.technical@sunmed.co.uk