

BULLETIN DETAILS

Product Group: Seating

Model Affected: J3 Backrest Range


Effective Date: Immediate effect

Serial No. Effective From: N/A

Bulletin Subject: Improvement to J3 Backrest Range

Subject of Bulletin:

In an ongoing effort to improve the design, functionality and performance of our JAY J3 backrest range, a further series of improvements have been introduced on the mounting hardware, comfort curve on Low Thoracic height shells and accessory bag. It also includes accessory discontinuations from some J3 options.

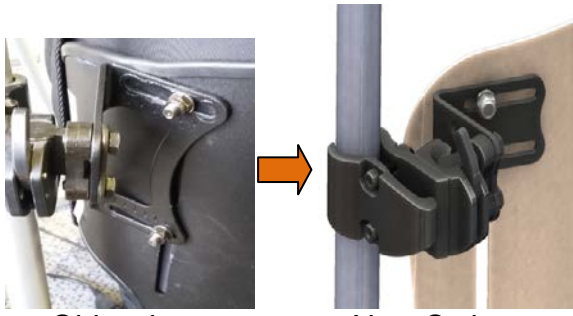
Feature	Benefit
<p>New Compact Mounting hardware</p>  <p>Old style New Style</p>	<p>Re-designed to meet the needs of active users:</p> <ul style="list-style-type: none"> ✓ 25% lighter than the current 2-point JAY J3 quick release hardware ✓ The lightest backrest mounting hardware available in the market ✓ Easier to fit and adjust ✓ Adjustable in width (5cm), depth (3cm), angle (+/- 22°), and height (by movement of brackets attached to shell or movement of backposts clamps) ✓ Easily fits majority of backpost tubes ✓ Longer clamp bolts – prevents clamp separation when fitting to wider back post tubes (20 and 30mm long bolts are supplied as standard) ✓ 136kg max user weight ✓ Transit tested and approved <p>The new compact mounting brackets will be included with every new Low thoracic and Mid thoracic J3 backrest but also available as a spare part.</p> <p>The old style hardware (code J3EUSH) will still be available as a spare part for the next 5 years.</p> <p>Please refer to http://parts.sunrisemedical.eu</p>
<p>New Heavy Duty Compact Mounting Hardware</p>	<p>The old 4-point mounting hardware and extended depth hardware will be discontinued and available only as spare parts</p> <p>Backrest depth extension brackets - code J3EUAHB 4 point hardware - code J3EU4PTHDSH</p> <p>They will be replaced by new Heavy Duty Compact hardware.</p>



Old style



New Style

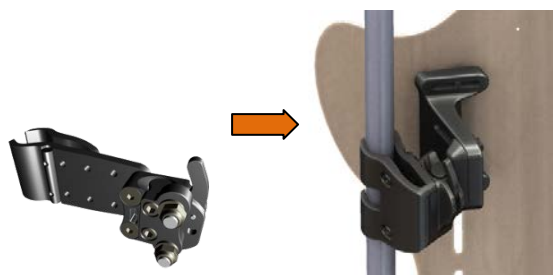

Benefits include:

- ✓ Designed for more passive or heavier users
- ✓ Supplied as standard on all **Upper Thoracic and Shoulder High** backrests
- ✓ **20% lighter** than old 4-point style but just as stable
- ✓ **Significantly less play and movement** due to re-design of clamp, cone and receiver
- ✓ Easier to fit and adjust
- ✓ No loss of angle, depth, width and height adjustment
- ✓ Improved styling
- ✓ **Longer clamp bolts** – prevents clamp separation when fitting to wider back post tubes (25mm and 30mm long bolts are both supplied as standard with all J3 backs)
- ✓ Higher maximum user weight of **225kg**
- ✓ Transit tested and approved

New comfort bend on top edge




The top edge of all low thoracic shells ordered from this point onwards will feature a comfort bend.

Benefit:

- The rearward flared edge reduces the pressure and discomfort at the top edge when leaning back to dress, function or alter position.

New Shell mounting plate on JAY J3 Carbon backrest

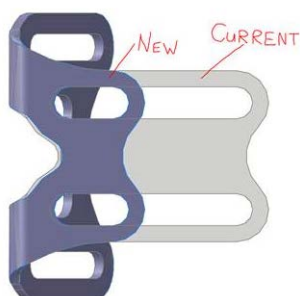

The shell mounting plate on the JAY J3 Carbon backrests will be replaced by a smaller, lighter bracket.

Benefits:

- ✓ 5% lighter compared to the current J3 Carbon quick release and fixed mounting brackets
- ✓ Reduced depth adjustment from 4.5cm to 2cm for a more compact, stylish look
- ✓ Retains the key range for adjustment of width, depth, angle and rotation

New Accessory Bag

The old style JAY J3 grey 'Vanity Flap' will be discontinued and replaced with the 'Accessory bag' from the JAY J3 Carbon backrest range.

Benefits:

- ✓ Central slot enables easy fitting
- ✓ Zipped compartment for storing essential items



Old Style

New Style

- ✓ Heavy weaved nylon for greater durability
- ✓ Shaped to ensure it is discreet and stylish

Please refer to www.sunrisemedical.co.uk, or contact your area sales manager for more information

If you require any further technical assistance please contact:

Matthew Eveleigh, JAY Clinical Specialist on matthew.eveleigh@sunmed.co.uk

Literature Affected:

Sunparts: Available from <http://parts.sunrisemedical.eu>

Order Forms: Updated and available via Customer Services

Website: Available on www.sunrisemedical.co.uk

Owners Manual: Available on www.sunrisemedical.co.uk

If you have any questions regarding this bulletin, please do not hesitate to contact:

Customer Services	Customer Support	Technical Service Centre	E-mail
Telephone:	+44 (0)845 605 66 88 Option 2	+44 (0)845 605 66 88 Option 3	enquiries@sunmed.co.uk
Fax:	+44 (0)845 605 66 89	+44 (0)845 605 66 89	help.technical@sunmed.co.uk