

Q700[®]

The Ultimate Experience



Q Series
Powerchairs

Created individually for whatever your clinical and lifestyle needs are.



Q700 M




Q700-UP M

THE ULTIMATE EXPERIENCE.



Q700 M

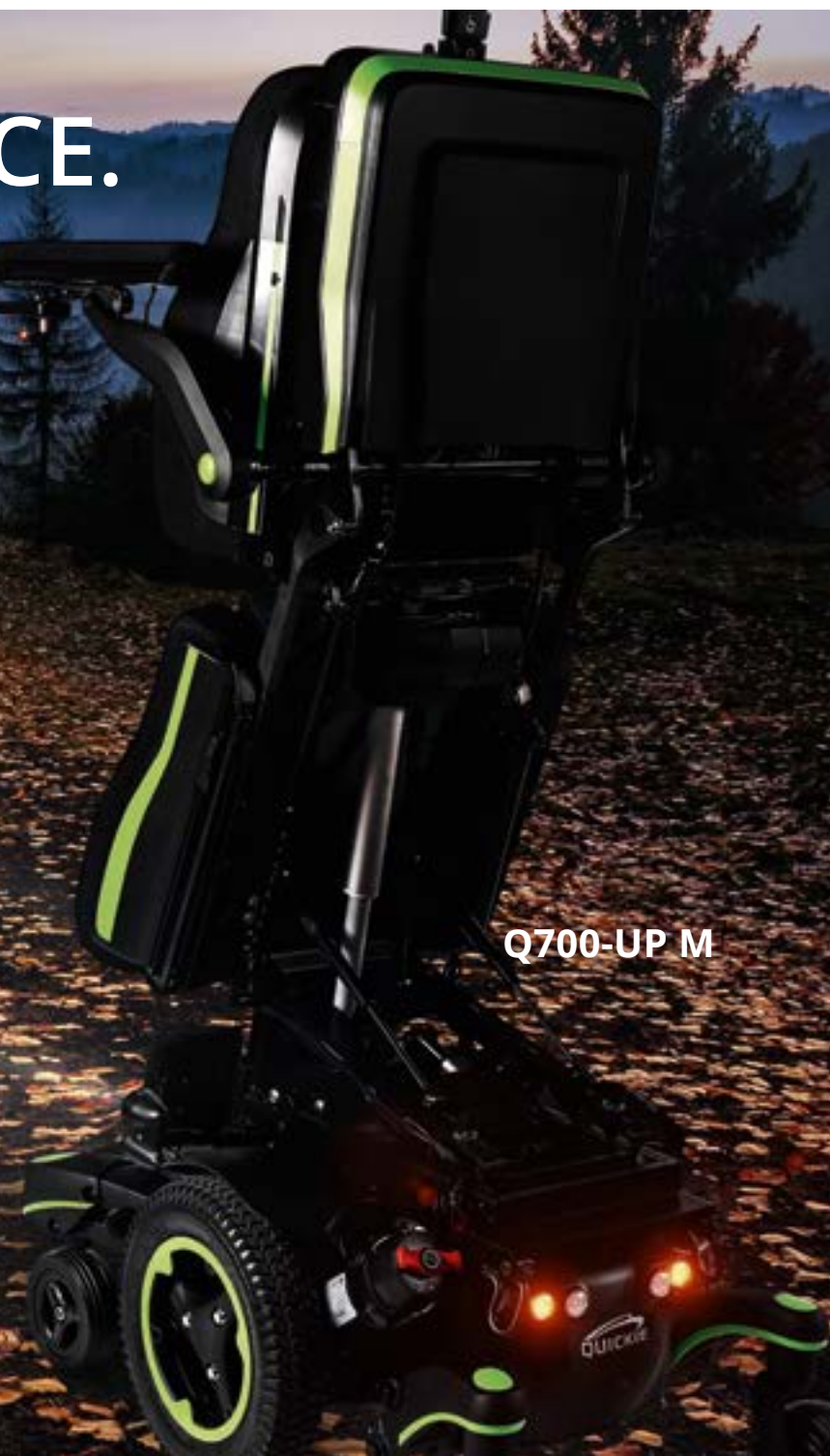

Q700-UP M

The QUICKIE brand delivers the world's most comprehensive range of innovative, technically advanced premium power chairs, created individually for whatever your clinical and lifestyle needs.

Resulting from years of experience with market leading drive base, seating and advanced electronics, the QUICKIE® Q700® M and Q700®-UP M offer the highest performance and technology in providing you with the ultimate experience!

HIGH-PERFORMANCE BASE

Advanced mid-wheel drive technology with the power and driving performance to get you where you want to go. Superb indoor agility combined with SpiderTrac™ 2.0 suspension and patented anti-pitch system for excellent outdoor performance and unbelievably smooth kerb climbing transitions.

INTELLIGENT SEATING

Mirroring your body's natural positions, the SEDEO ERGO's biomechanical seating ensures you'll never lose contact with your footplate, headrest or chin controls during seat transitions - maintaining good posture when tilted, reclined, elevated and in every position in between. Access up to 6 programmable seat positions all at the touch of a button.

INNOVATIVE CONTROLS

Seamlessly connect with your environment, your phones, computers and other electronics with a choice of industry-leading, intuitive and easy to operate controls, including Infrared and Bluetooth.

DISTINGUISHED STYLE

A streamlined, minimalistic design contrasted with sleek, black bodywork and striking colour inlays. Everything from the upholstery down to the wheels have been picked to turn heads.






MID-WHEEL DRIVE PERFORMANCE, EVOLVED.

The Q700 mid-wheel drive base sets the bar for powered wheelchair technology, inside and out. Travel further with its superior drive performance and feel secure over challenging terrain with SpiderTrac 2.0 suspension and its independent movement of all 6 wheels. It's a driving experience like never before.

TRAVEL WITH CONFIDENCE

With ultra-reliable high torque 4-pole motors powering speeds of 6, 8, 10 or 13 kph and a choice of 60 or 80 Ah batteries to maximise your range, you can always travel in confidence with the Q700 M.

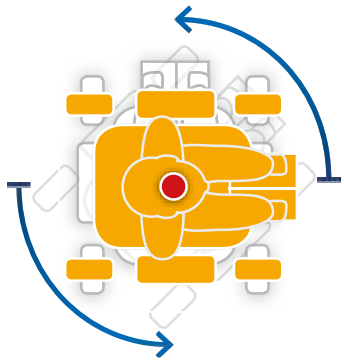

DESIGNED FOR INDOOR MANOEUVRABILITY

Turning on its own axis, enjoy intuitive driving and excellent manoeuvrability in tight compact spaces with this narrow (620 mm) base and a 1260 mm turning circle. Also, easily slide under tables and access Wheelchair Adaptive Vehicles (WAVs) with a low seat height from just 445 mm.


BUILT FOR OUTDOOR PERFORMANCE

Climb kerbs of up to 100 mm and enjoy a smoother, more responsive ride with the revolutionary SpiderTrac 2.0 suspension. With its independent movement of all six wheels, the front and rear wheels give the Q700 M the ability to crawl up and down steep obstacles whilst keeping all 6 wheels to the ground – even when ascending gradients of up to 10 degrees.




KEEPING YOU ON TRACK WITH GYRO CONTROL

The optional Gyro tracking system automatically detects and corrects the smallest deviation from your direction. By providing a pin-point accurate drive, it'll reduce the number of physical or cognitive actions required to keep you on-track – making for a less-demanding drive (especially if using special controls).

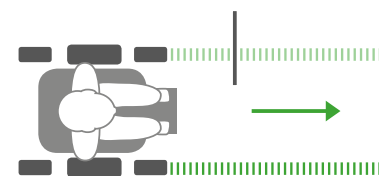

WITH GYRO CONTROL MODULE

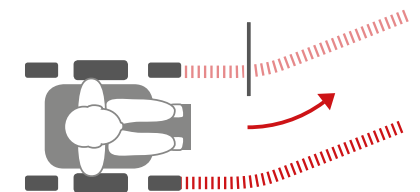

WITHOUT GYRO CONTROL MODULE



SEDEO ERGO

INTELLIGENT SEATING

It's the first-class seat that's designed specifically to fit you. The SEDEO ERGO with patented Biometric Repositioning is more effective than any other seating system when it comes to independent repositioning throughout the day.

EASY-ENTRY SYSTEM

The seat cushion is divided into 3 elements.

The two lateral side cushions can be swung-out and lifted off, ensuring a comfortable transfer. The middle part can also be substituted for a JAY cushion.

ADJUSTABLE TO FIT YOU

The premium headrest (1) can easily be adjusted or removed. All seat widths and depths (2) are infinitely adjustable with just one tool. Also, the seat and backrest width (3) can be adjusted independently of each other. Even the armrest height and angle-adjustment (4) are simple to adjust too.

PRESSURE RELIEF REMINDERS

Pressure relief reminders tell you when it's time to change position. A log can also be kept for clinicians to review.

PROGRAMMABLE SEATING

The SEDEO ERGO has up to 6 uniquely customisable seating positions that can be programmed to any position you desire. Whether you want to lie-down, cook or watch TV – every position can be completely customised to you.

FOOTPLATE TO THE GROUND

As the footplate goes all the way to the ground, you're able to do standing transfers in-and-out of your Q700-UP M.



BIOMETRIC REPOSITIONING

Moving out of a reclined position often requires use of the tilt function or to regain your proper seated posture.

With the patented Biometric Repositioning, when you move from reclined to seated, the SEDEO ERGO will first pass through an intermediary tilted position, using gravity to realign your pelvis in one smooth, effortless motion.

ANTI-SHEAR BACKREST

Transitioning in and out of recline doesn't need to be uncomfortable. The SEDEO ERGO's power recline is supported by 180 mm of anti-shear to better maintain the positioning of the headrest, and even chin control, during recline.



Q700 M

SEDEO ERGO SEATING SYSTEM

Live life on-the-go with the Q700 M SEDEO ERGO. Combining all of the benefits of the SEDEO ERGO seating system with the ability to travel at your pace when elevated with C-Me.



LIFT UP AND DRIVE WITH C-ME

No more waiting around, or driving at reduced speed when lifted – now you can keep up with friends as you rise up to stand out from the crowd with C-Me. Enabling you to rise up to 300 mm whilst travelling at speeds of up to 5 kph, you can also lift or descend at the touch of a button with no need to stop whilst transitioning. Just rise up and go at the same time!



TILT



RELAX



LIFTED (C-ME)



LIE DOWN

MEMORY POSITIONS

Get the most out of the Q700 M's powered option with Memory Seating. Save up to 4*programmable positions and transition to them with a simple press of a button.

*4 programmable positions by default, 2 more as an option.



Q700-UP M

SEDEO ERGO UP SEATING SYSTEM

For those who want to experience life from new heights, discover SEDEO ERGO UP. Combining all of the benefits of the SEDEO ERGO seating system with the power to stand whenever you want.

THE BENEFITS OF STANDING

With a standing angle of up to 85°, the SEDEO ERGO UP allows to stand whenever you want and offers many health and lifestyle benefits:

- Improved pressure relief for those who may have limited ability to independently weight shift or have to sit for long periods
- Enhanced function and ability to reach and access your environment
- Improved range of motion and a reduction of the risk of contractures
- Supporting vital organ capacity, including pulmonary, bowel and bladder function and improves circulation
- Promoting bone health and reducing the occurrence of skeletal deformities by encouraging standing
- Reducing abnormal muscle tone and spasticity

90° KNEE ANGLE


The new centre-mount footrest gives a true 90° knee angle in seated position. The ergonomic knee pad is also depth adjustable.




LIE DOWN




RELAX

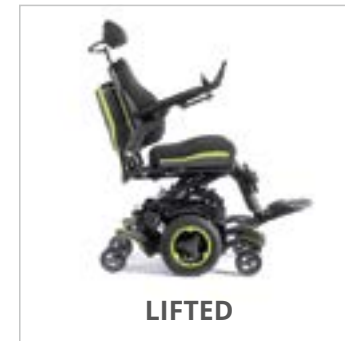


LIFTED

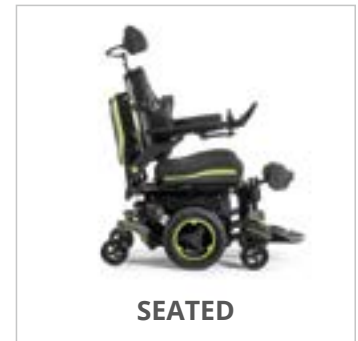

SEATED

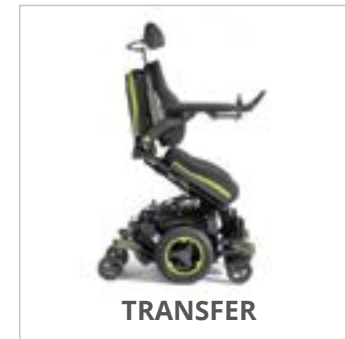

TRANSFER

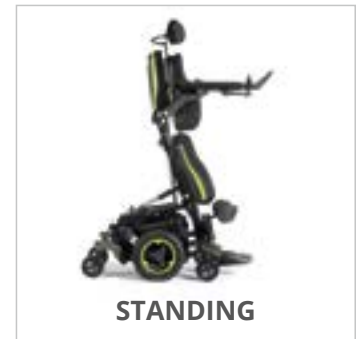

STANDING

MEMORY POSITIONS

Get the most out of the Q700-UP M's powered option with Memory Seating. Save up to 6 programmable positions and transition to them with a simple press of a button.

JAY® & WHITMYER® POSITIONING OPTIONS

A wide range of WHITMYER headrests and JAY cushions and backs are available with the SEDEO ERGO and SEDEO ERGO UP seating systems to meet every user's clinical needs and ensure we offer the best seating possible.



JAY CUSHIONS & BACKS

Designed for pelvic stability, skin protection and comfort, JAY cushions and backs are low maintenance, highly customisable and are available in a large variety to meet your needs.



WHITMYER HEADRESTS

Created to address a range of user needs from basic support to multiple-point control. Whitmyer headrests come in a variety of pad shapes and sizes for superior head positioning and comfort.





SWITCH-IT REMOTE SEATING APP

The Switch-It Remote Seating App* measures the seating angles and tracks the duration the client remains in the position in order to improve his/her overall health through pressure management.

- Bluetooth technology seamlessly connects to your wheelchair
- Live reports are created on the angle and duration of tilt, recline, leg elevation
- Reports can be shared with therapists and physicians via email
- Helps create pressure relief plans
- Alerts user to adjust seating



* Included as standard with the optional Ctrl+ 5 box.

A person wearing a blue and black jacket is shown from the side, holding a handheld electronic device. The device has a small screen displaying a graph and several buttons. The background is dark.

CONTROL AT YOUR FINGERTIPS

Be the master of your environment with the Q700 M Series and its advanced controller technology. Not only can you control your Q700 M the way you want, but also your household appliances, phones and computers.

R-NET

ADVANCED JOYSTICK

Every Q700 M Series comes with the superior R-Net advanced joystick as standard. A large, easy-to-read HD display and simple to control paddle switches allow for intuitive operation.

CONNECT

Seamlessly take control of your environment using the integrated Infrared functionality. Control your home entertainment system, TVs, DVDs, docking stations and other home entertainment systems simply through your joystick.

BLUETOOTH & IDRIVE

Control your phone, computer, tablet and any other Bluetooth and i-devices with this useful optional extra. Connect and control up to four devices.

QUICKIE'S PATENTED ASSIGNABLE BUTTONS

Change the function on any button on your joystick to any command you desire, thanks to QUICKIE's patented assignable buttons. You can also assign a second (long press) function to buttons, giving you quicker, simpler and smoother access to functions you commonly use without having to use your joystick.



SWITCH-IT

INNOVATIVE CONTROLS

Conveniently access your 5 programmable positions with the Ctrl+5 – a perfect combination with Memory Seating. It allows you to go straight to the function you need at the touch of a button.

ALTERNATIVE DRIVE CONTROLS

If driving with a traditional joystick is not the best option, alternative drive controls allow you to fully experience your chair's driving capabilities and benefits.



LINK-IT MOUNTING SYSTEM

Need to use additional switches or special controls? No problem, easily connect them using the universal Link-It Mounting System.

By using a combination of available modules and rods, you can create an endless number of combinations to accommodate what you need. Not only is it easy to position, Link-It also ensures the minimalistic, clean design of your wheelchair isn't compromised by additional controls.

So whether it's an extra controller on the armrest, a push button on the backrest or a switch on the headrest, just link it.





A close-up, low-angle shot of a black Quickie motor vehicle. The vehicle features a prominent white 'QUICKIE' logo on its side panel, which is accented with a red swoosh. Above the logo, two circular headlights are visible, each containing a red and a clear lens. The vehicle's body is black with red decorative stripes. A black tire with a multi-spoke wheel is in the foreground, and another tire is partially visible on the right. The background is dark and out of focus.

QUICKIE



DISTINGUISHED **STYLE**

With the Q700 M series, we've taken style to the next level. With a matt-black chassis finish, hidden wiring and components combined with sharp-tailored upholstery, it's a sleek and contemporary design.

The Q700 M also comes with a choice of 6 contemporary colour inlays to perfectly complement its striking matt-black bodywork. Fancy a change? Simply replace the colour inlays. Styling your powered wheelchair has never been so simple.



TECHNICAL DATA

	Q700 M	Q700-UP M		Q700 M & Q700-UP M
Max. user weight:	160 kg	120 kg	Max. reversing width:	1000 mm
Seat width:	400 – 560 mm	400 – 500 mm	Max. safe slope:	10°
Seat depth:	400 – 560 mm	400 – 540 mm	Max. kerb climbing:	100 mm
Seat height:	430 – 480 mm	450 – 480 mm	Wheel options:	14" drive wheels; 6" castor wheels
Backrest height:	510 mm – 610 mm	560 mm – 610 mm	Colour inlay choices	6 colours: red, blue, green, orange, white, black
Backrest angle:	85° to 175° (powered)	85° to 175° (powered)	Max. range (ISO 7176-4)	40 km (with 80 Ah batteries & 6 kph motors; dependant)
Total width:	620 mm	620 mm	Wheelchair weight:	from 160 kg (including 60 Ah batteries)
Total length:	1070 mm	1070 mm	Crash tested:	Yes - approved for forward facing transport in a vehicle (ISO 7176-19). Attachment: 4 chair fixing points with 2 double straps rear & 2 single straps front. No anchored belt restraints allowed. The ease of access and manoeuvrability in vehicles depends on the wheelchair size.
Total height:	1100 mm	1100 mm		
Seat memory functions:	Lie-down, seated, relax, tilted, lifted, free programmable	Standing, seated, transfer, lie-down, lifted, relax, tilted, free programmable	Transport without user:	Fix the wheelchair safe in the vehicle (e.g. car or plane)
Max. speed:	6 & 8 kph (standard), 10 & 13 kph (optional)	6 & 8 kph (standard), 10 & 13 kph (optional)	Intended use and environment:	For people with limited mobility; self-propelling or being pushed by an assistant, indoor and outdoor use. Wheelchair type class B
Battery size:	60 Ah; 80 Ah (optional)	60 Ah; 80 Ah (optional)		
Electronics:	R-Net controller; lights & indicators	R-Net controller; lights & indicators		



Sunrise Medical GmbH

Kahlbachring 2-4

69254 Malsch/HD

Tel.: +49 (0) 72 53 / 9 80-0

Fax: +49 (0) 72 53 / 9 80-222

E-Mail: international@sunrisemedical.de

www.SunriseMedical.eu

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.eu.

