

BULLETIN DETAILS

Product Group:	MANUAL WHEELCHAIRS		
Models Affected:	BREEZY RELAX2		
Effective Date:	JULY 2018	Serial No. Effective From:	
Bulletin Subject:	IMPROVED TILT/RECLINE RELEASE AND ELR WITH PROTECTIVE COVER		

Description:

In an ongoing effort to improve the design, functionality and performance of our Breezy Relax2, we have undergone an important update in the manufacturing process of the tilt and recline lever and also a further improved design of the cover for ELRs.

We have had feedback of the occasional breakages of the tilt and recline lever due to heavy use and wanted to improve the design. We have made the following improvements to the material and design:

- Lever made of Bayer PA6 without glass
- Lever reinforced (increased material thickness)
- Change the colour from black to grey (for easy differentiation)



We have also separated the lever from the drum brake. This is to ensure that breakages are reduced due to incorrect operation; We have found there have been cases of application of too much force, due to mistaking the correct lever.



Our new ELR cover has been designed for increased safety, the entrapment improvement was generated via a risk assessment only. The cover is clicked on and can also be retrofitted to existing ELRs that are currently in the field.




Old



New



Old

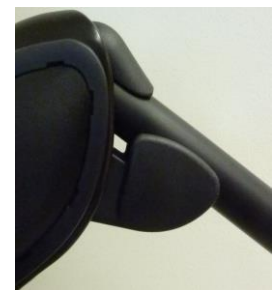

New

Important Note: This update is imminent and all chairs will now be supplied following this manufacturing process. Thank you for your product feedback and allowing us to continue to improve our products to better suit your user's needs.

Benefits now given to the end user:

- Improved reliability of the tilt and recline lever (we do not expect to receive any further breakages of this part).
- Increased safety for ELRs and ensuring that fingers are not accidentally caught.

Please refer to www.sunrisemedical.co.uk, or contact your Area Sales Manager for more information

Literature Affected:

Website:	No
Order forms:	No
Owner's Manual:	No
Parts Manual:	No

If you have any questions regarding this bulletin, please do not hesitate to contact:

Customer Services	Customer Support	Technical Service Centre	E-mail
Telephone:	+44 (0)845 605 66 88 Option 2	+44 (0)845 605 66 88 Option 3	enquiries@sunmed.co.uk
Fax:	+44 (0)845 605 66 89	+44 (0)845 605 66 89	help.technical@sunmed.co.uk