



UNFOLD YOUR FREEDOM

QS5x



Attempting the action shown in the image may cause injury or equipment damage. We recommend the use of anti-tips. The manufacturer does not bear any responsibility for injury or damage resulting from disregard of the instructions described in the owners manual.



5 GOOD REASONS TO TRY THE QS5 X[®]

1 FREEFOLD - OVER 40% ENERGY SAVED WITH EVERY FOLD - The super-precise FreeFold design delivers folding forces that require more than 40% less force than its closest competitor.

2 LIFTING IS A BREEZE - At a remarkably low transport weight from 7,3 kg for the Fixed Front and 7,5 kg for the Swing Away frame, it is the lightest aluminium box frame wheelchair on the market setting a new standard.

3 RIDE LONGER AND PROPEL MORE EFFORTLESSLY - Boasting not only the lightest but also one of the stiffest frames on the market thanks to proprietary manufacturing processes and technologies signature to Sunrise Medical's high-quality products - such as specially heat treated 7000 series aluminum, ovalized tubes and forged parts - the QS5 X frame design takes the driving experience to the next level.

4 STYLISH BUT ROBUST - With a sturdy box frame construction and the ovalized 4-arm cross brace, you can have peace of mind knowing your chair can handle whatever comes your way. Its sleek, modern lines give it a stylish and sophisticated look.

5 ADAPTABLE TO INDIVIDUAL NEEDS - Two frame versions, various frame angles and a tailored package of clinical-critical components, ensure it can adapt to individual needs, no matter if these may change over time - from the first-time user to the experienced wheelchair user.



Fixed-Front



Swing-Away

UNRIVALED LIGHTWEIGHT AND SUPERIOR RIGIDITY

- The QS5 X offers easy propulsion and lifting thanks to its lightweight aluminum frame. Made from heat-treated 7000 series aluminum for strength without added weight, it provides effortless maneuverability and efficient transfer of push to the wheels.

The **OVALIZED TUBING** of the cross-brace makes the QS5 X less likely to flex under pressure providing the superior rigidity and reduced flex.

Bright **LED LIGHTS** provide a great feature to light your way.



ALWAYS STRAIGHT TRACKING thanks to 3-way adjustable, **patented caster-link** avoiding inefficient counter steering.



Choose from **24 VIBRANT FRAME COLORS** and 6 cast-er fork and caster rim colors to create the perfect look!



INNOVATIVE FREEFOLD DESIGN - With more than **40% less force required** than other folding wheelchairs in the market, the QS5 X sets the new standard for ease and efficiency of folding.



The QS5 X features **TOUCH POINTS** on the caster link, knee-bend, and back frame that provide a **comfortable grip** and prevent slipping.

ROBUST BOX FRAME ensures the QS5 X can handle anything.

Featuring a reduced design with hidden adjustments that gives the QS5 X a **STYLISH** and sophisticated appearance.



PERFECT COMPANION FOR POWER ASSIST DEVICES - The front frame features high-quality forged parts - designed to withstand the high load power assist devices apply to the frame and deliver premium performance.

FITS YOUR LIFESTYLE

The **3D VISUALIZER** lets you see how your QS5 X will look like and what it will weigh.


Broad range of options to build the QS5 X that fits YOUR individual needs.

Easily upgrade your QS5 X with various **JAY®** cushions and backrest.

Want something unique? Our **BUILT-4-ME** team can help you create a truly individual QS5 X that's out of this world! Ask your dealer for more information.

JAY

built4me

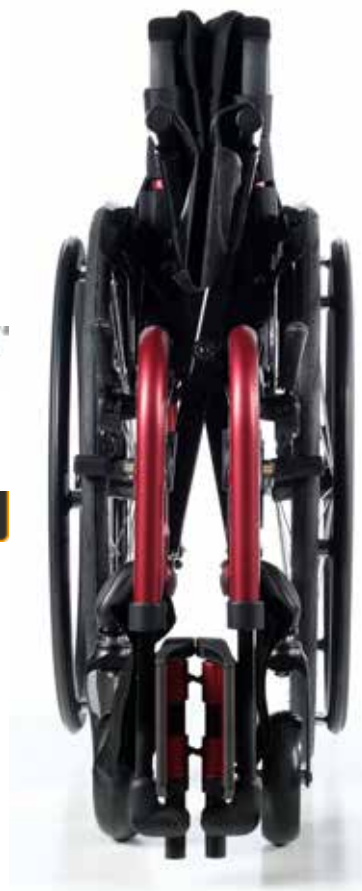

TECHNICAL DATA	QS5 X
Max. user weight	136 kg
Seat width	13 - 20" (in 20 mm steps)
Seat depth	13 - 20"
Front seat height	14 - 22"
Rear seat height	14 - 20"
Backrest height	10 - 20"
Backrest angle	from -15° to +15°
Total width	SW +180/240/270 mm at 0°/2°/4°camber and 25" wheels The total width can exceed 800 mm in certain configurations.
Total length	1113 mm
Total height	1020 mm
Max. pivot width	840 mm
Max. safe slope for use of parking brakes	7°
Frame versions	Swing Away or Fixed Front
Frame angle	Swing Away: 70° / 80° / 90°, Fixed Front: 75° / 85°
Front wheel	4", 5", 6", 7"
Rear wheel	22", 24", 25", 26"
Wheelchair weight	Swing Away: from 10,1 kg, Fixed Front: from 9,9 kg
Dismanteled weight	From 7kg
Mass of heaviest part	2,2 kg
Removable parts	Rear wheels (for transportation)
Frame colours	24 colours
Crash tested	Yes (approved for transport in a vehicle according ISO 7176-19)
Transport without user	Fix the wheelchair firmly and securely in the vehicle (e.g. car or plane)
Intended use and environment	For people unable to walk and with limited mobility; inside and outside. Operation by user or assistant.

Sunrise Medical Australia Pty Ltd

11 Daniel Street, Wetherill Park, NSW 2164
Phone: (02) 9678 6600
E-Mail: enquiries@sunrisemedical.com.au
Website: www.sunrisemedical.com.au

For further information on the full specification, options, and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.