

# MANUAL WHEELCHAIR REFERENCE GUIDE

---

AUSTRALIA





## YOU DRIVE OUR DESIGN

QUICKIE wheelchairs are designed for you and by you. In development, we consider how you interact with your wheelchair to maintain maximum efficiency, performance, fit, and style.

## LEADING TECHNOLOGY

We are continually developing new products, improving our flagship products, and offering new options and solutions to enhance your lifestyle.

## QUALITY AND DURABILITY

QUICKIE wheelchairs undergo rigorous testing that surpasses industry standards, so you can be confident that your product will hold up and perform throughout its lifetime.

## ADJUSTABILITY AND VERSATILITY

Once you choose a QUICKIE, our goal is that it continues to be the right chair for you. We offer maximum adjustability in our designs, interchangeable parts, and the True Fit Growth Program\* to help accommodate changes in your condition, environment or lifestyle.

\*Available on select QUICKIE wheelchairs

## BUILT-4-ME PERSONALISED

Our Built-4-Me Team is dedicated to creating product user-matched beyond what is listed on our product order forms. Experts in the engineering and manufacturing of configured rehab products will listen to your specific requests and build a chair as individual as your fingerprint!

## RIGID FRAME

Nitrum	4
Nitrum Hybrid	5
Krypton R	6
5 <sup>R</sup>	7
GP Series	8
QRi	9
RIGID Comparison Guide	10

## FOLDING FRAME

Krypton F	12
Xenon <sup>2</sup>	13
QS5 X	14
Quickie 2	15
QXi	16
M6	17
FOLDING Comparison Guide	18

## TILT IN SPACE

Iris	20
Iris with Power Tilt	21
SR45	22
TILT IN SPACE Comparison Guide	23

## JAY® Seating

JAY® Seating Selection	24
------------------------	----



# NITRUM


## LIGHTEST IN CLASS

ENGINEERED  
TO EFFICIENTLY  
CONVERT YOUR  
ENERGY INTO MOTION

*Nitrum boasts the largest range of sizes and adjustments in its class, so you can achieve your own perfect fit*



### LIGHT AND BACK SUPPORT

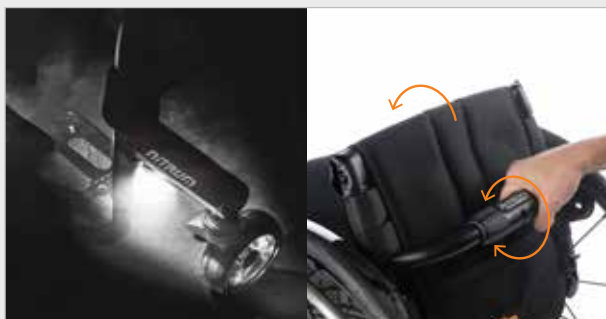

**LED Lights** mounted on each caster arm allow you to see and be seen.

**One-handed folding back support** offers you simple one-handed release and a clean look.

### PATENTED, ADJUSTABLE CASTER SYSTEM



Instead of several welded parts, we use a reinforced and patented caster mechanism in a single hollow- forged component.

#### Fine-Tuned Adjustment

- Angle
- Camber
- Height



# NITRUM

## HYBRID FRAME

MOST RIGID,  
STILL LIGHTWEIGHT

*Offers a reinforced front frame without compromising its proven lightweight construction.*



### FREESTYLE BACK OPTION



Freedom to move, ergonomic pushing.

### TENSION ADJUSTABLE SEAT SLING



Breathable fabric with block-and tackle tension-adjustment is easy to tighten, and its double-loop mechanisms will remain secure.

Helps position the pelvis lower than the thighs for additional postural stability.

# KRYPTON R

## CARBON FIBRE FRAME

REVOLUTIONARY  
CARBON FRAME  
WITH STAGGERING  
STRENGTH-TO-  
WEIGHT RATIO.

### RIGID FRAME

*Ultra-lightweight from  
3.65kgs in transit*



### CARBON BRAID TECHNOLOGY



This QUICKIE Krypton Range features a patented and fully automated BRAID process that makes it possible to build the perfect frame – strong, light and manufactured with a precise and reliable automated process.

### RIGIDITY & RAZOR-SHARP AGILITY



It's not only the exquisite material or the famously seductive styling, but also the athleticism and agility of a unique rigid wheelchair. Infinitely adjustable centre of gravity, fine-tuning of the front and rear seat heights and backrest angle adjustment ensures the Krypton R is designed for something truly special.



# QUICKIE 5<sup>R</sup>

## HIGHLY ADJUSTABLE FRAME

HIGHLY  
ADJUSTABLE  
TO ELEVATE YOUR  
MOBILITY

### Frame Styles:

#### CLASSIC

*Traditional frame length  
for maximum stability*

#### PERFORMANCE

*Shortened length for  
maximum responsiveness*



### AXLE PLATE DESIGN



The 5R's swept-back axle plate design isn't just for looks.

It is also incredibly quick to adjust and is available in two lengths to ensure you are getting the adjustments you need without adding extra weight.

### XTR SUSPENSION



The XTR Suspension Frame helps you fearlessly conquer curbs, bumps, and tough terrain. Add Mountain wheels, Knobby tires, wide casters, and Frog Legs suspension caster forks for the ultimate off-road wheelchair experience.

# GP SERIES

## RIGID WHEELCHAIR

### A CLASSIC QUICKIE RIGID WHEELCHAIR

#### **Available in:**

**GP** - Standard Style

**GPV** - Tapered for  
increased maneuverability

**GPSA** - Swing Away  
for easy transfers



#### THE ORIGINAL RIGID FRAME



The GP's classic box frame and one piece construction exceed our superior standard for quality and deliver an unmatched rigid feel.

#### SPORTS OPTION



The tapered GPV front frame, optional Quickie Performance wheels, and up to 12° of wheel camber combine for exceptional maneuverability on and off court.



## QRi RIGID MANUAL FRAME

NO-NONSENSE  
UTILITY AND  
UNBEATABLE VALUE



### COMMUTER FRIENDLY DESIGN





- Ultra lightweight and durable aluminum frame, WC-19 Transit approved up to 120kgs
- Open frame and fluid caster housing design for a stylish look
- Small footprint for maneuverability and easy transport

### COMFORTABLE FIT



The QRi comes standard with 0" or 1" of frame inset, a highly adjustable axle plate for perfecting your Center-of-Gravity setting, and multiple front frame angle options.



		
	Nitrum	Nitrum Hybrid
SEAT WIDTH	12 - 20"	12 - 20"
SEAT DEPTH	12 - 20"	12 - 20"
FRONT SEAT HEIGHT	16 - 21"	16 - 21"
REAR SEAT HEIGHT	13 - 20"	13 - 20"
FRAME STYLE	7020 Oval Aluminium	7020 Oval Aluminium
BACK HEIGHT	8 - 20"	8 - 20"
BACKREST ANGLE	-18° - 10°	-18° - 10°
FRAME / HANGER ANGLE	75° / 80° / 85°	75° / 80° / 85°
INSET	0" / 1" / 2"	0" / 1" / 2"
CAMBER	0° / 3° / 6°	0° / 3° / 6°
MAX. USER WEIGHT	125 kg	136 kg
PRODUCT WEIGHT	4.9 kg transit / 8.3 kg	5.8 kg transit / 9.2 kg
TRANSIT APPROVED	Yes*	Yes*
FRAME COLOUR	29 colours	29 colours
WARRANTY	Limited Lifetime	Limited Lifetime
ARTG	100167	100167

# RIGID FRAME SPECIFICATIONS

				
Krypton R	5 <sup>R</sup>	GP SA	GP/GPV	QRi
32 - 46 cm	12 - 20"	12 - 20"	12 - 22"	12 - 20"
34 - 48 cm	12 - 22"	12 - 22"	12 - 22"	12 - 20"
43 - 55 cm	16 - 21"	16.5 - 23"	15.5 - 23"	15 - 20"
38 - 50 cm	13 - 20"	16 - 20.5"	14 - 20.5"	14 - 20"
Carbon Fibre Open Frame	Open Series 5000 Aluminium	Traditional	Traditional	Open Series 5000 Aluminium
25 - 45 cm	10 - 20"	8 - 20"	8 - 20"	8 - 20"
-31° - 15°	-17° - 8°	-4° - 8°	-4° - 8°	-8° - 17°
92° / 100°	75° / 80° / 85°	60° / 70° / 90°	60° / 70° / 75° / 80°	75° / 80° / 85°
0" / 1"	0" / 1" / 2"	N/A	N/A on GP, 1" V taper on GPV	0, 1" V inset
0° / 3°	0°, 3°, 6°	(Adjustable) Set at 3°	(Adjustable) Set at 3°	0° / 3°
125 kg	120 kg	113 kg	113 kg	120 kg
3.65 kg transit / 6.2kg	6.4 kg transit / 9.8 kg	8.6 kg transit / 12 kg	7.6 kg transit / 11 kg GP	6.6 kg transit / 10 kg
No	Yes*	ISO Only - thru frame	N/A on GP, ISO Only - thru frame on GPV	Yes*
5 Finishes	28 colours	28 colours	28 colours	24 colours
Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime
100934	100167	100167	100167	100167

*\*Conforms to ANSI/RESNA WC-19 Standards. Please refer to the Quickie & Zippie WC-19 Transit Securement Supplement*



# KRYPTON F

## CARBON FIBRE FRAME

REVOLUTIONARY  
CARBON FRAME  
WITH STAGGERING  
STRENGTH-TO-  
WEIGHT RATIO

### FOLDING FRAME

*Ingeniously designed cross-brace with an open frame feel!*



### CARBON BRAID TECHNOLOGY



This QUICKIE Krypton Range features a patented and fully automated BRAID process that makes it possible to build the perfect frame – strong, light and manufactured with a precise and reliable automated process.

### LIGHTEST FOLDING WHEELCHAIR



From just 5.8kgs in transit, the QUICKIE Krypton F is our lightest folding wheelchair yet! From every angle, the Krypton F is conceived to portray beauty. The heart of the Krypton F is its ingeniously designed cross-brace that fits so neatly under the seat it's barely noticeable giving a very minimalistic, open-frame feel.

# XENON<sup>2</sup>

## FOLDING FRAME

ONE OF THE  
LIGHTEST FOLDING  
WHEELCHAIRS IN  
THE WORLD!

### **Available in:**


**FF** - Open fixed front frame

**DUAL** - Reinforced  
fixed front frame

**SA** - Swing Away for easy  
transfers



### EASY TO FOLD MECHANISM

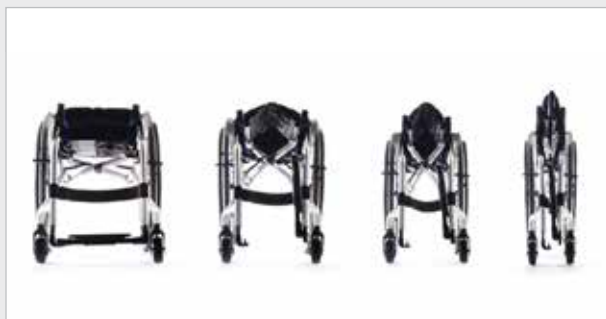

The aluminium latches on the cross-brace will also keep your wheels in parallel when folding and makes sure the chair stays folded when lifting.

### ADJUSTABLE TO SUIT YOU



You want to change your chair driving characteristic from passive to more active or the opposite way? Just change the centre of gravity setting by simply repositioning the axle stem. And don't forget to adjust the back angle for your ergonomic posture.

# QS5 X

## FOLDING FRAME

EASY PROPELLING  
AND LIFTING  
WITH THE  
FEATHERWEIGHT  
QS5 X


**Available in:**

**FF** - Open fixed front frame

**SA** - Swing Away for easy transfers



### UNRIVALED LIGHTWEIGHT

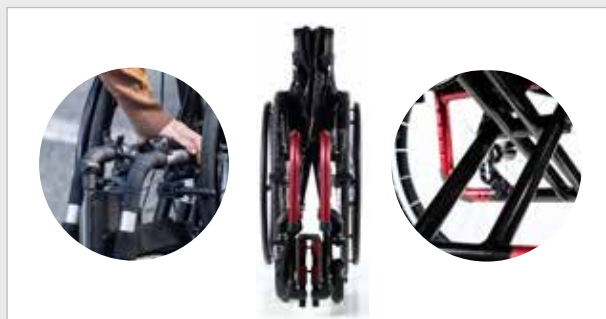

The QS5 X offers easy propulsion and lifting thanks to its lightweight aluminum frame. Made from heat-treated 7000 series aluminum for strength without added weight, it provides effortless maneuverability and efficient transfer of push to the wheels.

### SLEEK AND MODERN LINES



The QS5 X features a minimal design with hidden adjustments and seamlessly integrated caster link tubes. Few visible screws or holes to compromise aesthetics - this gives your QS5 X a stylish and sophisticated appearance.

# QUICKIE 2

## FOLDING WHEELCHAIR

VERSATILE AND  
MODULAR.  
THE INDUSTRY'S  
FAVORITE FOR OVER  
25 YEARS!

**Available in:**

**Q2** - Standard Modular  
frame

**Q2 HP** - High Performance



### MODULAR FRAME DESIGN



Featuring a modular design, the QUICKIE 2 offers the broadest range of possible configurations in the industry: over 100,000 different complex rehab options and counting. Designed to meet all of your unique needs without compromising fit or performance, only one wheelchair can live up to the name QUICKIE.

# QXI

## FOLDING WHEELCHAIR

LOW  
MAINTENANCE  
& SIMPLICITY



### DURABLE FRAME DESIGN



- Ultra lightweight and durable folding frame with a 120kgs weight capacity
- One-piece side frames reduce moving parts
- Three arm cross brace sits inside frame for rigid ride

### SIMPLE ADJUSTMENT & SERVICING



Most adjustments can be made in less than one minute, using only one tool - Center-of-Gravity (axle plate), Seat-to-floor height, Camber, Axle spacing and Back height.

# M6

## FOLDING BARIATRIC

IDEAL FOR  
BARIATRIC CLIENTS  
WEIGHING UP TO  
295 KG & UPTO  
75CM WIDE



### MODULAR & HIGH-STRENGTH FRAME



The M6's modular frame, features high strength steel for increased durability. Its interchangeable, heavy duty components accommodate changing conditions with simple parts replacements.

### LOCKING CROSS BRACE & STABILISER BAR



The M6 utilises a double locking cross brace and stabiliser bar to minimize frame flex and increase durability. The frame will hold an active user upto 295 kg.



**Krypton F**



**Xenon² SA®**



**Xenon² FF®**

	Krypton F	Xenon² SA®	Xenon² FF®
SEAT WIDTH	32 - 46 cm	32 - 50 cm	32 - 46 cm
SEAT DEPTH	34 - 48 cm	34 - 50 cm	34 - 50 cm
FRONT SEAT HEIGHT	43- 55 cm	43 - 56 cm	43 - 55 cm
REAR SEAT HEIGHT	37 - 50 cm	37 - 50 cm	37 - 50 cm
FRAME STYLE	Carbon Fibre Open Frame	Reinforced with swing-away front	Open frame with fixed front
BACK HEIGHT	25 - 47.5 cm	25 - 47.5 cm	25 - 47.5 cm
BACKREST ANGLE	75° to 103°	75° to 103° in 4° increments	75° to 103° in 4° increments
FRAME / HANGER ANGLE	92° / 100°	70° / 80°	92° / 100°
INSET	0 - 1 cm	N/A	0 - 1 cm
CAMBER	0° / 2° / 4°	0° / 2° / 4°	0° / 2° / 4°
MAX. USER WEIGHT	110 kg	125 kg	110 kg
PRODUCT WEIGHT	5.8 kg transit / 8.3 kg	7 kg transit / 10.3 kg	6.4 kg transit / 8.8 kg
TRANSIT APPROVED	No	Yes	Yes
FRAME COLOUR	5 Finishes	32 colours	32 colours
WARRANTY	Limited Lifetime	Limited Lifetime	Limited Lifetime
ARTG	100934	100934	100934

# FOLDING FRAME SPECIFICATIONS

				
Xenon² Dual®	QS5 X	Quickie 2	QXi	M6
32 - 50 cm	13 - 20"	12 - 24"	14 - 20"	22 - 30"
34 - 50 cm	13 - 20"	12 - 22"	14 - 20"	18 - 22"
43 - 55 cm	14 - 22"	13.5 to 22"	12.5 - 19.5"	17 - 20"
37 - 50 cm	14 - 20"	13 to 22"	12 - 19"	17 - 20"
Reinforced with fixed front	Traditional - Fixed or SA front	Traditional- Fixed or SA front	Traditional	Heavy Duty
25 - 47.5 cm	10 - 20"	9 - 21"	13 - 20"	17 - 20"
75° to 103° in 4° increments	-15° to +15° (optional)	-15° - +3° (Optional)	-15° - 3°	N/A
92° / 100°	SA 70°/80°/90° - Fixed 75°/85°	60°/ 70°/ 80°/ 90°	60°/ 70° / 80°	N/A
0 - 1 cm	0 / 1" each side	N/A	N/A	N/A
0° / 2° / 4°	0° / 2°	(Adjustable) Set at 3°	(Adjustable) Set at 3°	0°
140 kg	136 kg	136 kg Std / 158 kg HD	136 kg	295 kg
7.2 kg transit / 9.5 kg	7 kg transit / 9 kg	8.2 kg transit / 11.5 kg	9.9 kg transit / 13.25 kg	24 kg
Yes	Yes	Yes	Yes	No
32 colours	24 colours	24 colours	24 colours	24 colours
Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime
100934	100167	100167	100167	100167

# IRIS

## TILT IN SPACE

### INTELLIGENT ROTATION IN SPACE TECHNOLOGY

*One of the best and lightest  
tilt-in-space wheelchairs on  
the market*



### INTELLIGENT ROTATION



The Quickie IRIS features up to 55° of Intelligent Rotation in Space Technology to achieve superior positioning and is designed to facilitate feeding and respiratory function, reduce pressures and improve visual alignment.

### MONO BACKREST SYSTEM



The MONO Backrest System's single back post offers limitless seating and positioning capabilities. By mounting the backrest to a single central post, the backrest and chair widths may be selected independently to accommodate users with varying torso and hip widths.



# IRIS

## WITH POWER TILT

HAVE MORE  
CONTROL OVER THE  
POSITION OF YOUR  
BODY IN SPACE

### **Available in:**

#### **IRIS**

#### **IRIS Heavy Duty**

- improves comfort by assisting with sitting tolerance
- provides more independence to participate in daily activities



### ATTENDANT CONTROL



Allows the carer to give control of the tilt function to the user. The battery gauge is displayed on the top of the control so you never have to guess when it's time to charge. Ability to conveniently hang the attendant control anywhere on the chair.

### TOGGLE SWITCH



Can be mounted vertically or horizontally depending on needs. Charging port is conveniently located on the toggle switch for easy access.

# SR45

## TILT IN SPACE

PRE-CONFIGURED  
FOR EASE OF  
ORDERING



### 45° OF ROTATION-IN-SPACE



The Quickie SR45's seat frame rotates around the client's centre of gravity providing the shortest possible wheelbase without compromising stability. When the centre of gravity and centre of rotation are properly aligned, tilting is smooth, easy and controlled.

### STANDARD SMART PACKAGE



**SR45 comes standard with** - Non folding back canes, Height adjustable back with adjustable stroller handles, Foot tilt, Transit tie downs, Dual post height adjustable flip back with full length arm pads, 70" Swing in/ out hangers, solid seat pan and rear anti tips and composite foot plates.

# TILT IN SPACE SPECIFICATIONS



**IRIS**

**SR45**

SEAT WIDTH	14 - 22"	16 x 16"
SEAT DEPTH	14 - 22"	18 x 18"
FRONT SEAT HEIGHT	12.5 - 19"	20 x 18"
REAR SEAT HEIGHT	12.5 - 19"	17"
FRAME STYLE	Rotation In Space	Rotation In Space
BACK HEIGHT	15 - 24"	21"
BACKREST ANGLE	-5° to +30°	-5° to +30°
FRAME / HANGER ANGLE	60° / 70° / 80° / 90° / Elevating	70°
INSET	N/A	N/A
CAMBER	0°	0°
MAX. USER WEIGHT	136 kg / HD 181 kg	120 kg
PRODUCT WEIGHT	17.6 kg	17.6 kg
TRANSIT APPROVED	Yes	Yes
FRAME COLOUR	25 colours	Black
WARRANTY	Limited Lifetime	Limited Lifetime
ARTG	100167	100167



# SEATING SELECTION

## JAY® CUSHIONS



LowZone ECO



LowZone Switch



LowZone Gel Switch



ClimaZone



X2™



Easy™



Ion™



Union™



Fusion™



Cryo Fusion™



J2™



J2™ Deep Contour

## JAY® BACKS



Jay 3 Carbon



J3™ Plus



Focus Point



J3™ Backs

To learn about personalising JAY seating products, visit the JAY Your Way product page on [www.SunriseMedical.com.au](http://www.SunriseMedical.com.au)



Sunrise Medical Australia Pty Ltd  
11 Daniel Street, Wetherill Park, NSW 2164  
Phone: (02) 9678 6600  
Email: [enquiries@sunrisemedical.com.au](mailto:enquiries@sunrisemedical.com.au)  
Website: [www.sunrisemedical.com.au](http://www.sunrisemedical.com.au)