

# SR45™



The smart choice for intelligent rotation,  
adjustability and durability

# SR45™

The best just got better. The Quickie® SR45™ is the newest member of the acclaimed Quickie tilt-in-space family and features 45 degree of patented and proven IRIS tilt technology and a targeted offering of clinically effective options.



## Seat Frame Set-up Options

When ordering the SR45™ you have a choice of frame set up by Quickie or selecting a specific depth range.

### Frame Set-up by Quickie:

Client anatomical measurements are used to determine the frame configuration necessary to optimize the center-of-gravity to center-of-rotation relationship.

### Selection of Specific Depth Range:

If you choose to select one of the 3 depth range options, the frame can be adjusted to 5 different width x depth configurations.

# Nothing beats a Quickie!



## Rotation-in-Space Technology

The new Quickie SR45™ seat frame rotates around the client's center of gravity providing the shortest possible wheelbase without compromising stability. When the center of gravity and center of rotation are properly aligned, tilting is smooth, easy and controlled.

## Smart Option Mix

The portfolio of options and accessories available on the SR45™ meets the needs of a wide range of clients. Newly redesigned Quickie compatible options include a swing-in/swing-out forward mount latch hanger and low profile dual post armrest.




Low Profile Dual Post Armrest



Swing-in/Swing-out Hanger




Steel Base

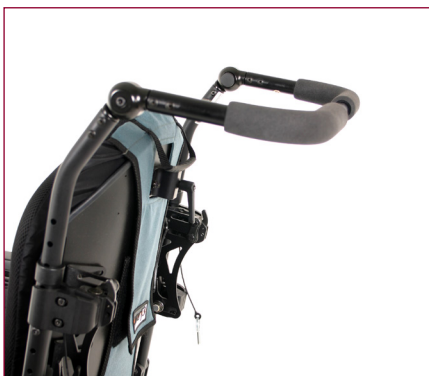

Adjustable Stroller Handle




Foot Release Tilt Actuator

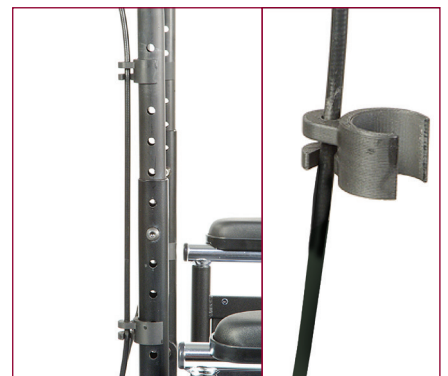

Cable Management



Short, small, light-easy to tilt.  
Unsurpassed maneuverability.

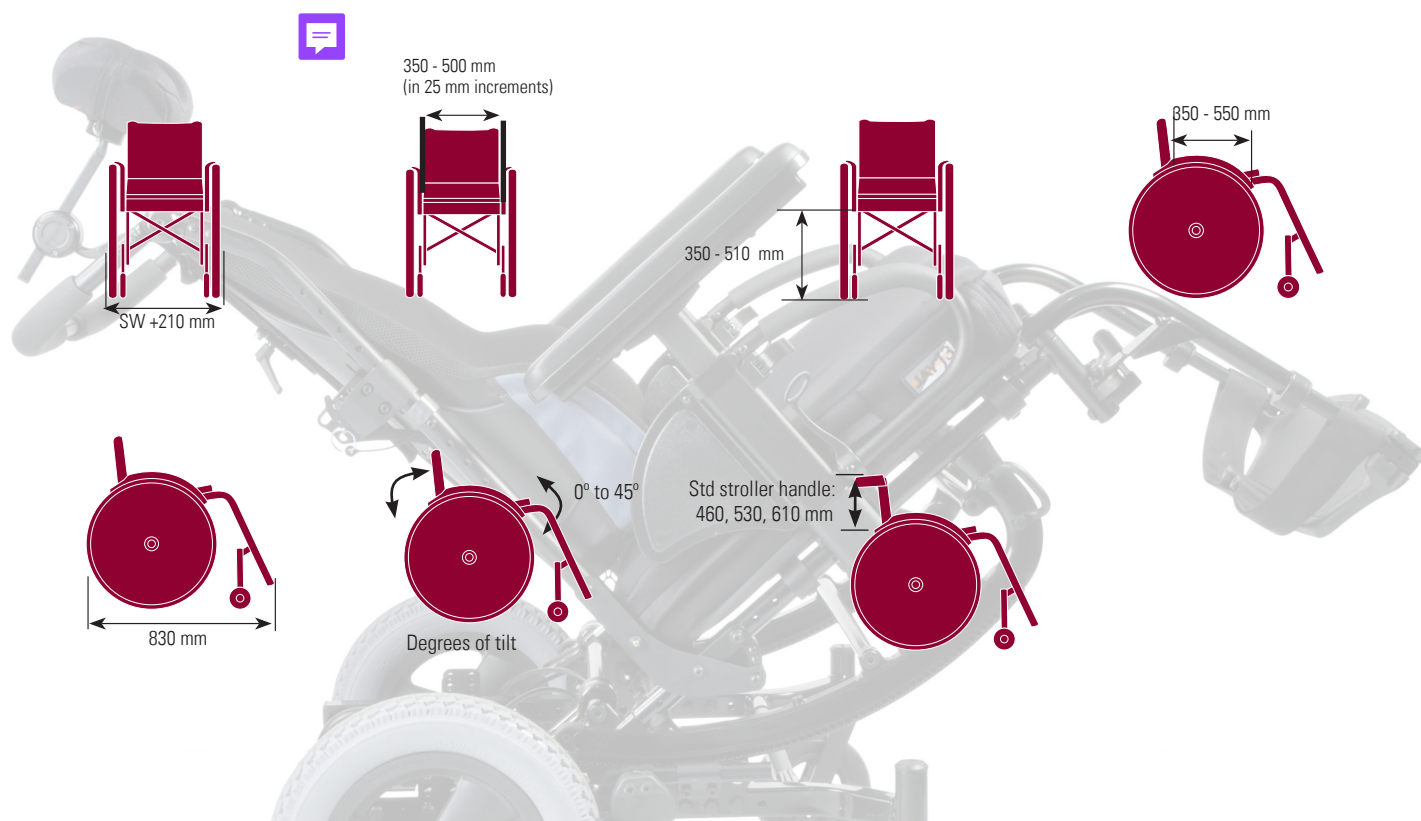

## Technical data sheet

Overall weight	from 17,6 kg	Frame style:	Tilt-in-space
Transport weight	from 13,6 kg	Back Options:	Stroller handle Back Post ( Low, Medium, Tall)
Max. user weight	120 kg	Hanger Options:	70°, 80° ( Swing Away and Heavy Duty) ELR

### IMPORTANT NOTE

Due to a recent update, the colours shown may not reflect current availability. For the most up-to-date selection, please refer to the latest order form or contact your local Sunrise Medical representative.

#### Sunrise Medical GmbH

Kahlbachring 2-4  
69254 Malsch/HD  
Tel.: +49 (0) 72 53 / 9 80-0  
Fax: +49 (0) 72 53 / 9 80-222  
E-Mail: [international@sunrisemedical.de](mailto:international@sunrisemedical.de)

For further information on the full specification, options, and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.



If you are visually impaired, this document can be viewed in PDF format at [www.SunriseMedical.eu](http://www.SunriseMedical.eu).



04.2025



[www.SunriseMedical.eu](http://www.SunriseMedical.eu)