

ComforFit®



Creating the right seating solution for your client often requires secondary supports. The ComforFit™ product line addresses a variety of positioning needs while providing exceptional comfort.

POSTURE SUPPORTS

POSITIONING BELTS

FOOT & ANKLE SUPPORTS


ComforFit®



A division of Sunrise Medical

Adaptive Equipment Systems | Kid Kart | RehabiliTech | ShurShape | Whitmyer Biomechanix

ComforFit Anterior Trunk Supports

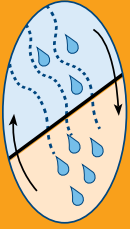

BREATHABLE FABRIC

Most ComforFit Trunk supports are constructed of a unique open-cell fabric which allows air to circulate through and moisture to escape, ensuring your client stays cooler and drier.

STRAP DESIGN

The design of the ComforFit upper D-ring adjustment straps keeps the buckles away from the user's face, head, and neck.

All ComforFit supports are made of an easy-to-care for fabric that is machine washable (air dry).



CLASSIC STYLE

The Classic style is designed with a large padded surface in the thoracic area for maximum support.

BREATHABLE



CONTOUR STYLE

Contour style features a narrowed sternum section especially suited for G-tube users and female clients, or to simply cover a smaller surface area.

BREATHABLE



BACKPACK STYLE

Allows control of forward flexion and allows more intimate shoulder retraction—can be used aggressively for upper torso control.

BREATHABLE



Y STYLE

Ideal for independent shoulder retraction and to assist in positioning where trunk rotation is an issue. Y-style supports allow for aggressive positioning while maintaining comfort and eliminates potential interference with feeding tubes.

NOT BREATHABLE



ZIPPER OPEN

Center zipper makes transfers quick and easy. Ideal for end-users with limited dexterity.

BREATHABLE

CENTER-OPENING W/ BUCKLE

Buckle closure is adjustable to allow for changes in the end-user's attire, or a limited amount of growth. Available in Structured only.

NOT BREATHABLE



ADJUSTABLE STRETCH STYLE

Great for when a client requires a varying amount of control throughout the day.

Lengthen or shorten the webbing straps that run down the front of the support to allow for more or less control as needed.

BREATHABLE



Foot & Ankle Support

Proper foot and ankle control is sometimes overlooked in the positioning process. There are several techniques available to the rehab professional depending upon the client's specific tone, pressure sensitivity, skin condition and lifestyle. ComforFit Foot & Ankle supports are available in a wide range of sizes to suit every client's need.



ANKLE POSITIONER


The Ankle Positioning System incorporates a soft, breathable, three-point ankle support that attaches to an ABS base plate. The base plate is secured to the footrest with the provided bolts (drilling the footrest may be required.) Adjustments may be made around the ankle or at the heel.



FOOT & ANKLE POSITIONER

The Foot & Ankle Positioning System incorporates all the features of the Ankle Positioning System, plus adds a ComforFit padded toe strap to ensure total control of the client's foot and ankle.

Toe strap is easily replaced—no sewing required!

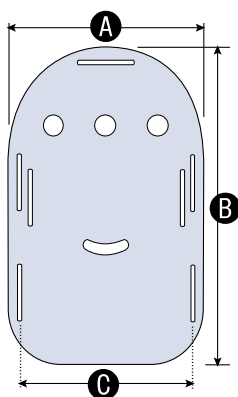

			BASE PLATE DIMENSIONS			ANKLE STRAP CIRCUMFERENCE
			A	B	C	
	ANKLE POSITIONER	FOOT & ANKLE POSITIONER				
XS	TST6501	—	4.5	6	4	5–7.5
S	TST6511	TST6510	5	7.75	4.5	6.5–9
M	TST6521	TST6520	5	7.75	4.5	8–10.5
L	TST6541	TST6540	6	9.75	5.5	9.5–12
XL	TST6561	TST6560	6.5	10.75	6	11–13.5




Shoe Holders

ComforFit Shoe Holders are designed to capture the complete shoe, or stockinged foot.

All Shoe Holders come standard with adjustable straps that have a curved buckle for added comfort. Removing and replacing straps is easy — no sewing is required for attachment.



Heel cup cutout provides clearance for Achilles tendon

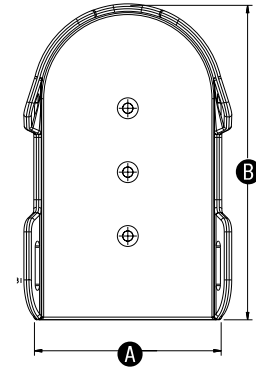

SHOE HOLDER W/ TENDON RELIEF

The Shoe Holder with Tendon Relief is constructed of an injection molded plastic in a single wall thickness.

High heel walls allow for more aggressive positioning, while a cutout at the back of the heel cup provides clearance for the client's Achilles tendon in plantar flexion.


Available only with padded straps.

	DYNAMIC	A	B
S	TST6061	3.5	5.25
M	TST6161	4	6.75
L	TST6261	4.75	8
XL	TST6361	5	10

SHOE HOLDER

Standard style Shoe Holders feature a low profile, double-wall construction with a rolled edge to eliminate hard edges.

Available with either padded or non-padded straps.

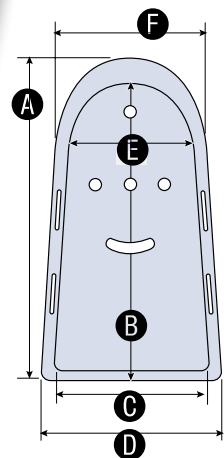

	PADDED STRAPS	NON-PADDED STRAPS
S	TST6001	TST6000
M	TST6101	TST6100
L	TST6201	TST6200
XL	TST6301	TST6300


A	B	C	D	E	F
5.75	5.25	3.25	4.375	3.25	4.25
7.75	7.25	3.5	4.75	3.25	4.25
9.75	9.25	4.5	5.75	4	5.25
12.25	11.75	6	7	5.25	6.25

DYNAMIC MODEL

Dynamic models provide excellent postural control, increased mobility and a soft fit to the body. The client moves easily within a limited range, yet gently returns to a neutral position when relaxed.

STRUCTURED MODEL


Structured models are as soft and comfortable as our other ComforFit products but with increased stabilizing support. They are perfect for clients needing positive control with little movement of the upper body.

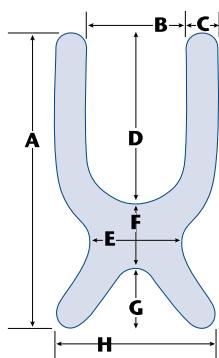

CLASSIC STYLE

	DYNAMIC	STRUCTURED
EI	ATS8107	—
XS	ATS8100	ATS8300
S	ATS8101	ATS8301
M	ATS8102	ATS8302
L	ATS8103	ATS8303
XL	ATS8104	ATS8304


A	B	C	D	E	F	G	H
9.75	3	2	5.25	3	3	1	6.5
12.5	3.5	2	7.5	3.5	3	2	7.5
15	4.5	2	9	4	4	2	8.5
19	5.5	2	11	5	5	3	9.5
22.5	6.5	2	13	6	6	3.5	10.5
25.5	6.5	2.5	14	7	7	4.5	11.5

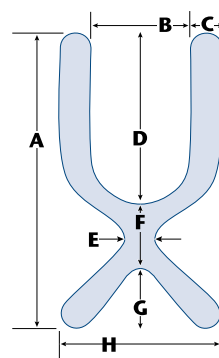

CONTOUR STYLE

	DYNAMIC	STRUCTURED
EI	ATS9107	—
XS	ATS9100	ATS9300
S	ATS9101	ATS9301
M	ATS9102	ATS9302
L	ATS9103	ATS9303
XL	ATS9104	ATS9304

A	B	C	D	E	F	G	H
9.75	3	1.75	5.25	1	3.5	1.5	6.5
12.5	3.5	2	7.5	1	3.5	1.5	7.5
15	4.5	2	9	1	4	2	8.5
18.5	5.5	2	11	1	5	2.5	9.5
22.5	6.5	2	13	1.5	6	3.5	10.5
25.5	6.5	2.5	14	1.5	7	4.5	11.5

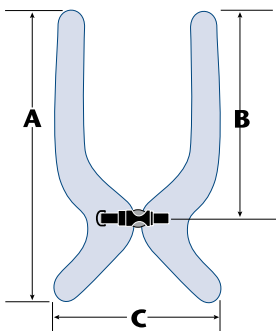


CENTER-OPENING WITH BUCKLE

	STRUCTURED	A	B	C
EI	ATS9407	9.75	7	5
XS	ATS9400	11	8	5.75
S	ATS9401	12.25	9	6.5
M	ATS9402	14.25	11	9.75
L	ATS9403	15.5	11.5	9.75
XL	ATS9404	17.5	12.75	10


NOTE: Dimensions are approximate and may vary due to binding finish. All measurements are listed in inches.

ComforFit Anterior Trunk Supports

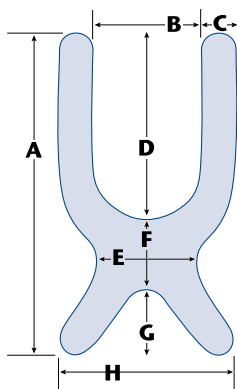


ZIPPER OPEN

	DYNAMIC	A	B	C	D	E	F	G	H
S	ATS8701	15	4.5	2	9	4	4	2	8.5
M	ATS8702	19	5.5	2	11	5	5	3	9.5
L	ATS8703	22.5	6.5	2	13	6	6	3.5	10.5
XL	ATS8704	25.5	6.5	2.5	14	7	7	4.5	11.5

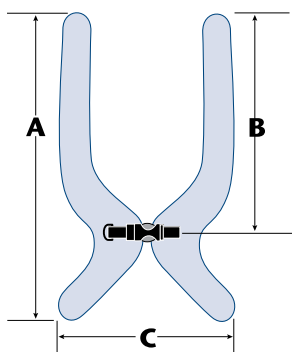


ADJUSTABLE STRETCH

	DYNAMIC	A	B	C
EI	ATS8807	9.75	7	5
XS	ATS8800	11	8	5.75
S	ATS8801	12.25	9	6.5
M	ATS8802	14.25	11	9.75
L	ATS8803	15.5	11.5	9.75
XL	ATS8804	17.5	12.75	10

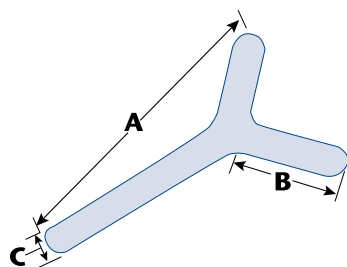

Y STYLE

	DYNAMIC	STRUCTURED	A	B	C
EI	ATS9167	ATS9367	15	7	1.75
S	ATS9161	ATS9361	19	8	2
M	ATS9162	ATS9362	22	9	2.25
L	ATS9163	ATS9363	25	11	2.5

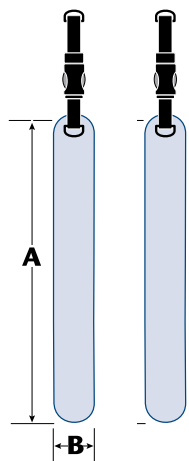

BACK PACK STYLE

	DYNAMIC	STRUCTURED	A	B
12"	ATS9601	ATS9701	12	2
18"	ATS9602	ATS9703	18	2

Pelvic Positioning Belts

ComforFit Pelvic Positioning Belts feature an ergonomic pad design. The small indent in the pad fits the ASIS for increased comfort while ensuring the positioning you intend to achieve.

Choose from a wide range of buckles and pad sizes. Metal anchors are included for mounting.



2-POINT BELTS W/ SIDE RELEASE BUCKLE

BUCKLE	STRAP	NON-PADDED	5.5" PAD	7.5" PAD
1" SR	Single Pull	APB6110	APB6111	APB6112
	Dual Pull	—	APB6121	APB6122

BUCKLE	STRAP	NON-PADDED	7.5" PAD	9.5" PAD
1.5" SR	Single Pull	APB6210	APB6212	APB6213
	Dual Pull	—	APB6222	APB6223

BUCKLE	STRAP	NON-PADDED	9.5" PAD	11.5" PAD
2" SR	Single Pull	APB6310	APB6313	APB6314

2-POINT BELTS W/ METAL PUSH-BUTTON BUCKLE

BUCKLE	STRAP	NON-PADDED	4.5" PAD	5.5" PAD	7.5" PAD
Li'l Kiddos™ 1" PB	Single Pull	APB5110	APB5115	APB5116	APB5117
	Dual Pull	—	APB5125	APB5126	APB5127

BUCKLE	STRAP	NON-PADDED	7.5" PAD	9.5" PAD
1.5" PB	Single Pull	APB7210	APB7212	APB7213
	Dual Pull	—	APB7222	APB7223

BUCKLE	STRAP	NON-PADDED	9.5" PAD	11.5" PAD
2" PB	Single Pull	APB7310	APB7313	APB7314

MEDIAL 4-POINT BELTS WITH METAL PUSH BUTTON BUCKLE

BUCKLE	3.5" PAD	4.5" PAD	5.5" PAD	6.5" PAD	7.5" PAD	9.5" PAD
Li'l Kiddos 1" PB	APB5145	APB5146	—	—	—	—
1.5" PB	—	APB7242	APB7243	APB7244	APB7245	APB7246

SR = Side Release
PB = Push Button (all metal)



To view other custom seating products in our range,
please visit: www.sunrisemedical.com.au or
please contact : marketing@sunrisemedical.com.au


JAY FIT SYSTEM

**Modular Adjustable
Contour Seating**

For Katie and her team of rehab professionals, the Jay Fit System is an effective way to give her the customized seating and positioning options she needs.

The Jay Fit System can keep everyone, including cowgirls like Katie, sitting tall in the saddle.



JAY 

JAY **SUREFIT**

Simply Custom Seating

