

# J3 CARBON

## FEEL THE DIFFERENCE



---

## IT'S SO LIGHT..

Carbon is lightweight and strong – therefore the ideal material for this backrest. Offering the same strength as aluminium (but 60 % lighter), the J3 Carbon back reduces the weight that a wheelchair user has to push, day-by-day. You will feel the difference!

## ... AND STYLISH

---

The J3 Carbon back will also impress you with its unique and beautiful finish. An active and clean design combined with an elegant appearance – the stage is truly yours.



# THE SUPPORT YOU NEED TO BE ACTIVE

---

The J3 Carbon back provides great stability and comfort. Through the right pelvis support it offers firm stabilisation of the upper trunk and encourages an optimal spinal curve. This enhances the stability you need for a very active life, reduces back pain and provides a great comfort. Plus, the force that's generated when propelling your wheelchair is transferred more effectively and efficiently through a solid back into the frame. That means you're generating greater movement of your wheelchair with less effort.

Just try it...



## For all types of chairs

Whether you use a rigid or folding frame, the J3 Carbon is your back. The fixed or quick release version fit on all types of chairs.

### → Quick & easy to remove


With the optional quick release hardware the backrest is quick and easy to release – ideal for folding chairs.



## Where are my keys?

You now have somewhere to store those easy to misplace items with the new accessory bag. Not only will it keep all of your valuables in one place, but it also smartly covers the gap between your back and cushion. You'll wonder what you ever did without it.

The usefully stylish accessory bag

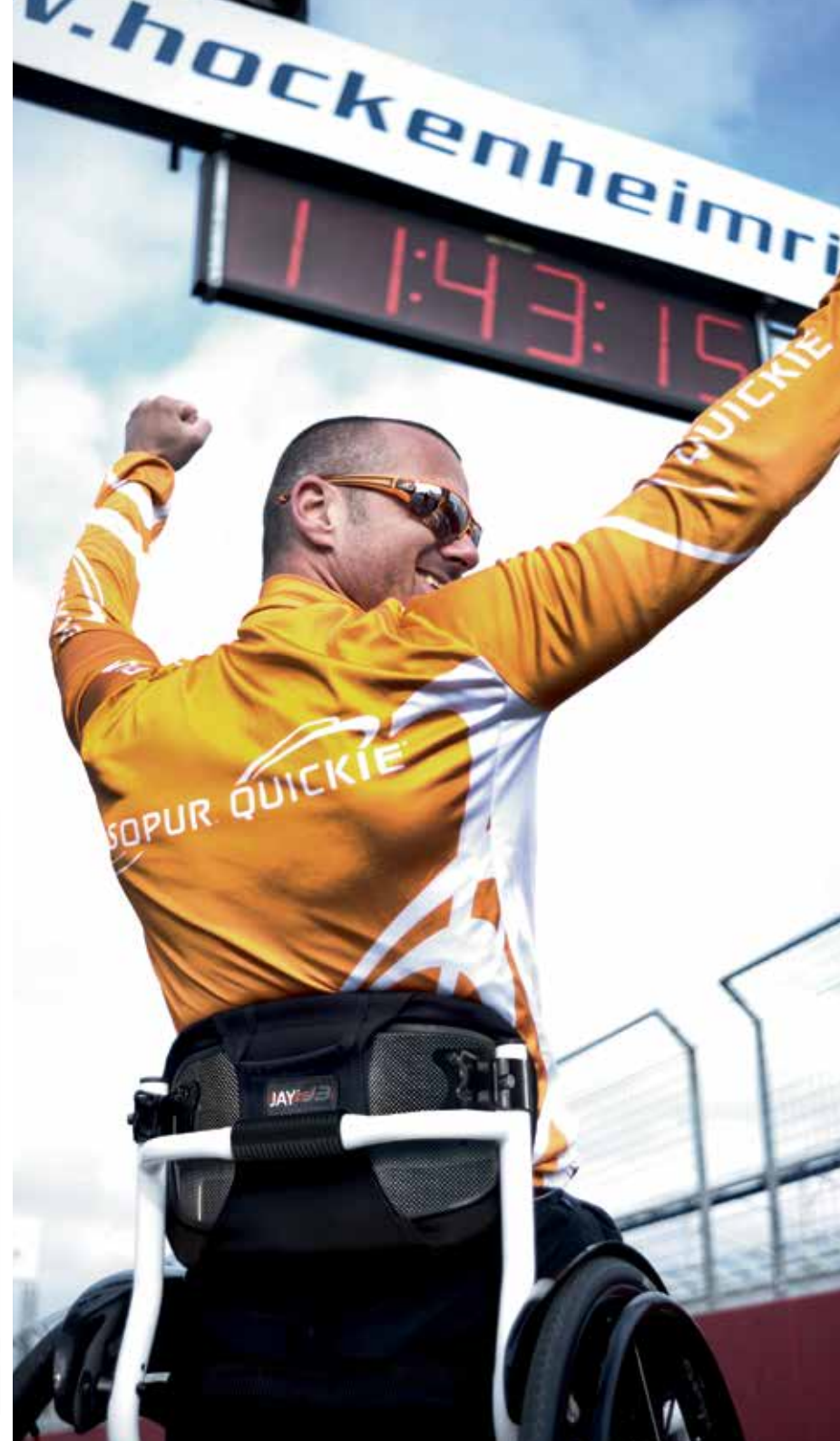

# ADJUSTABLE TO YOUR NEEDS


---

Two different mounting systems for folding  
and rigid chairs available

The J3 Carbon back is available with two different mounting systems. The fixed system offers rigid wheelchair users a mounting option that's extremely weight optimised. Alternatively, the quick-release system allows the backrest to be easily removed, ideal for mounting on a folding wheelchair. Regardless of choice, both mounting systems are designed as easy-to-handle, 2-point mounting systems and offer the same adjustments in angle, depth, height and width to help adapt the J3 Carbon back to your individual needs..



The quick release  
hardware  
ideal for folding  
wheelchairs

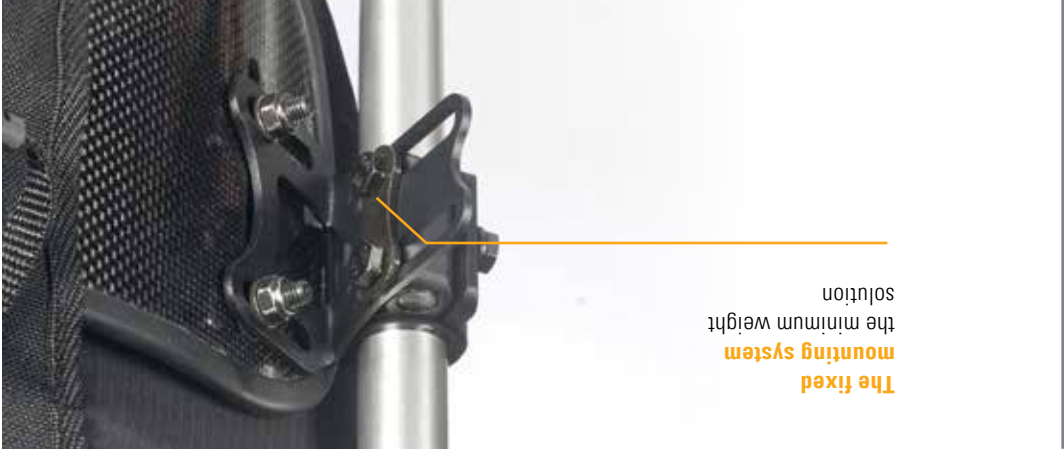

The fixed  
mounting system  
the minimum weight  
solution





→ **Active contour**  
Mild lateral trunk support  
through 5 cm contour depth



→ **Different support levels**  
Always the right support  
(four heights from 17 to 38 cm)

## Choose your back...!

Different needs demand adequate solutions. The four back heights (from 17 to 38 cm) focus on the active user with a need for low to moderate posterior trunk support. Mild lateral trunk support is also given by the 5 cm contour depth.



## Look cool and stay cool with the microclimatic cover and soft foam underlay

The J3 Carbon back's cover uses 3DX™ spacer fabric, which vents both heat and moisture and increases air flow through the back. Machine washable at 60°C and quick to dry, the cover overlaps the shell to provide you with extra protection. You can even adapt the shape of the cover to your individual needs with lumbar support pads.

→ **J3 Carbon cover**  
The 3DX™ technology offers a comfortable microclimate.



# TECHNICAL DATA

Max. user weight:	125 kg
Backrest width:	32 / 36 / 41 / 46 cm
Backrest height:	17 / 24 / 30 / 38 cm
Contour depth:	5 cm
Back angle:	± 20°
Backrest weight:	From 690 g
Hardware	Fixed and quick-release version
Accessories	Accessory bag and lumbar supports

For further information on the full specification, options and accessories, please refer to the order form.

All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at [www.SunriseMedical.eu](http://www.SunriseMedical.eu) or alternatively is available on request in large text.



**Sunrise Medical**

Kahlbachring 2-4

D-69254 Malsch/HD

Tel.: +49 (0) 72 53 / 9 80-0

Fax: +49 (0) 72 53 / 9 80-222

E-Mail: [international@sunrisemedical.de](mailto:international@sunrisemedical.de)

[www.SunriseMedical.eu](http://www.SunriseMedical.eu)

