

Patient Lifts

Sunlift



The SUNLIFT range consists of 4 models: MICRO, MINI, MIDI and MAJOR. All models feature exceptional rigidity thanks to the advanced design of all of the moving connections, eliminating play and offering a rigid feel with the hoist fully loaded. Actuators incorporate a double security button, light indicator of the battery charge level and an acoustic warning before total discharge happens. The battery capacity (2.9 Amph) provides a range of use at 50% higher than usual, allowing up to 60 double uses (up and down).



Sunlift Range




Technical Features

(common for the 4 models)

Batteries: 2 x 12v de gel
Batteries capacity: 2,9 Amph
Charger input: 230 v / 50 Hz 0.12 A
Charger output: 28 v 0.5 Amph

Measurements

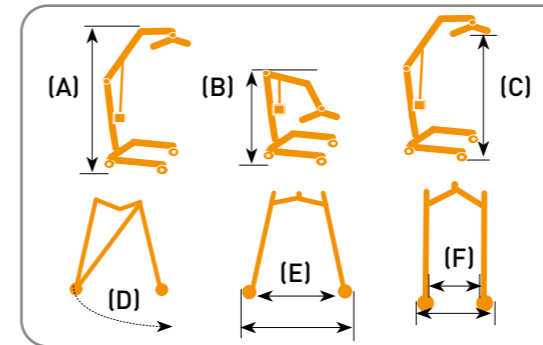

150 kg	175 kg
1.225 mm	1.340 mm
1.160 mm	1.275 mm
640 mm	650 mm
1.940 mm	2.006 mm
1.285 mm	1.435 mm
1.688 mm	1.813 mm
555 mm	575 mm
1.380 mm	1.500 mm
1.075 mm	1.185 mm
980 mm	1.105 mm
640 mm	650 mm
540 mm	520 mm
100 mm	110 mm
20 mm	15 mm
75 mm	75 mm
100 mm	100 mm
17,6 kg	22,75 kg
15,1 kg	18,1 kg
2,76 kg	2,76 kg

Common features



Legs can be opened and closed by pressing the ergonomic pedal.



Hand control classified with IPX4 (Water resistant) and ergonomic design.



Ergonomic push handles allowing elbows position at 90° for easier handling.

	SUNLIFT Micro	SUNLIFT Mini
Safe working load	130 kg	130 kg
Max total length (with wheels facing outwards)	1.010 mm	1.130 mm
Min total length (with wheels facing inwards)	970 mm	1.065 mm
Max width with legs closed	545 mm	545 mm
Max height (without load) (A)	1.725 mm	1.870 mm
Max height (without load) (B)	1.130 mm	1.280 mm
Max spreader bar height (without load) (C)	1.483 mm	1.628 mm
Min spreader bar height (without load)	460 mm	513 mm
Turning radius with legs open (D)	1.130 mm	1.240 mm
Legs open: (E)		
Max width (external)	870 mm	940 mm
Max width (internal)	740 mm	810 mm
Legs closed: (F)		
Max width (external)	545 mm	540 mm
Max width (internal)	465 mm	470 mm
Total leg height	100 mm	100 mm
Ground clearance	20 mm	20 mm
Castor wheels	75 mm	75 mm
Rear wheels	100 mm	100 mm
Mast, boom and electric components assembly weight	16,6 kg	17,54 kg
Base assembly weight	13,34 kg	13,6 kg
Batteries weight	2,76 kg	2,76 kg

Slings



	Quick fit sling	Quick fit sling deluxe	Access sling
Recommended use	85 % of patients	95% of patients	For easier access to toilet
Material	Polyester, Padded or Net	Padded	Polyester
Sizes	Medium or large	Medium	Medium
Independence level *	1 and 2	1 and 2	2 and 3
Others	-	With or without head support	-

	Long seat sling
Recommended use	Amputee patients
Material	Polyester
Sizes	Medium
Independence level *	1 and 2
Others	-



Selection of the most appropriate sling for each case:

	Toilet	Commode chair to bath	From chair to chair	From bed to chair	From floor to bed	From chair to commode chair
Quick fit sling	■	■	■	■	■	
Quick fit sling deluxe	■	■	■	■	■	
Access sling	■		■			■
Long seat sling	■			■	■	

*Level of independence 1 = Total independence 7 = Total independence



Sunrise Medical S.L

Polígono Bakiola 41
48498 Arrankudiaga
Vizcaya, España
Tel: 902 14 24 34
Fax: 94 648 15 75
www.SunriseMedical.es

All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have. If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.com or alternatively is available on request in large text

Slings