

A woman with long dark hair, wearing a brown suede jacket and blue jeans, is sitting in a black and red Quickie wheelchair on a city sidewalk. She is leaning against a wall on the left and looking towards the camera. The background is a blurred city street with other pedestrians and buildings.

**QUICKIE**  
Live without limits.

**NEON<sup>2</sup>**

**LIFESTYLE MEETS  
PERFORMANCE**



# NEON<sup>2</sup>

## MODERN. SMART. INDIVIDUAL.

---

Bursting with personality, the lightweight and affordable Neon<sup>2</sup> provides an energy-efficient ride that you'll never want to be without. The sleek, effortlessly foldable open-frame fits easily in the smallest car – and this is a wheelchair you'll never want to leave behind!

You are looking for individuality? Paint your personality all over the Neon<sup>2</sup> with a beautiful selection of colours for the frame, upholstery trims, wheels and forks. When you're done customising, it's not just the wheels that turn, heads will too.

Take a closer look at the new QUICKIE Neon<sup>2</sup> and we think you'll agree – you don't need to break the bank to enjoy a first class ride.




**// FINE-TUNE YOUR BACK UPHOLSTERY**

It's not just comfortable – Stand out from the crowd with a versatile pallet of 5 upholstery trims.



**// LIFE IS COLOURFUL**

Tailor your Neon<sup>2</sup> with a selection of 32 eye-catching colours.

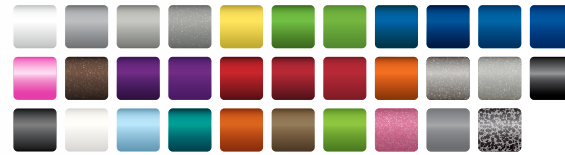


**// HEAD-TURNING LOOK**

There's also 5 distinctive styles to choose from for the castor fork and wheel.



**// THE SPECIAL ONE**

Accentuate your colour scheme further by matching the rim and hub of your 24" rear wheel with 3 anodized colour selections.

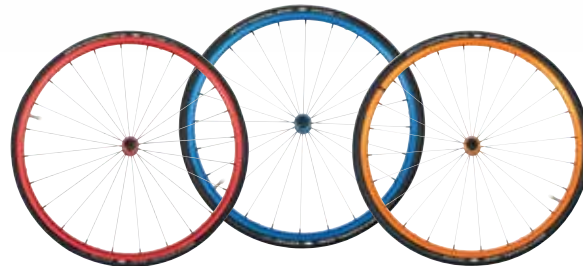

# NOTED FOR ITS **VERSATILITY**

---


At the heart of the Neon<sup>2</sup> is a reinforced open-frame. The firmness of the cross-brace and seat tubes provide an rigid-feel driving performance by reducing the flex of the wheelchair when in motion – all without sacrificing the benefits of a folding frame.

Speaking of folding, the unique open-rear of the frame allows for very compact space saving dimensions. The swing-away legrests not only ensure easy transfers in and out of your wheelchair, they're removable too. Combine it with the optional half-folding back and what you have is an extremely compact package that will fit in the smallest of cars.

The versatility of the Neon<sup>2</sup> doesn't end there. The durable aluminium axle-stem allows for simple adjustments to the rear seat height and centre-of-gravity. Need to adjust the backrest angle? It's just as easy via the aluminium backrest bracket – all without the need of being without your wheelchair.





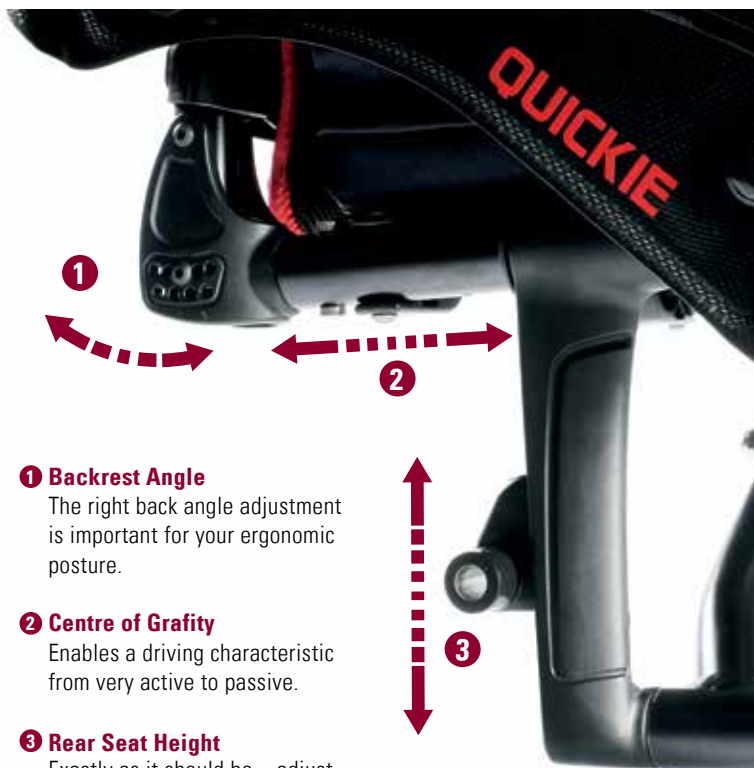

**1 Backrest Angle**

The right back angle adjustment is important for your ergonomic posture.

**2 Centre of Gravity**

Enables a driving characteristic from very active to passive.

**3 Rear Seat Height**

Exactly as it should be – adjust the height to your personal need.



**// PRACTICAL PORTABILITY WITH SWING AWAY FOOTRESTS**

They support easy standing transfers in and out of your wheelchair and offer very compact folded dimensions for easy portability by car.



**NEON<sup>2</sup>**

WHEN COMPACT  
DIMENSIONS AND EASY  
TRANSFERS COUNT.





## OPTIONS AND FEATURES OF THE NEON<sup>2</sup>

Configure your wheelchair to your individual needs. For further information on product options please check the Neon<sup>2</sup> order form.



### // ADD SOME COLOUR TO YOUR UPHOLSTERY

This very comfortable and smart backrest upholstery can be tailored by 5 different colour trims.



### // SEVERAL PUSH HANDLES

Choose the pushhandle which fits best for you. Decide between the integrated one (as shown) or the fold-down or height adjustable one.



### // ANGLE ADJUSTABLE BACK

Designer backrest solution integrated into the follow of the frame. Backrest angle from 75° to 103° in 4° increments.



### // ONE ARM FORK

This one arm fork especially in combination with the coloured aluminum forks gives your chair the special finish.



### // CARBON PLATFORM FOOTPLATE

Lightweight and stylish flip-up footplate (sideways) with angle and depth adjustment.



### // CARBON SIDEGUARD

Very individual and lightweight sideguard with a unique look.



### // FOLD-DOWN BACKREST

This half-folding backrest folds to the front and achieves a small wheelchair size for portability.



### // SURGE® LT HANDRIMS

Its Gripton rubber strip on top of the ergonomic oval shaped rim offers maximum traction.

# TECHNICAL DATA

NEON <sup>2</sup> // FOLDING WHEELCHAIR			
<b>Max. user weight:</b>	140 kg	<b>Backrest angle:</b>	75° to 103° in 4° increments
<b>Frame Version:</b>	Reinforced open frame with swing-away front	<b>Centre of gravity (COG):</b>	2 – 12 cm
<b>Seat width:</b>	30 – 50 cm	<b>Rear wheels:</b>	22" / 24" / 25" / 26"
<b>Seat depth:</b>	34 – 50 cm	<b>Camber:</b>	0° / 2° / 4°
<b>Seat height front:</b>	35 – 57 cm	<b>Chair weight:</b>	from 11,5 kg
<b>Seat height rear:</b>	33 – 53 cm	<b>Frame colours:</b>	32
<b>Backrest height:</b>	25 – 47,5 cm	<b>Crash tested:</b>	yes

For further information on the full specification, options and accessories, please refer to the order form.

All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at [www.SunriseMedical.eu](http://www.SunriseMedical.eu) or alternatively is available on request in large text.



## Sunrise Medical GmbH

Kahlbachring 2-4

69254 Malsch/HD

Tel.: +49 (0) 72 53 / 9 80-0

Fax: +49 (0) 72 53 / 9 80-222

E-Mail: [international@sunrisemedical.de](mailto:international@sunrisemedical.de)

[www.SunriseMedical.eu](http://www.SunriseMedical.eu)