

Gemino

Techi

100% ROLLATORS



GEMINO – 100% ROLLATOR FOR YOU

Are you looking for a rollator that fits your active lifestyle and allows you to walk comfortably and confidently, indoors and out? Then a Gemino rollator is the perfect choice for you. Our award-winning lightweight rollators set new standards in quality, comfort, safety and design. And let you be yourself!

The Gemino comes in various models, sizes and colours. Naturally, there's also a special/medical Gemino rollator for you, if you need something extra to walk safely. So there's always a Gemino that fits your needs - and personal style. Whatever your choice, you'll always get a 100% rollator. 100% comfort, safety and design.

Great when you can rely on something 100%





100% COMFORT



ENJOYABLE WALKING EXPERIENCE

Our Gemino rollators have been designed with your comfort and enjoyment in mind. A lightweight aluminium frame, high-quality materials and maximum stability ensure a great walking experience, indoors and out. Due to the height adjustment with memory function the ergonomic handles are always in the right position. And if it's time for a break, you can sit safely on the spacious seat. Naturally, you can also take your rollator with you when travelling further — it easily folds to a compact size.

Try the Gemino of your choice, and discover how easy and enjoyable walking can be!

// CHOOSE YOUR FRAME SIZE

Selecting the right size will not only help you walk comfortably and safely. Each frame size also has a different seat height - important for maintaining a correct posture when sitting (For more information, please see Technical Data on page 22).

// SMART PUSH HANDLE ADJUSTMENT

Find out how simple it is to adjust the ergonomic push handles to fit your height. They safely lock into place and due to the special memory function, they're always in the right position for you.

// EASY FOLDING WITH SELF-LOCKING SYSTEM

The Gemino's unique 'fold-click system' makes it very easy to fold and compact to store. Simply put the push handles in the lowest position to reduce the size further. Practical for small spaces and portability by car. By the way, did you know you can also safely manoeuvre in tight places when your rollator is (partially) folded?

// PRACTICAL EVERYDAY USE

With your Gemino you'll enjoy comfortable walking - and comfortable seating, if it's time for a break. Combine with a comfortable back support for additional support. Lots of things to carry? Simply use the practical basket – or choose a stylish bag.



WHEREVER YOU GO

Wherever you go, you'll feel confident with a Gemino rollator. Its sturdy, stable frame design ensures maximum safety and security in all settings. And the ergonomically designed push handles provide excellent support. But there's more. Many clever features are proof that we always go the extra mile for your safety. Gemino – enjoy your active live, in total safety!

// EXCELLENT, ULTRA-LIGHT BRAKING

Ergonomically designed, you'll find that operating the brake handles requires very little effort. The smart brake design is easy on the wheels and ensures optimal grip, even when the wheels show signs of wear.

// EASY NEGOTIATION OF OBSTACLES Kick tubes make it easy to climb kerbs or thresholds, while the protective side bumpers help you avoid obstacles.

// WELL VISIBLE IN THE DARK

Going out at night? 6 integrated reflectors on all sides make sure you're well visible in the dark. For even more safety, simply add the optional light.

// NEW STEP TIP

The stepping area is closer to the ground makes stepping up and down easier, reducing strain on your joints. The lightweight design ensures less effort when maneuvering, enhancing the overall mobility of your rollator.







GEMINO 20 GREAT FUNCTIONALITY AT AN ATTRACTIVE PRICE

Are you looking for a high-quality rollator that offers great value for money? Then the Gemino 20 is just right for you. Simple and safe to use, and it looks good too. With this lightweight rollator you'll enjoy easy and stable walking, while the 'fold-click system' allows for compact storage and easy transport. It's also fitted with a practical basket and spacious seat – great if you want to take a break.

The 'basic' Gemino 20 comes in 3 frame sizes (colour Grey). Choose the model that fits your height and enjoy 100% quality and functionality – at an attractive price.

// 100% GEMINO - SMART CHOICE FOR MOBILITY

Expect a real Gemino walking experience from your Gemino 20 rollator. It has all the technical features of the award-winning Gemino 30, and can be combined with lots of practical accessories.

// EASY STORAGE AND PORTABILITY

Your Gemino 20 simply folds to a compact package due to the unique 'fold-click system'. For practical storage and easy portability by car.

The new forks are essential for cutting down the rollator's overall weight and simplifying the wheel assembly. This makes changing wheels a breeze. GeMINOZO



GEMINO 30 ENJOY LIFE AND EXPRESS YOURSELF

Looking for a lightweight rollator that stands out? Gemino 30 - the original 100% rollator — will exceed your expectations when it comes to comfort, safety and design. Its modern design styling and high finishing express quality and class. Naturally, it allows for full personalization to suit your active life — and personal style!

With the award-winning Gemino 30, you'll get 'the best of the best'. Balancing high manoeuvrability with maximum stability, this design rollator will truly enhance your independence. Personalize with your choice of colour and accessories – and enjoy life to the full!

// AWARD-WINNING DESIGN

Our original 100% rollator is recognized by international design institutes like red dot and iF and German consumer magazine Öko-Test.

The new forks are essential for cutting down the rollator's overall weight and simplifying the wheel assembly. This makes changing wheels a breeze.



GEMINO 30 carbon

GEMINO 30 CARBON THE PREMIUM CARBON ROLLATOR FOR AN ACTIVE LIFESTYLE

May it be something special? Equipped with all the advantages of the Gemino 30, like Gemino's unique 'fold-click-system' and the huge range of accessories, the Gemino 30 Carbon is the choice for those who are looking for something special. Particularly light, particularly elegant and particularly maneuverable are the features that caracterizes our premium rollator and make it really stands out.

// SUPER LIGHTWEIGHT

The low weight of 5.4 (M)/5.6 kg makes the Gemino 30 Carbon easy to lift, easy to transport and store, and also super maneuverable in tight spaces.

// FREEDOM OF LEG MOVEMENT Parallel front tubes enable an increased seating space because the legs do not hit the tubes.

// IN FRAME INTEGRATED STEBTUPE The steptube was unobstrusively integrated into the carbon frame and allows to easily overcome curbs.







// GEMINO 30 CARBON COMFORT -COMFORT WHEELS OPTIONAL AVAILALBE

With soft wheels that offer premium comfort and superior damping to reduce impact and vibration - even on uneven grounds. With a modern design, they provide a stylish, updated look while ensuring a smoother ride.



IV.

// ANATOMICAL SHAPED HANDGRIPS

Ergonomically designed hand grips provide extra comfort and a secure feel and support while walking due to the extra hand support.



// CARBON FRAME

Carbon is characterized by extreme strength and stiffness combined with low weight and this makes the Gemino 30 Carbon the perfect companion for an active lifestyle



// ADJUSTABLE BACK SUPPORT

The adjustable back support allows to lean back and relax. The built-in seat together with the back support allow for a comfortable rest when needed.



// EASY HANDLING OF PURCHASES

A lot to carry? The advanced net bag can be easily opened with one hand, offers plenty of space for purchases and carries up to 5 kg.







GEMINO 30 COMFORT HIGH INDOOR MANOEUVRABILITY WITH OUTDOOR PERFORMANCE

Are you looking for a lightweight rollator that allows you to walk comfortably and confidently, in- and outdoors? Then a Gemino 30 Comfort rollator is the perfect choice for you. Compact with large soft wheels and superb rolling characteristics.

// WHEELS THAT MAKE THE DIFFERENCE

Strolling the cobblestone streets or walking down the forest path. The-large, soft, shock absorbing wheels allows you to walk comfortable – on any surface.

// INDOOR AND OUTDOOR USE

The compact frame gives you high inndoor manouevrability and the large soft wheels allows you to negotiate obstacles wiht ease whether it is indoor tresholds or outdoor curbs.

// EASILY HANDLING EXTRA WEIGHT Lots of things to carry? The extra strong basket carries up to 10 kg.









Gemino 50

// DESIGNED FOR SUPERIOR STABILITY

The Gemino 60's sturdy, extra wide and long frame together with the large wheels ensure excellent stability and safety, also on rough terrain.

// WHEELS THAT MAKE THE DIFFERENCE

The extra-large, soft wheels allow you to negotiate virtually any obstacle with ease. And ensure ultimate walking comfort – on any surface.



// BASKET WITH LID

Unlike net alternatives, the new weather-resistant basket guarantees no issues with rain or dampness. The basket can be fully closed preventing unauthorized access and adding an extra layer of security to your belongings.



// ENJOY MORE SPACE

The large wheelbase makes for maximum stability, whilst offering extra walking space and a very spacious seat – also available with adjustable comfort back support. Great if you want to take a break.



// EASILY HANDLING EXTRA WEIGHT

All Gemino 60 models are suitable for heavy-duty use and take a maximum user weight of 150 kg. Lots of things to carry? The extra strong basket carries up to 10 kg.



GEMINO 60 ULTIMATE COMFORT AND SAFETY – ON ANY TERRAIN

Let's go! Gemino 60 is the perfect companion for your active lifestyle – lightweight yet rock-solid. With its extra-large, soft wheels this design rollator allows you to walk comfortably and securely, also on uneven terrain. Strolling the cobblestone streets in the old town, enjoying a walk in the forest or simply negotiating thresholds in your home? It'll all be so much easier with the robust Gemino 60!

The Gemino 60 is designed for active outdoor use and easily copes with extra weight. At the same time, it quickly folds together – for easy storage and portability by car. Gemino 60 - stay active and enjoy ultimate comfort, wherever you go!

// GEMINO 60 ATR - ALL TERRAIN WHEELS

The all terrain wheels (ATR) are designed for use on uneven ground, making it ideal for winter conditions and for usage on any terrain. Enjoy increased comfort and safety with enhanced grip and damping.







The new forks are essential for cutting down the rollator's overall weight and simplifying the wheel assembly. This makes changing wheels a breeze.

EXTRA SUPPORT AND STABILITY. ALSO GREAT FOR GAIT TRAINING

If you need extra support when walking, the Gemino 30 Walker is the perfect choice for you. This lightweight rollator is fitted with fully adjustable forearm supports, enabling you to walk comfortably and securely. The integrated push handles and brakes are easy to reach and very light to operate. Combining the benefits of the Gemino 30 with enhanced support and stability makes the Gemino 30 Walker also ideal for gait

training.

// ADVANCED BRAKING OPTIONS FOR EVEN MORE SECURITY

Looking for even more security? Complement your Gemino 30 Walker with a slowdown brake or innovative SpeedControl braking system. With SpeedControl your rollator will brake automatically when you start to accelerate. You'll keep an upright posture and walk confidently and comfortably close to your rollator.

// FOREARM SUPPORTS FOR ENHANCED STABILITY AND SAFETY

You can easily adjust the height, depth and angle of the forearm supports to suit your needs. The brakes and push handles are fully integrated – easy to reach and light to operate. For comfortable and safe walking.



// NEW LASER FUNCTION WALKER

The innovative laser function is now available for our Walker, specifically designed for patients with Parkinson's. This feature provides visual cues to help initiate movement, improving gait and reducing the risk of freezing episodes. When freezing occurs, simply activate the laser line with a button on the push handles. By focusing on the line that appears on the ground, you can step over it and resume moving with ease.





OPTIMAL SUPPORT AND STABILITY – INDOORS AND OUT

Going out a lot? With the Gemino 60 Walker you get the best of both worlds: The comfortable rolling characteristics and stability of the Gemino 60 and the additional support and safety you'd expect from a 'walker'. Adjust the forearm supports and push handles to fit your needs, and enjoy more freedom indoors - and outdoors. Naturally, you can also use this lightweight rollator for extended gait training!

// ADJUSTABLE FOREARM SUPPORTS FOR ENHANCED STABILITY

The height and angle adjustable forearm supports with integrated push handles and brakes are designed for optimal support and stability, and easy use. Adjust to suit your height and enjoy a comfortable walking position.

// COMFORTABLE AND SAFE OUTDOOR USE

The sturdy frame design coupled with extra-large soft wheels allow you to go out in total comfort and safety. Cobbled streets or thresholds in your home? With the Gemino 60 Walker you'll overcome any obstacle with ease.





GEMINO 30 SpeedControl

GEMINO

// ALSO A PERFECT COMPANION FOR CHILDREN

The Gemino 30 S SpeedControl, the smallest model, offers children every opportunity to discover the world. Thanks to the dynamic braking system they enjoy total safety when walking – making gait training fun!





the rollator's overall weight and simplifying the wheel assembly. This makes changing wheels a breeze.

GEMINO 30 SPEEDCONTROL THE ROLLATOR INNOVATION

You're looking for a rollator that offers maximum security when walking? Then the Gemino 30 SpeedControl with its innovative dynamic braking system is just what you need. You'll feel completely secure and enjoy a smooth gait and better, upright position. SpeedControl helps prevent falls and is also great for gait training – for adults and children.

// 3-LEVEL DYNAMIC BRAKING FOR MAXIMUM SECURITY

The innovative SpeedControl braking system features 3-level breaking and can be easily adjusted to your needs. When you start to accelerate, the centrifugal brakes in the rear wheels will automatically be activated. Naturally, you'll experience no resistance when walking at low speed or up slopes. SpeedControl promotes a steady pace, ensures that you walk comfortably, close to the rollator and prevents the rollator from 'running away' from you.

// GREAT AID FOR GAIT TRAINING

It's the Gemino 30 SpeedControl's unique ability to smooth out your gait through offering controlled resistance that also makes it a great training partner if you need to learn to walk again. You can concentrate fully on the sequence of movements, while your rollator promotes an upright posture and a steady pace. And use the different braking levels for adjusting the rollator to the progress you are making in your rehabilitation.

// NEW STEP TIP

The stepping area is closer to the ground makes stepping up and down easier, reducing strain on your joints. The lightweight design ensures less effort when maneuvering, enhancing the overall mobility of your rollator.





SPECIALLY DESIGNED FOR PEOPLE WITH PARKINSON'S

Some people need a little extra help to be able to walk safely. If you've Parkinson's disease, you may be one of those people. Our Gemino 30 Parkinson rollator has been developed with you in mind and offers features that are fully tailored to your particular needs. Get out and about safely – and with confidence!

// SLOWDOWN BRAKE PROMOTES A STEADY PACE

The Gemino 30 Parkinson is fitted with a slowdown brake, which provides additional resistance and helps you walk at a steady pace. Adjust the level of resistance to your specific needs.

// INCREASED SAFETY WITH REVERSE BRAKE

Due to the reverse braking system the rollator immediately stops moving once the brake is released. This increases safety if you experience 'freezing' while walking - when it feels like your feet are glued to the ground.

// LASER DEVICE HELPS OVERCOME FREEZING

When 'freezing' occurs, you can activate a laser line with a button on the push handles. Through concentrating on the line that appears on the ground, you can then step over it and get moving again.

// NEW STEP TIP

The stepping area is closer to the ground makes stepping up and down easier, reducing strain on your joints. The lightweight design ensures less effort when maneuvering, enhancing the overall mobility of your rollator.

The new forks are essential for cutting down the rollator's overall weight and simplifying the wheel assembly. This makes changing wheels a breeze.

Gemino 30 comfort parkinson

Left or right mount

ENHANCED MOBILITY FOR PARKINSON'S PATIENT

Discover the Parkinson's version of our popular Gemino 30 Comfort rollator – maintaining its maneuverability while still offering compactness and lightness. Experience superior stability and smooth rolling characteristics with large, soft wheels that provide a comfortable, secure ride even on uneven terrain. Plus, the features that give you that little bit of extra support to walk safely when you need it.

// 20% MORE COMFORT THAN USUAL WHEELS!

Walk with ease on any surface with our Gemino 30 Comfort Parkinson rollator. Equipped with large, soft, and shock-absorbing wheels, it ensures a comfortable stroll and effortless navigation over obstacles.

// ADJUSTABLE SLOWDOWN BRAKE PROMOTES A STEADY WALKING PACE

The lightweight GEMINO 30 Comfort Parkinson rollator is fitted with a slowdown brake, which provides additional resistance and helps you walk at a steady pace. Need more, or less, resistance? The level of resistance can easily be adjusted to your specific needs.

// LASER DEVICE HELPS OVERCOME FREEZING

In moments of freezing, the Gemino 30 Comfort Parkinson rollator supports you with a built-in laser line. A simple press of a button on the push handles activates the laser line, allowing you to focus on it as it appears on the ground. Step over with confidence and get moving again effortlessly. The new forks are essential for cutting down the rollator's overall weight and simplifying the wheel assembly. This makes changing wheels a breeze.





// INCREASED SAFETY WITH REVERSE BRAKING SYSTEM

The reverse braking system is another great feature of the GEMINO 30 Comfor Parkinson: The rollator immediately stops moving once the brake is released. This increases safety for you if you experience ,freezing' while walking - when it feels like your feet are glued to the ground.

// NEW STEP TIP

The stepping area is closer to the ground makes stepping up and down easier, reducing strain on your joints. The lightweight design ensures less effort when maneuvering, enhancing the overall mobility of your rollator.

100% CHOICE PRACTICAL ACCESSORIES FOR YOU

Gemino rollators are not just easily adjustable. They also come with a full range of attractive accessories to make everyday life easier. A stylish bag, a practical crutch holder, a comfortable back support or a light for extra safety? Choose the options that fit your needs and lifestyle!

// COLOURFUL ADVANCED BAGS Compliment your Gemino with a stylish, closable bag in your favourite colour. With handles and shoulder strap for easy carrying.

Not available for Gemino 60 M





Graphite grey





Pink



Mint green







The Shop-n-Go trolley bag is the perfect accessory for your rollator, making daily errands and trips a breeze. Easily store and protect your groceries in the trolley bag, take them home safely with your rollator, and remove the bag effortlessly with a single movement thanks to the clever fastening mechanism. With smooth-running wheels and an extendable handle, you can pull the bag effortlessly, perfect for negotiating single stairs or when the rollator is not to be taken into the living area. The Shop-n-Go is compatible with the latest models Gemino 20, Gemino 30 and Gemino 60 and can also be retrofitted to older models.







Red

// BAGS FOR GEMINO 30 S





Alcantara black

Grey, black print



// UMBRELLA****

Rainy weather conditions? With this height-adjustable umbrella you can go out safely, as you can keep both hands on the push handles.



// BACK SUPPORT* Provides you support to sit comfortably and safely when it's time to take a break. No need to remove when folding the rollator.



// ADJUSTABLE COMFORT BACK SUPPORT**

For additional comfort, choose this large, softly padded back support. Quick to adjust, easily folds with the rollator.



// TABLE TRAY Practical tray with cup holder

allows you to transport food and beverages safely and with ease.



// CRUTCH HOLDER Keeps your crutch or walking

cane securely fastened when walking. Easy to mount, easy to access.



// BASKET FOR OXYGEN CYLINDER***/**** Carrier bag with shoulder strap for safe transport of oxygen cylinders.



// SLOWDOWN BRAKE**** Provides adjustable resistance for steady walking. Easy to turn on/off.



// UNIVERSAL BRACKET Suitable for fitting your choice of personal accessories.



// LIGHT Ensures you're well visible in the dark – for additional safety.



// INFUSION HOLDER**** Height and angle adjustable drip holder.





Automatic braking at times of sudden accelaration

// SPEEDCONTROL*

With the innovative SpeedControl braking system – available with all Gemino 30 models – you'll feel completely secure when walking. Your rollator will brake automatically when you start to accelerate, ensuring smooth walking and a comfortable, upright position at all times. Naturally, you'll experience no resistance when rolling at normal pace or up slopes. SpeedControl helps prevent falls and is also great for gait training - for adults and children!



// ONE HAND BRAKE Comfortable and safe usage with one hand.



// REVERSE BRAKE KIT The rollator immediately stops moving when the brake is released.



// STABILIZER BAR More safe feeling, easy removable for transport. Only available for Gemino 30 and Gemino 60 Walker

* check availability per model


** not available for Gemino 30 / 30 M Walker

*** not available for Gemino 60 M

**** not available for Gemino 30 Carbon / 30 M Carbon



GEMINO – SELECTION GUIDE

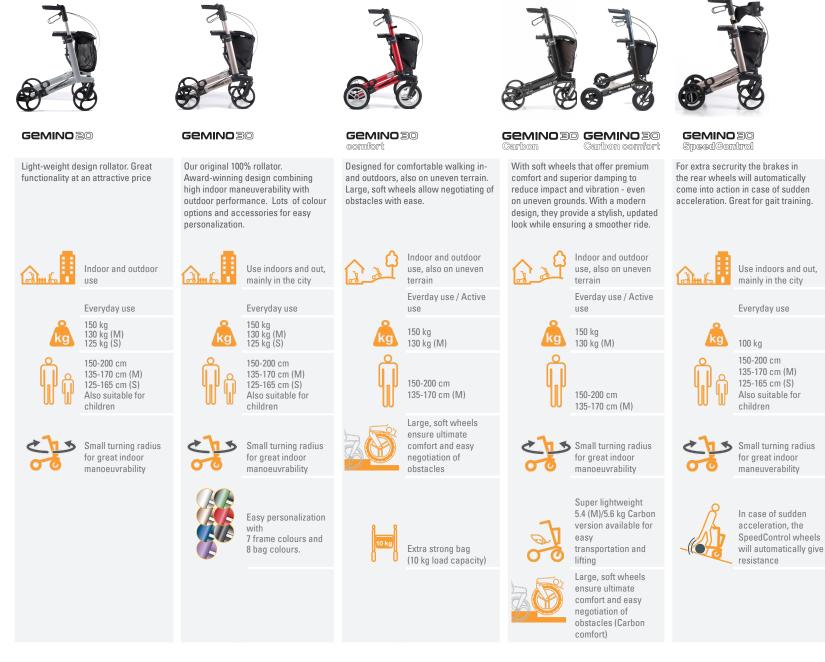

CHOOSE THE GEMINO ROLLATOR THAT SUITS YOU



GEMINO30 parkinson

Specially designed for people with Parkinon's. Featuring a slowdown brake promoting a steady pace, and a reverse braking system and laser light to help with 'freezing'.





Laser light can be activated to help overcome 'freezing'



GEMINO30 comfort parkinson

Specially designed for people with Parkinson's, this innovative solution features a slowdown brake that promotes a steady pace, a reverse braking system, and a laser light to combat 'freezing.' It's also designed for comfortable walking on uneven terrain.



150-200 cm 135-170 cm (M)

Large, soft wheels ensure ultimate comfort and easy negotiation of obstacles

Reverse braking system: the rollator immediately stops moving once the (P)brake is released



Laser light can be activated to help overcome 'freezing'



Gemino 50

Designed for any kind of terrain and heavy-duty use. Lightweight, yet rocksolid. All the benefits of the awardwinning Gemino 30 and more.

> Mostly outdoor use, also on rugged terrain, occasionally indoors

> > Active / heavy-duty use

GEMINO 60

The all terrain wheels (ATR) are

designed for use on uneven ground,

making it ideal for winter conditions

and for usage on any terrain. Enjoy

increased comfort and safety with

150 kg

150 kg (M)

150-200 cm 135-170 cm (M)

Mostly outdoor use,

also on rugged terrain,

Active / heavy-duty use

Extra wide and long

Extra-large and soft

comfort and easy

negotiation of any

Extra wide seat

obstacle

wheels ensure ultimate

frame for ultimate

stability

occasionally indoors

enhanced grip and damping.

ATR



150-200 cm 135-170 cm (M)



Extra-large and soft wheels ensure ultimate comfort and easy negotiation of any obstacle



Extra strong bag (10 kg load capacity)



GEMINO30 walker

Lightweight walker with adjustable forearm supports. The perfect companion if you need extra support and stability indoors. Also great for gait training.






GEMINO 50 walker

Optimal support and stability – indoors and out. Combines the comfortable rolling characteristics and stability of the Gemino 60 with the additional support and safety you'd expect from a walker.





GEMINO TECHNICAL DATA

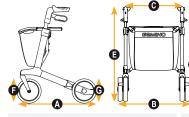








Gemino30 GEMINO 30 Carbon COMINIC 20 Carbon comfort **GEMINO 30** Comfort

					GEMINO 30 Carbon comfort				
		20	20 M	205	30	EOM	305	30	BOM
DIMENSIONS							not Carbon		
Maximum user weight	kg	150	130	125	150	130	125	150	130
Recommended user height	cm	150-200	135-170	125-165	150-200	135-170	125-165	150-200	135-170
Maximum load capacity in basket	kg	5	5	5	5	5	5	10	10
Total length \Lambda	cm	65	65	68	65	65	68	68	68
Total width 🚯	cm	60	60	55	60	60	55	62	62
Total weight without basket (with basket)	kg	7.1 (7.3)	7.0 (7.2)	6.8 (6.6)	7.1 (7.4) 5.6 (5.9) Carbon	7.0 (7.3) 5.4 (5.7) Carbon	6.6 (6.8)	7.7 (8)	7.5 (7.8)
Width between centre of push handles O	cm	47	47	42	47	47	42	47	47
Seat height D	cm	62	55	48	62	55	48	64	57
Total height 🕒	cm	78 - 100	69 - 88	65 - 77	78 - 100	69 - 88 73 - 88 Carbon	65 - 77	80 - 102	71 - 90
Total height folded	cm	80	76	69	80	76	69	82	77
Total lenght folded	cm	65	65	68	65	65	68	68	68
Total width folded	cm	23	23	23	23	23	23	23	23
Diameter front wheels 🕞	cm	20 x 3.6	20 x 3.6	20 x 3.6	20 x 3.6	20 x 3.6	20 x 3.6	23 x 4.2	23 x 4.2
Diameter rear wheels G	cm	20 x 3.6	20 x 3.6	20 x 3.6	20 x 3.6	20 x 3.6	20 x 3.6	23 x 4.2	23 x 4.2
PERFORMANCE									
Turning radius	cm	84	84	79	84	84	79	88	88
COLOUR OPTIONS									
Black						for CarbonM only			
Silver grey					√ not Carbon	√ not Carbon	\checkmark	\checkmark	\checkmark
Champagne					√ not Carbon	√ not Carbon			
Sapphire blue					√ not Carbon	√ not Carbon	\checkmark		
Pink					√ not Carbon		\checkmark		
Apple green					√ not Carbon	√ not Carbon			
Scarlett red					√ not Carbon				
CERTIFICATES									
CE		√	√	V	√	√	√	√	V
EN 12182		√	1	V	√	√	√	√	V
ISO 11199-2		V			\checkmark		\checkmark		V
ISO 11199-3									

CE Gemino Walker rollators: EN 12182, ISO 11199-3 - All other Gemino rollators: EN 12182, ISO 11199-2

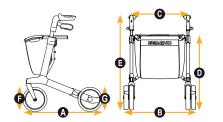











GEMINO30



GEMINO 30

Walker Walker SpeedControl comfort parkinson **GEMINO** GOATR Parkinson 50 50M IO EOM 50 50M IC 30 M 305 IO 30 M IO 30**M** Maximum user weight 150 150 130 130 130 100 100 100 130 130 150 130 kg 150 Recommended user height cm 150-200 135-170 150-200 135-170 150-200 135-170 150-200 135-170 125-165 150-200 135-150-200 135-170 Maximum load capacity in basket kg 10 Total length \Lambda 68 cm 74 74 66 68 Total width **B** 64 62 62 cm 64 60 60 64 64 60 60 55 60 60 Total weight without basket (with basket) kg 8.5 (8.8) 8.2 (8.5) 10.6 (10.9) 10.4 (10.7) 11.8 (12.1) 11.6 (11.9) 8.7 (8.9) 8.6 (8.8) 8.2 (8.4) 7.5 (7.8) 7.4 (7.7) 8.1 (8.4) 7.9 (8.1) 47 Width between centre of push handles O cm 51 51 41 41 45 45 47 47 42 47 47 47 Seat height D 57 cm 62 53 62 55 62 53 62 55 62 55 64 48 Total height 🕒 76-95 73 - 83 98 - 115 84 - 100 98 - 115 89 - 98 78 - 100 65 - 77 78 - 100 69 - 88 80 - 102 71 - 90 69 - 88 cm 77 Total height folded 80 80 82 cm 84 97 89 109 101 76 69 80 76 74 74 65 74 74 65 65 68 68 Total lenght folded 65 68 65 cm 23 23 23 cm 23 23 35.5 35.5 35.5 35.5 23 23 23 23 Total width folded Diameter front wheels cm 28 x 4.2 28 x 4.2 20 x 3.6 20 x 3.6 28 x 4.2 28 x 4.2 20 x 3.6 20 x 3.6 20 x 3.6 23 x 4.2 23 x 4.2 20 x 3.6 20 x 3.6 Diameter rear wheels **G** 23 x 4.2 23 x 4.2 20 x 3.6 20 x 3.6 23 x 4.2 23 x 4.2 20 x 3.6 23 x 4.2 23 x 4.2 cm PERFORMANCE Turning radius cm 84 84 95 84 79 8/ 84 Black Silver grey Champagne Sapphire blue Pink Apple green Scarlett red CE EN 12182 ISO 11199-2 ISO 11199-3

For further information on the full specification, options and accessories, please refer to the order form.

All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.





Sunrise Medical GmbH

Kahlbachring 2-4 69254 Malsch/HD Tel.: +49 (0) 72 53 / 9 80-0 Fax: +49 (0) 72 53 / 9 80-222 E-Mail: international@sunrisemedical.de www.SunriseMedical.eu For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have. If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.eu.

