



Explore Sunrise Medical
mobility bases and Leckey
modular seating systems to
find the perfect combination

Supporting Children

As part of the Sunrise Medical family, Leckey is delighted to be able to offer a range of Leckey seating systems on Sunrise Medical mobility bases. From early intervention to older children with more complex needs, the perfect combination can now be achieved.



Squiggles Seat with ZIPPIE Iris

Perfectly designed for the unique needs of young children, the much-loved **Squiggles Seat** helps with everyday functions at home, in the classroom or even out and about.

- Squiggles is not a scaled down version of another seat. Importantly, kids are not over-supported in Squiggles, helping them to reach their full developmental potential.
- The one-piece seat base and adjustable sacral support can be used to accommodate flexion or encourage extension which is pertinent at this early stage. The optional contoured shoulder support is angle adjustable and helps bring the hands to midline.
- Squiggles provides a snug, comfortable fit which kids love. This is achieved via the shape of the seat and back, the laterals with their fabric tip which wrap securely around the child and the way in which the 4-point pelvic harness integrates with the hip pad.





- With its 40° or 55° rotation in space ranges, **Sunrise Medical's ZIPPIE Iris** wheelchair helps children remain comfortable and secure whilst tilting, and is easy for parents and caregivers to control.
- An optional folding frame with patented XLOCK folding crossbrace for extra security allows you to easily transport the base in a compact size.
- The Iris is an intuitive, easy to adjust base with built in growth options to accommodate your growing child's needs.
- The Iris base also comes with an adjustable COG setting, allowing you to optimise the position of the base in relation to the Leckey seat. This ensures that you can always achieve maximum stability and ease of propulsion. It also allows you to adjust this setting so as the child and the seat grow, the base is still set up correctly.

Mygo Seat with QUICKIE Iris

Mygo Seat provides unrivalled postural support and comfort, helping more children with moderate to complex needs optimise their independence and opportunities for skill development.

- The highly adjustable sacral support enables fixed pelvic rotations and obliquities, as well as both anterior and posterior tilts, to be accommodated. In combination with a ramped seat base and 4-point pelvic harness accessory, Mygo Seat is able to mimic the support a therapist's hands would offer.
- Mygo Seat's innovatively designed split seat base enables the legs to move independently beyond the midline to accommodate windsweeping postures. A range of up to 13° abduction and 10° adduction may be achieved whilst offsetting the hip guides can aid in accommodating more extreme windsweeping by providing additional space. This ensures a more comfortable, stable and functional position for the child.
- A choice of laterals, the shoulder section accessory and a variety of harnesses gives Mygo Seat its exceptional ability to support the upper body and head which is vital for promoting upper limb function.




- **Sunrise Medical's QUICKIE Iris** offers up to 55° rotation in space ranges, whilst being extremely easy to manoeuvre.
- The dynamic back helps maintain proper positioning and protects the wheelchair over time.
- As the lightest wheelchair base around, the QUICKIE Iris ensures that caregivers and users have the smoothest experience and best range of adjustments.
- The Iris 'Rotation in Space' technology is achieved via a rocker system, as opposed to a traditional single pivot tilt. This means that the user's centre of mass is always over the centre of the chair, so no matter how much tilt is applied the base is stable and easily manoeuvrable. The rocker system also allows for a much smoother tilting mechanism, meaning it is less likely to spark a cognitive or tonal reaction from the user when compared to an often jerky and uncontrolled single pivot tilt.



The full range of possible seat and base combinations are listed below, giving you more versatility than ever before to allow children to explore the world in the best possible way. For more information, visit [leckey.com](https://www.leckey.com)

Squiggles Seat	Mygo Seat	Kit Seat
ZIPPIE Iris	ZIPPIE Iris	ZIPPIE RS
ZIPPIE TS	ZIPPIE TS	QUICKIE Iris
	QUICKIE TS	QUICKIE SR45
	QUICKIE SR45	QUICKIE Salsa M&M2
	QUICKIE Salsa M & M2	
	ZIPPIE RS	

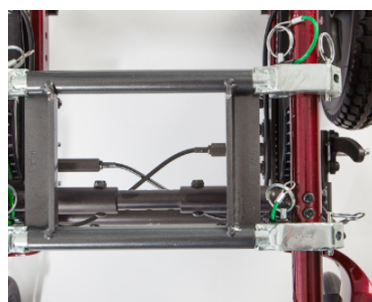

Innovative Interface

The Squiggles, Mygo and Kit seats have a new innovative interface with three different widths to fit any Sunrise Medical base quickly and easily. The aesthetically discreet combination ensures both function and style are maintained.

Quick & easy adjustability across seat & base options

Due to the fact the interface can be adjusted and widened as needed across seat sizes, you have the flexibility to try out different base options, future proofing your purchase and ensuring versatility.

Technical Information

The below table indicates the combinations that have been crash tested together, and those that have been tested as a seat and base separately. All listed combinations are commercially available in the market.

Unit	Base	Interface	Sales Code	Crash Tested? (Y/N)
SQUIGGLES SEAT	ZIPPIE Iris	Squiggles J	194-600	Y
SQUIGGLES SEAT	ZIPPIE TS	Squiggles J	194-600	Y
MYGO SEAT SIZE 1	ZIPPIE Iris	Mygo K	194-700	Y
MYGO SEAT SIZE 1	ZIPPIE TS Folding	Mygo A	117-775	Y - seat & base separately
MYGO SEAT SIZE 1	QUICKIE TS	Mygo A	117-775	Y - seat & base separately
MYGO SEAT SIZE 1	QUICKIE SR45	Mygo A	117-775	Y - seat & base separately
MYGO SEAT SIZE 2	ZIPPIE Iris	Mygo K	194-700	Y
MYGO SEAT SIZE 2	QUICKIE TS	Mygo A	117-775	Y- seat & base separately
MYGO SEAT SIZE 2	QUICKIE Salsa M	Salsa	137-860	Y - seat & base separately
MYGO SEAT SIZE 2	QUICKIE Salsa M2	Salsa	137-860	Y - seat & base separately
MYGO SEAT SIZE 2	QUICKIE SR45	Mygo A	117-775	Y - seat & base separately
MYGO SEAT SIZE 2	ZIPPIE TS Folding	Mygo A	117-775	Y - seat & base separately
MYGO SEAT SIZE 2	ZIPPIE RS	Mygo A	117-775	Y - seat & base separately
MYGO SEAT SIZE 3	QUICKIE Iris	Mygo A	117-775	Y - seat & base separately
MYGO SEAT SIZE 3	QUICKIE SR45	Mygo A	117-775	Y - seat & base separately
MYGO SEAT SIZE 3	ZIPPIE RS	Mygo A	117-775	Y - seat & base separately
MYGO SEAT SIZE 3	QUICKIE Salsa M	Salsa	137-860	Y - seat & base separately
MYGO SEAT SIZE 3	QUICKIE Salsa M2 Mini	Salsa	137-860	Y - seat & base separately
MYGO MAX	QUICKIE Salsa M	Salsa	137-873	Y - seat & base separately
MYGO MAX	QUICKIE Salsa M2	Salsa	137-873	Y - seat & base separately
KIT SEAT 1	ZIPPIE RS	Kit A	133-790	Y- seat & base separately
KIT SEAT 1	QUICKIE SR45	Kit L	194-800	Y
KIT SEAT 1	QUICKIE Iris	Kit L	194-800	Y
KIT SEAT 1	QUICKIE Salsa Mini	Salsa	137-873	Y
KIT SEAT 1	QUICKIE Salsa M2	Salsa	137-873	Y - seat & base separately
KIT SEAT 2	QUICKIE Iris	Kit L	194-800	Y
KIT SEAT 2	QUICKIE SR45	Kit L	194-800	Y
KIT SEAT 2	QUICKIE Salsa M2	Salsa	137-873	Y- seat & base separately

James Leckey Design Limited
19 Ballinderry Road
Lisburn BT28 2SA
Northern Ireland

(+44) 28 9260 0750
hello@leckey.com

leckey.com



LECKEY

James Leckey Design Limited has made every effort to ensure that the information contained in this document is accurate and up to date. At the time of going to print, the information presented is accurate to the best of the company's knowledge. James Leckey Design Limited however, cannot accept responsibility for changes to product information after going to print.