

Q400 M

NEXT GENERATION PERFORMANCE



Q-Series
Powerchairs

Advanced engineering. Customised for you.



Q400 M

ADVANCED ENGINEERING. CUSTOMISED FOR YOU.

QUICKIE, JAY, WHITMYER and SWITCH-IT bring you the most comprehensive range of premium, innovative and technically advanced powerchairs, each created individually for you whatever your clinical or lifestyle need.



THE POWER OF QUICKIE

The QUICKIE Q-series is the **most comprehensive range** of powered wheelchairs in the world! Encompassing the **most advanced drive wheel technology** and offering the very best drive traction and performance in each class, the Q-series utilises **award-winning seating** across the range to provide optimal postural support to suit the requirements of each and every user and always ensure maximum comfort.

The Q-series incorporates a comprehensive catalogue of options so that each powerchair can be tailor-made to suit the unique needs of its owner. This includes **cutting-edge electronics** and an integrated environmental interaction capability, as well as innovative special controls and mounts to provide solutions whatever the user's level of function.

Q-series - the future of powered mobility for everyone!

CLINICAL NEED



Q100

SEATING
SLING SEAT

BASES
RWD



Q200

SEATING
SLING SEAT

BASES
RWD



Q400

SEATING
SEDEO LITE

BASES
MWD



Q500

SEATING
SEDEO PRO

BASES
MWD
FWD
RWD



Q700

SEATING
SEDEO ERGO
SEDEO PRO

BASES
MWD
FWD
RWD



Q700-UP

SEATING
SEDEO ERGO UP

BASES
MWD
FWD

□ RWD: Rear-Wheel Drive
 □ MWD: Mid-Wheel Drive
 □ FWD: Front-Wheel Drive

BASE PERFORMANCE

Q400 M OVERVIEW



The Q400 M delivers a smooth, responsive and comfortable ride, all in the NEW updated chassis, providing a compact chair for indoor driving while offering great outdoor performance.

The Sedeo Lite seating system features a wide range of adjustments, all of which are quick and easy to make. With the ability to adjust to narrower and wider dimensions than ever before and with a multitude of options including JAY cushions and backs, Whitmyer headests and Switch-IT special controls, the Q400 M offers a solution for all.



ADVANCED OUTDOOR AND INDOOR PERFORMANCE

With the Q400 M base, the highly customisable SEDEO LITE seating (including JAY and WHITMYER) and a wide range of high end R-Net and SWITCH-IT control systems, you can tailor your Q400 M individually to your clinical and lifestyle needs.



PERFORMANCE BASE

Q400 M: Mid-Wheel Drive

SEATING SYSTEM

SEDEO LITE

JAY & WHITMYER

CONTROLS & CONNECTIVITY

VR2 control

R-NET controls

SWITCH-IT controls

LINK-IT mounting system

DISTINGUISHED STYLE

STYLISH DESIGN and choice of colours

DRIVE BASES

Q400 M MID-WHEEL DRIVE BASE

The new QUICKIE Q400 M base balances superb indoor manoeuvrability with outstanding outdoor performance. Its powerful motors combined with 6 wheel suspension and QUICKIE's anti-pitch technology ensure you have the most comfortable and safe driving experience.

SUPERB INDOOR MANOEUVRABILITY

Enjoy fantastic manoeuvrability from the mid-wheel drive technology with the ability to turn on its own axis, the Q400 M provides a compact turning circle of 112 cm making it easier than ever to get out of those tricky spots. With a 61 cm base width, navigating through tight doorways, hallways and cafés will be a breeze.



GREAT OUTDOOR PERFORMANCE



NEXT GENERATION SUSPENSION

The Q400 M boasts a new generation of suspension providing a responsive, smooth and comfortable ride. With the six wheel suspension it provides a soft and cushioned ride over even rough terrain and cambers.

Kerb climb up to 700 mm and have the best traction, even when ascending gradients of up to 6°.



ANTI-PITCH TECHNOLOGY

QUICKIE's anti-pitch technology eliminates pitching when decelerating, meaning you are always safe and in control. With all your weight located directly above the drive wheel to maximise traction and stability, you've never been in a safer driving position.



ENJOY YOUR TIME OUTSIDE

With reliable 2-pole motors powering speeds of 6 and 10 kph and a choice of 50, 60 or 80 Ah batteries to maximise your range, you can always enjoy your time outside with the Q400 M.



SEDEO LITE

The SEDEO LITE seating system combines an adjustable and configurable seating frame with the broadest choice of support and positioning options. Through the wide range of options available, the SEDEO LITE seating is designed to be individually tailored to your needs.

QUICK & EASY TO ADJUST

All adjustments to the SEDEO LITE can be done in a very quick and easy way thanks to its self-explanatory scales, and using just one tool.

It can be easily customised to suit individual requirements and quickly adapted should your needs change.

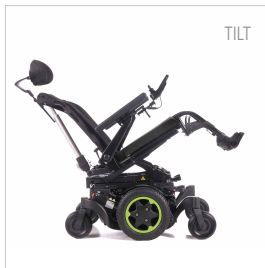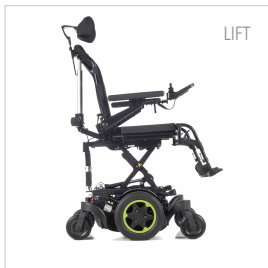

MULTI-ADJUSTABLE SEATING



SLING BACK WITH OPTIONAL COMFORT CUSHIONS

The SEDEO LITE offers height adjustable, breathable and comfortable sling backs with adjustable straps as standard. Integrated on the prescription form are NEW Sedeo Lite back cushions which offer a high level of comfort and lateral support.

JAY cushions and backs are also on the prescription form for simple ordering for clients who require a greater level of support and positioning.

16" SEAT TO FLOOR HEIGHT

With the SEDEO LITE we are able to achieve a low seat height of 415 mm with tilt or 425 mm with a lift and tilt, making accessing Wheelchair Adapted Vehicles and driving under a surface easier than ever.

SEAT RAIL TECHNOLOGY

The SEDEO LITE always comes with a seat rail making it easier, quicker and more precise than ever when mounting belts, harnesses, knee guides and hip guides. Simply attach to the rail and slide to exactly where you need them.

UNIVERSAL HEADREST

The SEDEO LITE utilises the same headrest stem as the entire range, making it easier than ever to fit the headrest that is just right for you.

FLIP BACK ARMRESTS

The easy to operate flip back armrest gives you free transfer space. You can choose from a removable or reclining version. Also available: single post removable armrest.

JAY & WHITMYER POSITIONING OPTIONS

A wide range of WHITMYER headrests and JAY cushions and backs are available with the SEDEO LITE seating systems to meet everyone's clinical needs and ensure we offer the best seating possible.



WHITMYER HEADRESTS

Created to address a range of user needs from basic support to multiple-point control. WHITMYER headrests come in a variety of pad shapes and sizes for superior head positioning and comfort.

JAY CUSHIONS & BACKS

Designed for pelvic stability, skin protection and comfort, JAY cushions and backs are low maintenance, highly customisable and are available in a large variety to meet your needs.



VR2 & R-NET

ADVANCED JOYSTICK

Every Q400M comes with a choice of VR2 control and R-Net joysticks depending on your requirements. The R-Net boasts the option of the LCD advanced joystick, which offers a large, easy-to-read HD display and simple to control paddle switches allow for intuitive operation.

CONNECT


Seamlessly take control of your environment using the integrated Infrared functionality. Control your home entertainment system, TVs, DVDs, docking stations and other home entertainment systems simply through your joystick.

BLUETOOTH & IDRIVE

Control your phone, computer, tablet and any other Bluetooth and i-devices with this useful optional extra. Connect and control up to four devices.

QUICKIE'S PATENTED ASSIGNABLE BUTTONS

Change the function on any button on your joystick to any command you desire, thanks to QUICKIE's patented assignable buttons. You can also assign a second (long press) function to buttons, giving you quicker, simpler and smoother access to functions you commonly use, without having to use your joystick.

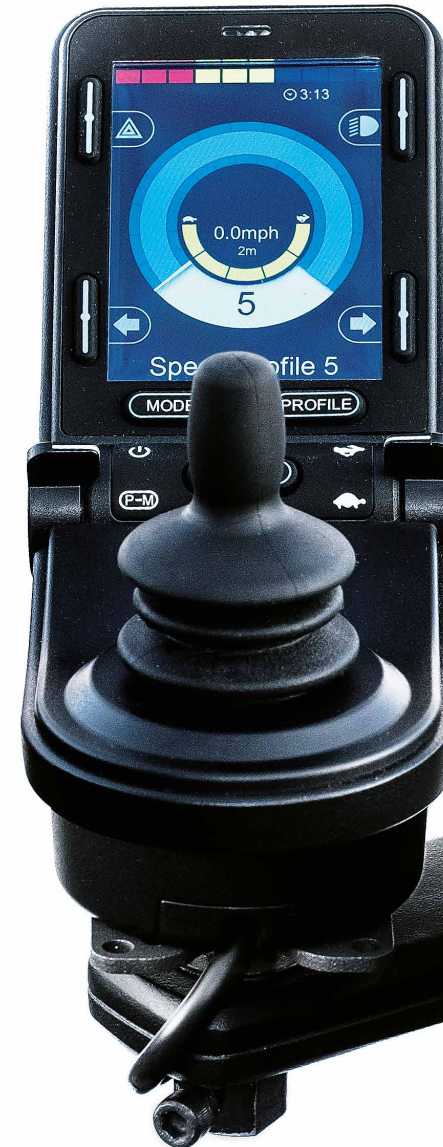


SWITCH-IT

INNOVATIVE CONTROLS

Conveniently access your 5 programmable positions with the Ctrl+5. It allows you to go straight to the function you need at the touch of a button.

SPECIAL CONTROLS


If driving with a traditional joystick isn't possible for you, we have an expansive range of alternative drive controls which can be tailored and fitted to work best for you. Allowing you to fully experience your chair's driving capabilities and benefits.

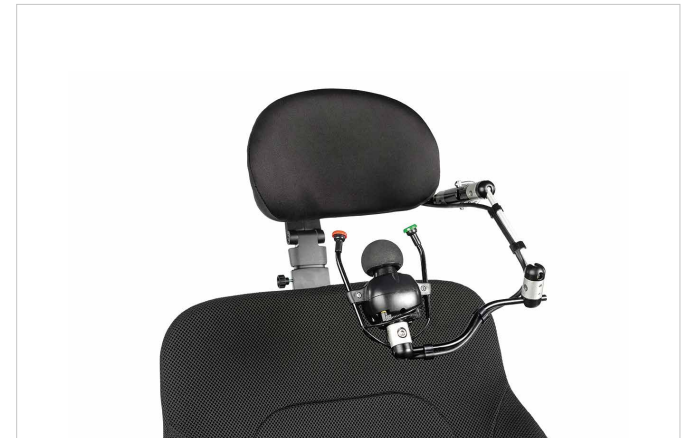

LINK-IT MOUNTING SYSTEM

Need to use additional switches or special controls? No problem, easily connect them using the universal Link-It Mounting System.

By using a combination of available modules and rods, you can create an endless number of combinations to accommodate what you need. Not only is it easy to position, Link-It also ensures the minimalistic, clean design of your wheelchair isn't compromised by additional controls. So whether it's an extra controller on the armrest, a push button on the backrest or a switch on the headrest, just link it.

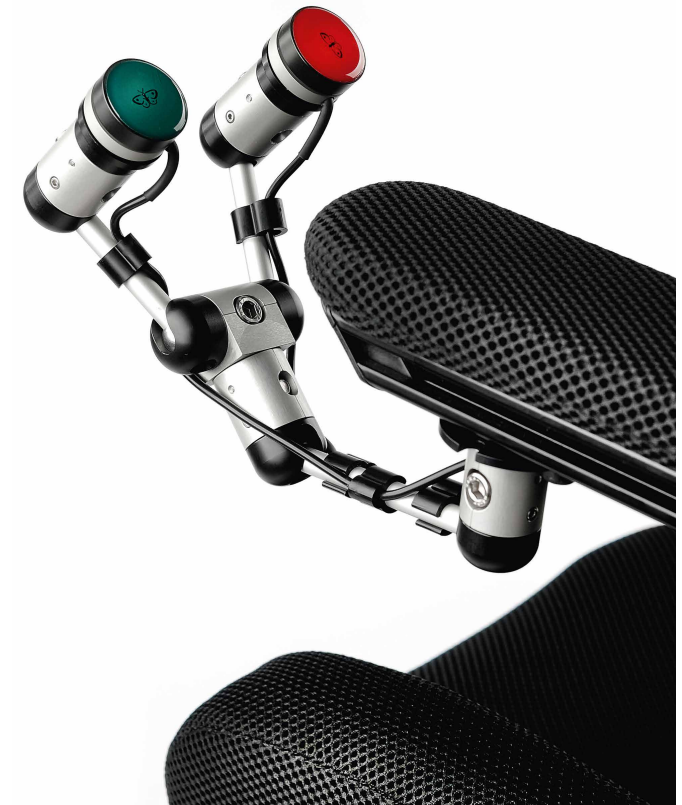






DESIGN AND APPEARANCE

DISTINGUISHED STYLE

With the Q400 M, we've taken style to the next level. With a matt-black chassis finish, hidden wiring and components, it's a sleek and contemporary design.

Fancy a change?

The Q400 M also comes with a choice of 4 contemporary colour inlays to perfectly complement its striking matt-black bodywork. Simply replace the colour inlays as desired. Styling your powered wheelchair has never been so simple.

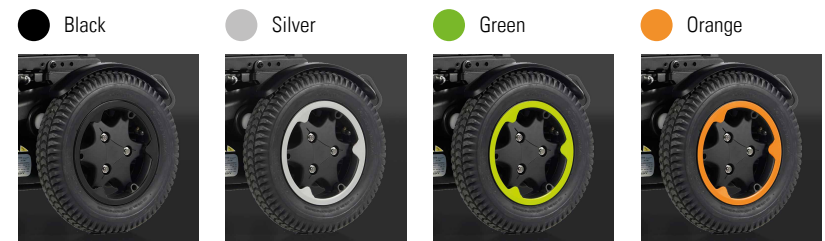

TECHNICAL DATA / Q400 M

SEDEO LITE		Q400 M	
Max. user weight	136 kg	Max. speed	10 kph
Seat width	385-510 mm	Battery size	50 Ah; 60 Ah; 80 Ah
Seat depth	410-560 mm	Electronics	VR2, R-Net
Seat height	415-490 mm	Turning radius	560-630 mm
Backrest height	460-510 mm	Max reversing width	1100 mm
Total height	Max. 1400 mm	Max. safe slope	6°
Powered recline	85°-120°	Max. kerb climbing	75 mm
Tilt angle	30°/50°	Drive wheels	13"
Seat functions	Tilt, Lift&Tilt, Recline	Castor wheel options	7"
		Colour inlay choices	4 Colour Accents
		Max. range (ISO 7176-4)	Max. 40 km
		Total width	610-620 mm
		Total length	Max. 1130 mm
		Wheelchair weight	From 130 kg



Sunrise Medical
 Thorns Road, Brierley Hill
 West Midlands, DY5 2LD
 Tel.: +44 (0) 845 605 66 88
 Fax: +44 (0) 845 605 66 89
www.SunriseMedical.co.uk

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at SunriseMedical.co.uk

