







TRY IT, CHANGE IT AND REARRANGE IT.



# FIT 2U BACK

For individuals who require a more complex seating solution, when an off-the-shelf back cannot accommodate an individual's orthopedic asymmetries, or when a molded back cannot accommodate changes. Combined with the JAY GS Cushion it makes an excellent solution

for growing and changing needs – ideally for pediatric clients.

#### MULTIDIMENSIONAL SUPPORT


To accommodate orthopedic asymmetries or a pronounced trunk rotation.

#### FLEXIBLE & QUICK TO CHANGE

Adjustments can be made upon delivery and over time.

#### **DOVETAIL INSERTS**

Add or remove interlocking inserts to create personalized contour.

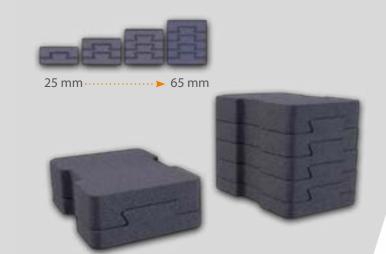





JAY Rotational Swing-Away Lateral Additional adjustment for thoracic support.

FOUR-POINT HARDWARE (FP)


JAY ROTATIONAL SWING-AWAY LATERAL

ZIPPIE IRIS° with JAY FIT 2U° BACK and JAY GS° Cushion



# **Dual Cover System**

The outer cover can easily be removed for cleaning while the inner cover keeps the inserts in place and dry.

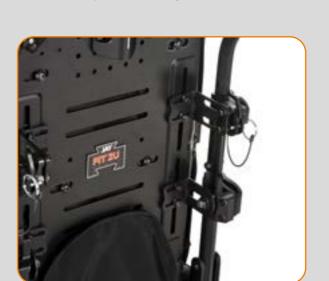

### Four-Point Hardware (FP)

Utilizes patented J3 hardware with quick release levers for fast back removal.
Two-point hardware is also available.



# Height and Width Adjustable


Five centimeters of independent height and width adjustment to handle changing postures with growth.



## **Positioning Options**

Mounting positions optimized for WHITMYER headrests with 50.8 mm lateral movement.

JAY FIT 2U Back	
Back Support Width	250 / 300 / 350 / 400 / 450 / 500 mm
Back Support Height	300 / 350 / 400 / 450 / 500 / 550 mm
Max. User Weight	350 mm wide and under: 75 kg 400 mm and above: 136 kg
Transport Weight (Without Hardware)	300 mm x 300 mm: 2.36 kg 400 mm x 500 mm: 4.13 kg
Cover Material	Stretch
Hardware	4-Point Angle: +/- 10° Width: 50 mm Depth: 35 mm or 2-Point Only available as spare part
Accessories	Large offering of JAY® and Whitmyer® seating and positioning options to meet individual needs

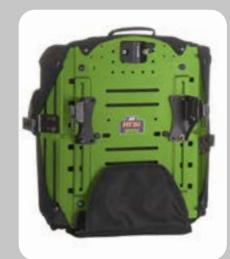

Personalize the back shell to match your chair or express your creativity.
Shown in Green Apple.



IMPORTANT NOTE:
Print processes cannot accurately represent paint finishes. Actual colors my vary.

For more information on the complete technical data, options and accessories can be found in the order form. The information is subject to change at any time. If you have any questions, please contact Sunrise Medical. If you are visually impaired, this document can be veiwed in PDF format at www.SunriseMedical.de.





Sunrise Medical GmbH Kahlbachring 2-4 69254 Malsch/HD

Tel.: +49 (0) 72 53 / 9 80-0

Fax: +49 (0) 72 53 / 9 80-222

E-Mail: export@sunrisemedical.de

www.SunriseMedical.eu