

Q300 M MINI WITH SEDEO PRO

■ IMPROVING PEOPLE'S LIVES I

INTRODUCTION



Product overview and key features

- Drive base
- Seating
- Options



Q300 M MINI



Q300 M Mini

Compact in size – BIG on performance

Available with Sedeo Lite, now adding Sedeo Pro

NEW

Q300 M Mini (Teens)

Exceptional drive performance for growing children and smaller adults

Zippie Q300 M Mini (Kids)

One of the narrowest pediatric power chairs on the market



Q300 M Mini Kids



Q300 M Mini Teens



0300 M Mini with SEDEO LITE



Q300 M Mini with SEDEO PRO

THE POWER OF QUICKIE



Part of the Quickie power range

- One modular power chair family
- Simple to prescribe and match product to clinical need
- The most efficient power chair fleet on the market



Q300 M MINI DRIVE BASE



Ultra slim 520 mm compact base (with 12" wheels)

Turning circle of just 1120 mm

• go where others can't

TRUE mid-wheel drive powerchair with all wheel suspension

- Keeps you safe and stable
- Super strong 4-pole motors
- 70 mm kerbs (14" drive wheels)
- Range of 36 km with 56 Ah batteries

Proven indoor and outdoor performance!





OUTDOOR PERFORMANCE



Anti pitch technology


You're kept safe and stable in your seat with Anti-Pitch technology, this prevents jolting you forwards when you come to a halt on a slope.

All-wheel independent suspension

Our PATENTED all-wheel independent suspension ensures all six wheels remain on the ground whilst climbing and transitioning obstacles. Any shocks or bumps are absorbed by the wheel's suspension, providing a soft and smooth ride.

Kerb climbing capability of 70 mm

With 14" drive wheels

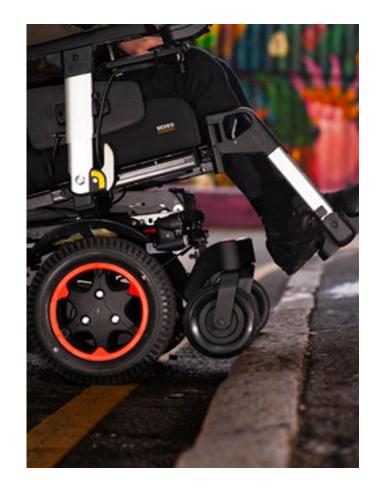

MOTORS



Sunrise Medical uses the industries best motors from Linix

Preserve power

Like a car on the motorway versus in the city, our motors are designed for optimal driving performance, so you get the power when you need it (and preserve it when you don't).

Power boost

Tricky incline? Our controller and motors work together to provide you with an extra 10 seconds of power boost when you need it the most.

Reduce thermal rollback

The larger motor cooling fins help heat to escape easier, reducing thermal rollback so you maintain a constant speed - even on slopes or uneven surfaces!

Designed to last

Unlike other powerchairs, our gears are polished metal (not plastic). Providing lower noise and higher efficiency, they're bred for tough conditions.

Options offered

•2-pole standard: 6 kph

•2-pole performance: 6 & 10 kph

•4-pole: 6 or 8 kph

• Cross-compatibility - Q300, Q400, Q500, Q700

LIGHTS




Stylish Rear/Front LED Lights with indicators

Fender lights (w/o indicators) NEW





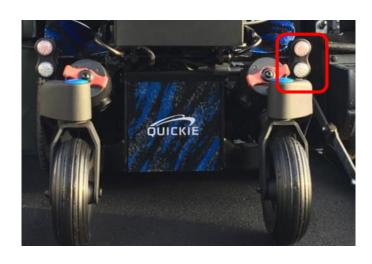

BASE OPTIONS



Base shroud colors

SHROUD:



















Color accents

ACCENT:

















SEDEO PRO SEATING



FOR CLIENTS
WITH MODERATE
TO COMPLEX
CLINICAL NEEDS

High-end fully adjustable and configurable seating frame with broadest range of options and positioning accessories

SEDEO PRO offers great flexibility to adapt to the individual user who have degenerative diseases or need to safeguard against occurrence of skeletal deformities with moderate to complex needs.

All existing Sedeo Pro options like

- Backrest cushions including flex back
- Seat cushions
- Armrests
- Legrests



POWERED SEATING OPTIONS



Low seat to floor height from 410 mm

With fixed tilt module

New lift – tilt 50°



- Also available lift tilt 30°
- Tilt in space 30°, 50° or -5° till 25°
- Helping with pressure relief by redistributing you weight from the cushion to the backrest

C-Me

- Drive elevated at 4,8 kph
- Our suspension doesn't lock-out so you get a smooth and comfortable ride



SEAT AND BACKREST FRAME



Easily adjustable and configurable seating frame

- Seat width 380 560 mm
- Additional armrest width +20 and -20 mm (for single post armrests)
- Seat depth 400 600 mm
- Integrated slotted seat rail

Backrest frames


Tool angle adjustable back, 90°-120°

- Stepless adjustment, easy to adjust with user in the chair
- Foldable backrest as standard

Powered recline, 90°-130°

Powered recline, 90°-150° with 100 mm (4") mechanical anti-shear

- Allows headrest, head and even chin control to stay in place during recline operation User adjustable manual recline back, 90°-130°
- For user who is able to shift body weight made easy with gas spring
- Operating lever for adjustment under the armrest

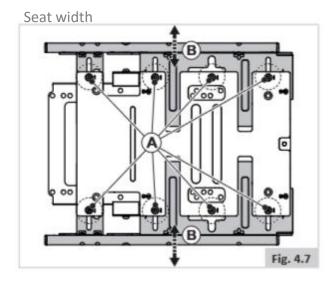





Manual recline, 90° - 130°

Powered recline, 90° - 130°

■ IMPROVING PEOPLE'S LIVES

POWER RECLINE WITH RELAX POSITION



Powered recline back 150°

- Antishear backrest
- With relax position (coming soon)
- Allows to adjust Anti-Shear recline and tilt at the same time and keeps headrest in place during recline.





SEAT AND BACK CUSHIONS



A – SEDEO PRO Back Cushion 🕥	
OMP030109	SEDEO PRO, Back Cushion, Profiled
OMP030113	SEDEO PRO, Back Cushion, Profiled, Adjustable Straps
OMP030114	SEDEO PRO, Back Cushion, Deep Contour, Passive
OMP030115	SEDEO PRO, Back Cushion, Flexible, Deep Contour, Passive, Adjustable Straps Only Available For NL
OMP030116	SEDEO PRO, Back Cushion, Flexible, Deep Contour, Active, Adjustable Straps
OMP030117	SEDEO PRO, Back Cushion, Flexible, Extra Deep Contour, Active, Adjustable Straps
C – SEDEO Seat Cushions 🕥	

SEDEO, Seat Cushion, Stabilizer, Height 90 mm


QMP020078

ARMRESTS



4 types of armrests

Single post flip back armrest Single post removable armrest Reclining flip back armrest Single post SEDEO LITE armrest

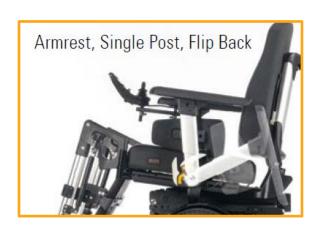







NEW ARMREST



Multi adjustable armrest NEW



- 2 extra hinges
- Height, width and angle adjustable
- Wrap around the user keeping the wheelchair's width as slim as possible





LEGRESTS



Centre mount legrest


- Angle adjustable
- Smaller footplate for small Q300 M Mini base


Power center mount legrest (small version – no tracks at the sides)

- Length compensating (260mm), user is kept in ideal position
- Separate footplates mounted inwards
- Option 'footplate to floor'
- Operate the two actuators separately
- Lower the footplate to the floor for easy forward transfers

Swing away legrests



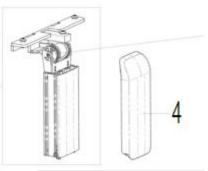

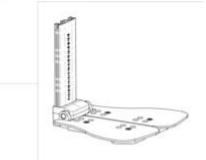





SUMMARY



Design Objective

Complete the Q300 M Mini Family by Integrating Sedeo Pro Seating onto the Q300 M Mini Base

Product Summary

- Existing Q300M Mini base now available with proven Sedeo Pro
- Improvements
- Now available with Lift tilt 50°
- Powered anti shear back recline back 150° with relax position (coming soon)
- Fender lights (w/o indicators)

